

At the high dose (1.2 µg), the glycoPEGylated rFSH had somewhat higher *in vivo* activity than the unPEGylated rFSH.

G-CSF

28. GlycoPEGylation of G-CSF produced in CHO cells

Preparation of Asialo-Granulocyte-Colony Stimulation Factor (G-CSF). G-CSF produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂ and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(*p*-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel were run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(α2,3)-Sialyl-PEG. Desialylated G-CSF was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis

according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(alpha2,8)-Sialyl-PEG. G-CSF produced in CHO cells, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of G-CSF-(alpha2,6)-Sialyl-PEG. G-CSF, containing only O-linked GalNAc, is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

G-CSF produced in CHO cells was treated with Arthrobacter sialidase and was then purified by size exclusion on Superdex75 and was treated with ST3Gal1 or ST3 Gal2 and then with CMP-SA-PEG 20Kda. The resulting molecule was purified by ion exchange and

gel filtration and analysis by SDS PAGE demonstrated that the PEGylation was complete. This is the first demonstration of glycoPEGylation of an O-linked glycan.

Glucocerebrosidase

29. Glucocerebrosidase-mannose-6-phosphate produced in CHO cells

This example sets forth the procedure to glycoconjugate mannose-6-phosphate to a peptide produced in CHO cells such as glucocerebrosidase.

Preparation of asialo-glucoceramide. Glucocerebrosidase produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL sialidase-agarose conjugate for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer, and once with 0.2 mL of the Tris-EDTA buffer. All supernatants are pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure

1). Asialo-glucocerebrosidase from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using

SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Glucocerebrosidase-SA-linker-Mannose-6-phosphate (procedure

2). Glucocerebrosidase, produced in CHO but incompletely sialylated, is dissolved at 2.5
5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated
with 1 mM CMP-sialic acid-linker-Man-6-phosphate and 0.1 U/mL of ST3Gal3 at 32°C for 2
days. To monitor the incorporation of sialic acid-linker-Man-6-phosphate, a small aliquot of
the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the
peptide is separated from the free label by gel filtration on a Toso Haas TSK-Gel-3000
10 analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the
peptide is quantitated using an in-line fluorescent detector. When the reaction is complete,
the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using
PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the
reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and
15 reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by
MALDI-TOF MS.

30. Glucocerebrosidase-transferrin

This example sets forth the procedures for the glycoconjugation of proteins, and in
20 particular, transferrin is glycoconjugated to glucocerebrosidase. The GlcNAc-ASN structures
are created on glucoceraminidase, and Transferrin-SA-Linker-Gal-UDP is conjugated to
GNDF GlcNAc-ASN structures using galactosyltransferase.

Preparation of GlcNAc-glucocerebrosidase (Cerezyme™). Cerezyme™

(glucocerebrosidase) produced in CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM
25 Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL Endo-H-agarose conjugate
for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with
the appropriate buffer and a IEF gel and SDS-PAGE performed according to Invitrogen
procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The
beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and
30 once with 0.2 mL of the Tris-EDTA buffer and all supernatants are pooled. The supernatant
is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice

more against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Transferrin-SA-Linker-Gal-glucocerebrosidase. Transferrin-SA-Linker-Gal-UDP from above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 2.5 mg/mL GlcNAc-glucocerebrosidase and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of glucocerebrosidase, the peptide is separated by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1) and the product detected by UV absorption. The reaction mixture is then purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

GM-CSF

31. Generation and PEGylation of GlcNAc-ASN Structures: GM-CSF produced in *Saccharomyces*

This example sets forth the preparation of Tissue-type Activator with PEGylated GlcNAc-Asn structures.

Recombinant GM-CSF expressed in yeast is expected to contain 2 N-linked and 2 O-linked glycans. The N-linked glycans should be of the branched mannan type. This recombinant glycoprotein is treated with an endoglycosidase from the group consisting of endoglycosidase H, endoglycosidase-F1, endoglycosidase-F2, endoglycosidase-F3, endoglycosidase-M either alone or in combination with mannosidases I, II and III to generate GlcNAc nubs on the asparagine (Asn) residues on the peptide/protein backbone.

The GlcNAc-Asn structures on the peptide/protein backbone is then be modified with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case the galactose-PEG is the terminal residue.

In the second case the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment the GlcNAc-Asn structures on the peptide/protein backbone can be galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an α 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

Herceptin™

32. Glycoconjugation of mithramycin to Herceptin™

This example sets forth the procedures to glycoconjugate a small molecule, such as mithramycin to Fc region glycans of an antibody molecule produced in mammalian cells. Here, the antibody Herceptin™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

Preparation of Herceptin™-Gal-linker-mithramycin. Herceptin™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-mithramycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the mithramycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Interferon α and Interferon β 33. GlycoPEGylation of Proteins expressed in Mammalian or Insect Systems:
EPO, Interferon α and Interferon β

This example sets forth the preparation of PEGylated peptides that are expressed in
5 mammalian and insect systems.

Preparation of acceptor from mammalian expression systems. The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have glycans terminating in galactose. Most peptides from mammalian expression systems will have terminal sialic acid that first
needs to be removed.

10 **Sialidase digestion.** The peptide is desialylated using a sialidase. A typical procedure involves incubating a 1 mg/mL solution of the peptide in Tris-buffered saline, pH 7.2, with 5 mM CaCl_2 added, with 0.2 U/mL immobilized sialidase from *Vibrio cholera* (Calbiochem) at 32°C for 24 hours. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The resin is then removed by centrifugation or
15 filtration, and then washed to recover entrapped peptide. At this point, EDTA may be added to the solution to inhibit any sialidase that has leached from the resin.

Preparation from insect expression systems. EPO, interferon-alpha, and interferon-beta may also be expressed in non-mammalian systems such as yeast, plants, or insect cells. The peptides to be glycoPEGylated using CMP-sialic acid PEG need to have
20 glycans terminating in galactose. The majority of the N-glycans on peptides expressed in insect cells, for example, are the trimannosyl core. These glycans are first built out to glycans terminating in galactose before they are acceptors for sialyltransferase.

Building acceptor glycans from trimannosyl core. Peptide (1 mg/mL) in Tris-buffered saline, pH 7.2, containing 5 mM MnCl_2 , 5 mM UDP-glcNAc, 0.05 U/mL
25 GLCNACT I, 0.05 U/mL GLCNACT II, is incubated at 32°C for 24 hours or until the reaction is substantially complete. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. After buffer exchange to remove UDP and other small molecules, UDP-galactose and MnCl_2 are each added to 5 mM, galactosyltransferase is added to 0.05 U/mL, and is incubated at 32°C for 24H or until the reaction is substantially
30 complete. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The peptides are then ready for glycoPEGylation.

Building O-linked glycans. A similar strategy may be employed for interferon alpha to produce enzymatically the desired O-glycan Gal-GalNAc. If necessary, GalNAc linked to serine or threonine can be added to the peptide using appropriate peptide GalNAc transferases (e.g. GalNAc T1, GalNAc T2, T3, T4, etc.) and UDP-GalNAc. Also, if needed, galactose can be added using galactosyltransferase and UDP-galactose.

GlycoPEGylation using sialyltransferase. The glycopeptides (1 mg/mL) bearing terminal galactose in Tris buffered saline + 0.02% sodium azide are incubated with CMP-SA-PEG (0.75 mM) and 0.4 U/mL sialyltransferase (ST3Gal3 or ST3Gal4 for N-glycans on EPO and interferon beta; ST3Gal4, or ST3Gal1 for O-glycans on interferon alpha) at 32°C for 24 hours. Other transferases that may work include the 2,6 sialyltransferase from *Photobacterium damsella*. The acceptor peptide concentration is most preferably in the range of 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA-PEG should be sufficient for there to be excess over the available sites, but not so high as to cause peptide solubility problems due to the PEG, and may range from 50 µM up to 5 mM, and the temperature may range from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH.

34. GlycoPEGylation of Interferon α produced in CHO cells

Preparation of Asialo-Interferon α . Interferon alpha produced from CHO cells is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, 5 mM CaCl₂ and concentrated to 500 µL in a Centricon Plus 20 centrifugal filter. The solution is incubated with 300 mU/mL Neuraminidase II (*Vibrio cholerae*) for 16 hours at 32 °C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed. The reaction mixture is then added to prewashed N-(p-aminophenyl)oxamic acid-agarose conjugate (800 µL/mL reaction volume) and the washed beads gently rotated for 24 hours at 4 °C. The mixture is centrifuged at 10,000 rpm and the supernatant was collected. The beads are washed 3 times with Tris-EDTA buffer, once with 0.4 mL Tris-EDTA buffer and once with 0.2 mL of the Tris-EDTA buffer and all supernatants were pooled. The supernatant is dialyzed at 4 °C against 50 mM Tris-HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris-HCl pH 7.4, 1 M

NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples of native and desialylated G-CSF are dialyzed against water and analyzed by MALDI-TOF MS.

5 **Preparation of Interferon-alpha-(alpha2,3)-Sialyl-PEG.** Desialylated interferon-alpha is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST3Gal1 at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide
10 is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF
15 analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and desialylated Interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Interferon-alpha-(alpha2,8)-Sialyl-PEG. Interferon-alpha produced in CHO, which contains an alpha2,3-sialylated O-linked glycan, is dissolved at 2.5
20 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of CST-II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction has CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free
25 label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis
30 according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon-alpha are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Interferon- α -(α 2,6)-Sialyl-PEG. Interferon- α , containing only O-linked GalNAc, was dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM CMP-sialic acid-PEG and 0.1 U/mL of ST6GalNAcI or II at 32°C for 2 days. To monitor the incorporation of sialic acid-PEG, a small aliquot of the reaction had CMP-SA-PEG-fluorescent ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The fluorescent label incorporation into the peptide is quantitated using an in-line fluorescent detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples of native and PEGylated interferon- α are dialyzed against water and analyzed by MALDI-TOF MS.

35. GlycoPEGylation of Interferon- β -1a with PEG (10 kDa) and PEG (20 kDa)

This example illustrates a procedure PEGylate Interferon- β with either PEG (10 kDa) or PEG (20 kDa).

Briefly, Interferon- β -1a (INF- β) was obtained from Biogen (Avonex™). The INF- β was first purified by Superdex-75 chromatography. The INF- β was then desialylated with *Vibrio cholerae* sialidase. The INF- β was then PEGylated with SA-PEG (10 kDa) or SA-PEG (20 kDa) and purified with Superdex-200 chromatography.

Superdex-75 chromatography purification. INF- β (150 μ g) was applied to a Superdex-75 column (Amersham Biosciences, Arlington Heights, IL) and eluted with PBS with 0.5 M NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol. The eluant was monitored for absorbance at 280 nm (Figure 172A and 172B) and fractions were collected. Peaks 4 and 5 were pooled, concentrated in an Amicon Ultra 15 spin filter (Millipore, Billerica, MA), and the buffer was exchanged to TBS with 5 mM CaCl₂, 0.02% Tween-20, 20 mM histidine and 10% glycerol.

Sialidase Reaction. The INF- β was then desialylated with *Vibrio cholera* sialidase (70 mU/ml, CALBIOCHEM®, EMD Biosciences, Inc., San Diego, CA) on agarose in TBS

with 5 mM CaCl₂, 0.02% Tween-20, 20 mM histidine and 10% glycerol. The reaction was carried out at 32°C for 18 hours. The INF- β was removed from the agarose with a 0.22 μ m Spin-X™ filter (Corning Technology, Inc., Norcross, GA). Figure 173A depicts the MALDI analysis of glycans released from native INF- β . The native INF- β has many glycoforms containing terminal sialic acid moieties. Figure 173B depicts the MALDI analysis of glycans released from desialylated INF- β . The desialylated INF- β has primarily one glycoform which is bi-antennary with terminal galactose moieties.

Lectin Dot-Blot Analysis of Sialylation. Samples of the INF- β from the desialidase reaction were dot-blotted onto nitrocellulose and then blocked with Tris buffered saline (TBS: 0.05M Tris, 0.15M NaCl, pH 7.5) and DIG kit (glycan differentiation kit available from Roche #1 210 238) blocking buffer. Some of the blots were incubated with *Maackia amurensis* agglutinin (MAA) labeled with digoxigenin (DIG) (Roche Applied Science, Indianapolis, IL) to detect α 2,3-sialylation of INF- β . These blots were washed with TBS then incubated with anti-digoxigenin antibody labeled with alkaline phosphatase, then washed again with TBS and developed with NBT/X-phosphate solution, wherein NBT is 4-nitro blue tetrazolium chloride and X-phosphate is 5-bromo-4-chloro-3-indoyl phosphate. The left side of Figure 174 depicts the results of the MAA blot of INF- β after the desialylation reaction. The INF- β is partially disialylated, as indicated by the decrease in dot development as compared to native INF- β in the desialylated samples.

Other blots were incubated with *Erthrina cristagalli* lectin (ECL) labeled with biotin (Vector Laboratories, Burlingame, CA) to detect exposed galactose residues on INF- β . After incubation with 2.5 μ g/ml ECL, the blots were washed in TBS and incubated with streptavidin labeled with alkaline phosphatase. The blots were then washed again and developed. The right side of Figure 174 depicts the ECL blot after development. The increased intensity of the dot of desialylated INF- β as compared to the native INF- β indicate more exposed galactose moieties and therefore extensive desialylation.

PEGylation of Desialylated INF- β with SA-PEG (10 kDa). Desialylated INF- β (0.05 mg/ml) was PEGylated with ST3Gal3 (50 mU/ml) and CMP-SA-PEG (10 kDa) (250 μ M) in an appropriate buffer of TBS + 5 mM CaCl₂, 0.02% Tween 20, 20 mM histidine, 10%

glycerol for 50 hours at 32°C. Figure 175 depicts the SDS-PAGE analysis of the reaction products showing PEGylated INF- β at approximately 98 kDa.

PEGylation of Desialylated INF- β with SA-PEG (20 kDa). Desialylated INF- β (0.5 mg/ml) was PEGylated with ST3Gal3 (170 mU/ml) and CMP-SA-PEG (20 kDa) in an appropriate buffer of TBS + 5 mM CaCl₂, 0.02% Tween 20, 20 mM histidine, 10% glycerol for 50 hours at 32°C. Figure 176 depicts the SDS-PAGE analysis the products of the PEGylation reaction. The PEGylated INF- β has many higher molecular weight bands not found in the unmodified INF- β indicating extensive PEGylation.

Superdex-200 Purification of INF- β PEGylated with PEG (10 kDa). The products of the PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1ml/min and 30 cm/hr flow. The eluant was monitored for absorbance at 280 nm (Figure 177) and fractions were collected. Peaks 3 and 4 were pooled and concentrated in an Amicon Ultra 15 spin filter.

Bioassay of INF- β PEGylated with PEG (10 kDa).

The test is inhibition of the proliferation of the lung carcinoma cell line, A549. The A549 cell line are lung carcinoma adherent cells growing in RPMI + 10% FBS at 37°C 5% CO₂. They can be obtained from ATCC # CCL-185. Wash the cells with 10 ml of PBS and remove the PBS. Add 5 ml of trypsin, incubate for 5 minutes at room temperature or 2 minutes at 37°C. When the cells are detached resuspend into 25 ml of media and count the cells. Dilute the cells at a concentration of 10000 cells/ml and add 200 μ l / well (96 wells plate). Incubate for 4 hours at 37°C 5% CO₂. Prepare 1 ml of IFN β at a concentration of 0.1 μ g/ml. Filter it under the hood with a 0.2 μ m filter. Add 100 μ l per well (8 replicates = 1 lane). Incubate for 3 days (do not let the cells go to confluence). Remove 200 μ l of media (only 100 μ l per well left). Add 25 μ l of MTT (Sigma) (5 mg/ml filtered 0.22 μ m). Incubate for 4 hours at 37°C and 5% CO₂. Aspirate the media gently and add 100 μ l of a mixture of isopropanol (100 ml and 6N HCl. Aspirate up and down to homogenize the crystal violet. Read OD 570nm (remove the background at 630 or 690 nm).

Figure 178 depicts the results of the bioassay of the peaks containing INF- β PEGylated with PEG (10 kDa) as eluted from the Superdex-200 column.

Superdex-200 Purification of INF- β PEGylated with PEG (20 kDa). The products of the PEG (20 kDa) PEGylation reaction were separated on a Superdex-200 column (Amersham Biosciences, Arlington Heights, IL) in PBS with 0.5 NaCl, 0.02 Tween-20, 20 mM histidine and 10% glycerol at 1 ml/min flow. The eluant was monitored for absorbance at 280 nm (Figure 179) and fractions were collected. Peak 3 contained most of the INF- β PEGylated with PEG (20 kDa).

Endotoxin test of INF- β PEGylated with PEG (20 kDa).

Limulus Lysate Test was performed, BioWhittaker # 50-647U

Table 24. Results of the endotoxin test of INF- β PEGylated with PEG (20 kDa).

	Concentration		
INF- β with PEG (20 kDa)	10 EU/ml	0.06 mg/ml	0.16 EU/ μ g
INF- β with PEG (20 kDa)	1 EU/ml	0.07 mg/ml	0.014 EU/ μ g
Native INF- β	40 EU/ml	0.1 mg/ml	0.4 EU/ μ g

RemicadeTM

36. GlycoPEGylation of RemicadeTM antibody

This example sets forth the procedure to glycoPEGylate a recombinant antibody molecule by introducing PEG molecules to the Fc region glycans. Here RemicadeTM, a TNF-R:IgG Fc region fusion protein, is the exemplary peptide.

Preparation of RemicadeTM-Gal-PEG (10 kDa). RemicadeTM is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-PEG (10 kDa) and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the PEG in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer

exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Rituxan™

37. Glycoconjugation of geldanamycin to Rituxan™

This example sets forth the glycoconjugation of a small molecule, such as geldanamycin, to the Fc region glycans of an antibody produced in CHO cells, such as Rituxan™. Here, the antibody Rituxan™ is used, but one of skill in the art will appreciate that the method can be used with many other antibodies.

Preparation of Rituxan™-Gal-linker-geldanamycin. Rituxan™ is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 1 mM UDP-galactose-linker-geldanamycin and 0.1 U/mL of galactosyltransferase at 32°C for 2 days to introduce the geldanamycin in the Fc region glycans. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

When the reaction is complete, the reaction mixture is purified using a Toso Haas TSK-Gel-3000 preparative column using PBS buffer (pH 7.1) and collecting fractions based on UV absorption. The fractions containing product are combined, concentrated, buffer exchanged and then freeze-dried. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Rnase

38. Remodeling high mannose N-glycans to hybrid and complex N-glycans: Bovine pancreatic RNase

This example sets forth the preparation of bovine pancreas RNase with hybrid or complex N-glycans. The high mannose N-linked glycans of the RNase are enzymatically

digested and elaborated to create hybrid N-linked glycans. Additionally, the high mannose N-linked glycans of the RNase are enzymatically digested and elaborated to create complex N-linked glycans.

High mannose structures of *N*-linked oligosaccharides in glycopeptides can be modified to hybrid or complex forms using the combination of α -mannosidases and glycosyltransferases. This example summarizes the results in such efforts using a simple *N*-Glycan as a model substrate.

Ribonuclease B (RNaseB) purified from bovine pancreas (Sigma) is a glycopeptide consisting of 124 amino acid residues. It has a single potential *N*-glycosylation site modified with high mannose structures. Due to its simplicity and low molecular weight (13.7 kDa to 15.5 kDa), ribonuclease B is a good candidate to demonstrate the feasibility of the *N*-Glycan remodeling from high mannose structures to hybrid or complex *N*-linked oligosaccharides. The MALDI-TOF spectrum of RNaseB (Figure 180A) and HPLC profile for the oligosaccharides cleaved from RNaseB by *N*-Glycanase (Figure 180B) indicated that, other than a small portion of the non-modified peptide, the majority of *N*-glycosylation sites of the peptide are modified with high mannose oligosaccharides consisting of 5 to 9 mannose residues.

Conversion of high mannose *N*-Glycans to hybrid *N*-Glycans. High mannose *N*-Glycans were converted to hybrid *N*-Glycans using the combination of α 1,2-mannosidase, GlcNAcT-I (β -1,2-*N*-acetyl glucosaminyl transferase), GalT-I (β 1,4-galactosyltransferase) and α 2,3-sialyltransferase /or α 2,6-sialyltransferase as shown in Figure 181.

As an example, high mannose structures in RNaseB were successfully converted to hybrid structures.

Man₅GlcNAc₂-R was obtained from Man₅₋₉GlcNAc₂-R catalyzed by a single α 1,2-mannosidase cloned from *Trichoderma reesei* (Figure 182). RNase B (1 g, about 67 μ mol) was incubated at 30°C for 45 hr with 15 mU of the recombinant *T. reesei* α 1,2-mannosidase in MES buffer (50 mM, pH 6.5) in a total volume of 10 mL. Man₆₋₉GlcNAc₂-protein structures have been successfully converted to Man₅GlcNAc₂-protein with high efficiency by the recombinant mannosidase.

Alternately, Man₅GlcNAc₂-R was obtained from Man_{5,9}GlcNAc₂-R catalyzed by a single α 1,2-mannosidase purified from *Aspergillus saitoi* (Figure 183). RNase B (40 μ g, about 2.7 nmol) was incubated at 37°C for 42.5 hr with 25 μ U of the commercial *A. saitoi* α 1,2-mannosidase (Glyko or CalBioChem) in NaOAc buffer (100 mM, pH 5.0) in a total
5 volume of 20 μ l. Man_{6,9}GlcNAc₂-protein structures were successfully converted to Man₅GlcNAc₂-protein by the commercially available mannosidase. However, a new peak corresponding to the GlcNAc-protein appears in the spectrum, indicating the possible contamination of endoglycosidase H in the preparation. Although several mammalian alpha-mannosidases were required to achieve this step, the fungal α 1,2-mannosidase was very
10 efficient to remove all α 1,2-linked mannose residues.

GlcNAcT-I then added a GlcNAc residue to the Man₅GlcNAc₂-R (Figure 184). The reaction mixture after the *T. reesei* α 1,2-mannosidase reaction containing RNase B (600 μ g, about 40 nmol) was incubated with non-purified recombinant GlcNAcT-I (34 mU) in MES buffer (50 mM, pH 6.5) containing MnCl₂ (20 mM) and UDP-GlcNAc (5 mM) in a total
15 volume of 400 μ l. at 37°C for 42 hr. A GlcNAc residue was quantitatively added to Man₅GlcNAc₂-protein by the recombinant GlcNAcT-I.

A Gal residue was then added using GalT 1 (Figure 185). The reaction mixture after the GnT-I reaction containing RNase B (120 μ g, about 8 nmol) was incubated at 37°C for 20 hr with 3.3 mU of the recombinant GalT-1 in Tris-HCl buffer (100 mM, pH 7.3) containing
20 UDP-Gal (7.5 mM) and MnCl₂ (20 mM) in a total volume of 100 μ l. A Gal residue was added to about 98% of the GlcNAc-Man₅GlcNAc₂-protein by the recombinant GalT 1.

The next step was the addition of a sialic acid using an α 2,3-sialyltransferase or an α 2,6-sialyltransferase (Figure 186). As an example, ST3Gal III, an α 2,3-sialyltransferase was used. The reaction mixture after the GalT-1 reaction containing RNase B (13 μ g, about 0.87
25 nmol) was incubated at 37°C for 16 hr with 8.9 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing CMP-Sialic acid (5 mM) and MnCl₂ (20 mM) in a total volume of 20 μ l. A sialic acid residue was added to about 90% of the Gal-GlcNAc-Man₅GlcNAc₂-protein by recombinant ST3Gal III using CMP-SA as the donor. The yield can be further improved by adjusting the reaction conditions.

For convenience, no purification or dialysis step was required after each reaction described above. More interesting, GalT 1 and ST3Gal III can be combined in a one-pot reaction. Similar yields were obtained as compared with the separate reactions. The reaction mixture after the GlcNAcT-I reaction containing RNase B (60 μ g, about 4 nmol) was
5 incubated at 37°C for 20 hr with 1.7 mU of recombinant GalT 1, 9.8 mU of recombinant ST3Gal III in Tris-HCl buffer (100 mM, pH 7.3) containing UDP-Gal (7.5 mM), CMP-sialic acid (5 mM) and MnCl₂ (20 mM) in a total volume of 60 μ l.

As shown in Figure 187, SA-PEG (10 kDa) was successfully added to the RNaseB. The reaction mixture after the GalT-1 reaction containing RNase B (6.7 μ g, about 0.45 nmol)
10 was dialyzed against H₂O for 1 hour at room temperature and incubated at 37°C for 15.5 hours with 55 mU of the recombinant ST3Gal III in Tris-HCl buffer (50 mM, pH 7.3) containing CMP-SA-PEG (10 kDa) (0.25 mM) and MnCl₂ (20 mM) in a total volume of 20 μ l. PEG-modified sialic acid residues were successfully added to the Gal-GlcNAc-Man₅GlcNAc₂-peptide by the recombinant ST3Gal III. The yield can be further improved by
15 adjusting the reaction conditions.

Conversion of high mannose N-Glycans to complex N-Glycans. To achieve this conversion, a GlcNAc β 1,2Man₃GlcNAc₂-peptide intermediate is obtained. As shown in Figure 188, there are at least four feasible routes to carry out the reaction from Man₅GlcNAc₂-peptide to this intermediate:

20 **Route I:** The Man₅GlcNAc₂-peptide produced by the fungal α 1,2 mannosidase is a substrate of GlcNAc transferase I (GlcNAcT-I, enzyme 2) which adds one GlcNAc. The terminal α 1,3- and α 1,6-linked mannose residues of GlcNAcMan₅GlcNAc₂-peptide is removed by Golgi α -mannosidase II (ManII, enzyme 5). This route is a part of the natural pathway for the processing of N-linked oligosaccharides carried out in higher organisms.

25 **Route II:** Two mannose residues are first removed by an α -mannosidase (enzyme 6), then a GlcNAc is added by GlcNAcT-I (enzyme 2). Other than its natural acceptor Man₅GlcNAc₂-R, GlcNAcT-I can also recognize Man₃GlcNAc₂-R as its substrate and add one GlcNAc to the mannose core structure to form GlcNAcMan₃GlcNAc₂-peptide.

Route III: The α 1,6-linked mannose is removed by an α 1,6-mannosidase, followed
30 by the addition of GlcNAc by GlcNAcT-I and removal of the terminal α 1,3-linked mannose

by an α 1,3-mannosidase. From the experimental data obtained, GlcNAcT-I can recognize this $\text{Man}_4\text{GlcNAc}_2$ -peptide as acceptor and add one GlcNAc residue to form $\text{GlcNAcMan}_4\text{GlcNAc}_2$ -peptide.

Route IV: Similar to Route III, α 1,3-linked mannose is removed by an α 1,3-mannosidase, followed by GlcNAcT-I reaction. Then the terminal α 1,6-linked mannose can be removed by an α 1,6-mannosidase.

After the function of GlcNAcT-I (responsible for the addition of the GlcNAc β 1,2-linked to the α 1,3-mannose on the mannose core) and GlcNAcT-II (responsible for the addition of a second GlcNAc β 1,2-linked to the α 1,6-mannose on the mannose core), the $\text{GlcNAc}_2\text{Man}_3\text{GlcNAc}_2$ -peptide can be processed by GalT 1 and sialyltransferase to form bi-antennary complex N- Glycans. Other GlcNAc transferases such as GlcNAcT-IV, GlcNAcT-V, and/or GlcNAcT-VI (Figure 188 and Figure 189) can also glycosylate the $\text{GlcNAc}_2\text{Man}_3\text{GlcNAc}_2$ -peptide. Additional glycosylation by the GalT 1 and sialyltransferases will form multi-antennary complex N-glycans. The enzyme GlcNAcT-III catalyzes the insertion of a bisecting GlcNAc, thus preventing the actions of ManII and subsequent action of transferases GlcNAcT-II, GlcNAcT-IV and GlcNAcT-V.

Tissue-Type Plasminogen Activator (TPA)

39. Fucosylation of TPA to create Sialyl Lewis X

This example sets forth the preparation of Tissue Tissue-type Plasminogen Activator (TPA) with N-linked sialyl Lewis X antigen.

Sialylation. TPA expressed in mammalian cells will often contain a majority of the glycans terminating in sialic acid, but to ensure complete sialylation, it would be beneficial to first perform an *in vitro* sialylation. TPA in a suitable buffer (most preferably between pH 5.5 and 9, for example Tris buffered saline, pH 7.2) is incubated with CMP sialic acid and sialyltransferase for a time sufficient to convert any glycans lacking sialic acid to sialylated species. Typical conditions would be 1 mg/mL TPA, 3 mM CMP sialic acid, 0.02 U/mL ST3Gal3, 32°C for 24 hours. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of CMP-SA should be

sufficient for there to be excess over the available sites, and might range from 50 μ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other sialyltransferases that may be capable of adding sialic acid in 2,3 linkage include ST3Gal4; microbial transferases could also be used.

Fucosylation. Typical conditions for fucosylation would be 1 mg/mL TPA, 3 mM GDP-fucose, 0.02 U/mL FTVI, 5 mM $MnCl_2$, 32°C for 24H in Tris buffered saline. Microbial growth can be halted either by sterile filtration or the inclusion of 0.02% sodium azide. The TPA concentration is most preferably in the range 0.1 mg/mL up to the solubility limit of the peptide. The concentration of GDP-fucose should be sufficient for there to be excess over the available sites, and might range from 50 μ M up to 50 mM, and the temperature from 2°C up to 40°C. The time required for complete reaction will depend on the temperature, the relative amounts of enzyme to acceptor substrate, the donor substrate concentration, and the pH. Other fucosyltransferases that may be capable of making sialyl Lewis x include FTVII, FTV, FTIII, as well as microbial transferases could also be used.

40. Trimming of high mannose to tri-mannose core structure: Tissue-type Plasminogen Activator produced in CHO

This example sets forth the preparation of Tissue-type Plasminogen Activator with a trimannose core by trimming back from a high mannose glycan.

Tissue-type plasminogen activator (TPA) is currently produced in Chinese Hamster Ovary (CHO) cells and contains a low amount of high mannose N-linked oligosaccharide. The mannoses can be trimmed down using a variety of the specific mannosidases. The first step is to generate $Man_5GlcNAc_2(Fuc0-1)$ from $Man_9GlcNAc_2(Fuc0-1)$. This can be done using mannosidase I. Then either $GlcNAcT1$ ($GlcNAc$ transferase I) is used to make $GlcNAc1Man_5GlcNAc_2(Fuc0-1)$ or Mannosidase III is used to make $Man_3GlcNAc_2(Fuc0-1)$. From $Man_3GlcNAc_2(Fuc0-1)$, $GlcNAc1Man_3GlcNAc_2(Fuc0-1)$ can be produced using $GlcNAcT1$ or from $GlcNAc1Man_5GlcNAc_2(Fuc0-1)$, $GlcNAc1Man_3GlcNAc_2(Fuc0-1)$ can be produced using Mannosidase II. $GlcNAc1Man_3GlcNAc_2(Fuc0-1)$ is then converted into $GlcNAc2Man_3GlcNAc_2(Fuc0-1)$ using $GlcNAcTransferase II$ ($GlcNAcTII$). The two

terminal GlcNAc residues are then galactosylated using GalTI and then sialylated with SA-PEG using ST3GalIII.

Conversely, TPA can be produce in yeast or fungal systems. Similar processing would be required for fungal derived material.

5

41. Generation and PEGylation of GlcNAc-Asn structures: TPA produced in Yeast

This example sets forth the preparation of PEGylated GlcNAc-Asn structures on a peptide such as TPA expressed in yeast.

10

Yeast expression is expected to result in a TPA which contains a single N-linked mannan-type structure. This recombinant glycoprotein is first treated with endoglycosidase H to generate GlcNAc structures on the asparagine (Asn) residues on the peptide.

The GlcNAc-Asn structures on the peptide/protein backbone are then be modified with galactose or galactose-PEG using UDP-galactose or UDP-galactose-6-PEG, respectively, and a galactosyltransferase such as GalT1. In one case, the galactose-PEG is the terminal residue. In the second case, the galactose is further modified with SA-PEG using a CMP-SA-PEG donor and a sialyltransferase such as ST3GalIII. In another embodiment, the GlcNAc-Asn structures on the peptide/protein backbone may be galactosylated and sialylated as described above, and then further sialylated using CMP-SA-PEG and an α 2,8-sialyltransferase such as the enzyme encoded by the *Campylobacter jejuni* cst-II gene.

15
20

Transferrin

42. GlycoPEGylation of Transferrin

This example sets forth the preparation of asialotransferrin and its sialylation with PEG-CMP-sialic acid.

25

Preparation of Asialo-transferrin. Human-derived holo-Transferrin, (10 mg) was dissolved in 500 μ L of 50 mM NaOAc, 5 mM CaCl₂, pH 5.5. To this solution was added 500 mU Neuraminidase II (*Vibrio cholerae*) and the reaction mixture was shaken gently for 20.5 hours at 37 °C. The reaction mixture was added to the prewashed N-(p-aminophenyl)oxamic acid-agarose conjugate (600 μ L) and the washed beads gently rotated

30

for 24 hours at 4 °C. The mixture was centrifuged at 10,000 rpm and the supernatant was collected. The reaction mixture was adjusted to 5 mM EDTA by addition of 100 µL of 30 mM EDTA to the washed beads, which were gently rotated for 20 hours at 4 °C. The suspension was centrifuged for 2 minutes at 10,000 rpm and the supernatant was collected.

- 5 The beads were washed 5 times with 0.35 mL of 50 mM NaOAc, 5 mM CaCl₂, 5 mM EDTA, pH 5.5 and all supernatants were pooled. The enzyme solution was dialyzed twice at 4 °C into 15 mM Tris-HCl, 1 M NaCl, pH 7.4. 0.3 mL of the transferrin solution (3.3 mL total) was removed and dialyzed twice against water. The remainder was dialyzed twice more at 4 °C against phosphate buffered saline. The dialyzed solution was stored at -20 °C.
- 10 Protein samples were analyzed by IEF Electrophoresis. Samples (9 µL, 25 µg) were diluted with 16 µL Tris buffer and mixed with 25 µL of the sample loading buffer and applied to Isoelectric Focusing Gels (pH 3-7). Gels were run and fixed using standard procedures. Gels were stained with Colloidal Blue Stain.

- Sialyl-PEGylation of asialo-Transferrin.** Desialylated transferrin (250 µg) and
- 15 CMP-sialic acid or CMP-SA-PEG (1 kDa or 10 kDa)(0.05 µmol) were dissolved in 69 µL 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2 in 1.5 mL plastic tubes. The tubes were vortexed briefly and 100 mU ST3Gal3 (90 µL) were added (total volume 250 µL). The tubes were vortexed again and mixed gently for 24 hours at 32 °C. The reactions were stopped by freezing at -80 °C. Novex Tris-Glycine 8-16% 1 mm gels were used for SDS
- 20 PAGE analysis (Figure 190). Samples (25 µL, 25 µg) were mixed with 25 µL of sample loading buffer and 0.4 µL of β-mercaptoethanol and heated for 6 minutes at 85 °C. Gels were run using standard conditions and stained with Colloidal Blue Stain. IEF gels were also performed as described above Figure 191). Samples were also dialyzed against water analyzed by MALDI-TOF.

- 25 **Results.** MALDI was also performed. Native transferrin (78729); asialotransferrin (78197); resialylated transferrin (79626/80703); with SA-PEG 1k (79037 (1); 80961 (2); 82535 (3); 84778 (4)); with SA-PEG 5k (90003 (2); 96117 (3); 96117 (4)); with SA-PEG 10k (100336 (2); 111421 (3); 122510 (4)).

43. Transferrin-GDNF

This example sets forth the procedures for the glycoconjugation of proteins, and in particular, transferrin is glycoconjugated to GDNF. Transferrin-SA-Linker-Gal-UDP is prepared from transferrin. The galactose residue is removed from GDNF glycans, and
5 Transferrin-SA-Linker-Gal-UDP is conjugated to GDNF glycans using a galactosyltransferase.

Preparation of agalacto-GDNF. GDNF produced in NSO cells (NSO murine myeloma cells) is dissolved at 2.5 mg/mL in 50 mM Tris 50 mM Tris-HCl pH 7.4, 0.15 M NaCl, and is incubated with 300 mU/mL beta-galactosidase-agarose conjugate for 16 hours at
10 32°C. To monitor the reaction a small aliquot of the reaction is diluted with the appropriate buffer and a IEF gel performed according to Invitrogen procedures. The mixture is centrifuged at 10,000 rpm and the supernatant is collected. The supernatant is dialyzed at 4 °C against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃ and then twice more against 50 mM Tris -HCl pH 7.4, 1 M NaCl, 0.05% NaN₃. The dialyzed solution is then concentrated
15 using a Centricon Plus 20 centrifugal filter and stored at -20 °C. The conditions for the IEF gel are run according to the procedures and reagents provided by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Transferrin-SA-Linker-Gal-UDP. Asialo-transferrin is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 0.05% NaN₃, pH 7.2. The solution is
20 incubated with CMP-sialic acid-linker-Gal-UDP (molar amount to add 1 molar equivalent of nucleotide sugar to transferrin) and 0.1 U/mL of ST3Gal3 at 32°C for 2 days. To monitor the incorporation of sialic acid, a small aliquot of the reaction has ¹⁴C-SA-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label
25 incorporation into the peptide is quantitated using an in-line radiation detector.

The solution is incubated with 5 mM CMP-sialic acid and 0.1 U/mL of ST3Gal3 (to cap any unreacted transferrin glycans) at 32°C for 2 days. The incorporation into the peptide is quantitated using an in-line UV detector. After 2 days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) and collecting
30 fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE

and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

Preparation of Transferrin-SA-Linker-Gal-GDNF. The transferrin-SA-Linker-Gal-UDP prepared as described above is dissolved at 2.5 mg/mL in 50 mM Tris-HCl, 0.15 M NaCl, 5 mM MnCl₂, 0.05% NaN₃, pH 7.2. The solution is incubated with 2.5 mg/mL agalacto-GDNF and 0.1 U/mL of galactosyltransferase at 32°C for 2 days. To monitor the incorporation of galactose, a small aliquot of the reaction has ¹⁴C-galactose-UDP ligand added; the label incorporated into the peptide is separated from the free label by gel filtration on a Toso Haas G3000SW analytical column using PBS buffer (pH 7.1). The radioactive label incorporation into the peptide is quantitated using an in-line radiation detector.

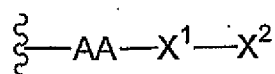
When the reaction is complete, the solution is incubated with 5 mM UDP-Gal and 0.1 U/mL of galactosyltransferase (to cap any unreacted transferrin glycans) at 32°C for 2 days followed by addition of 5 mM CMP-SA and 0.1 U/mL of ST3Gal3. After 2 additional days, the reaction mixture is purified using a Toso Haas G3000SW preparative column using PBS buffer (pH 7.1) collecting fractions based on UV absorption. The product of the reaction is analyzed using SDS-PAGE and IEF analysis according to the procedures and reagents supplied by Invitrogen. Samples are dialyzed against water and analyzed by MALDI-TOF MS.

The disclosures of each and every patent, patent application, and publication cited herein are hereby incorporated herein by reference in their entirety.

While this invention has been disclosed with reference to specific embodiments, it is apparent that other embodiments and variations of this invention may be devised by others skilled in the art without departing from the true spirit and scope of the invention. The appended claims are intended to be construed to include all such embodiments and equivalent variations.

What is claimed:

1. A cell-free, in vitro method of remodeling a peptide comprising poly(ethylene glycol), the peptide having the formula:



5 wherein

AA is a terminal or internal amino acid residue of the peptide;

$\text{X}^1 - \text{X}^2$ is a saccharide covalently linked to the AA, wherein

X^1 is a first glycosyl residue; and

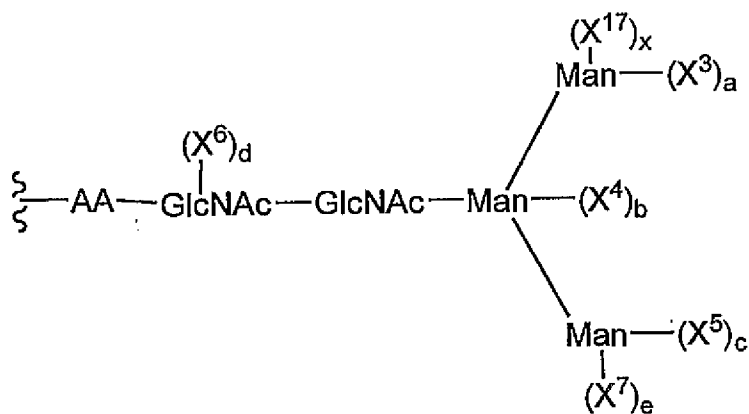
10 X^2 is a second glycosyl residue covalently linked to X^1 , wherein X^1 and X^2 are selected from monosaccharyl and oligosaccharyl residues;

the method comprising:

(a) removing X^2 or a saccharyl subunit thereof from the peptide, thereby forming a truncated glycan.

15 2. The method according to claim 1 wherein said truncated glycan is formed by removing a Sia residue.

3. The method according to claim 1 wherein said peptide has the formula:




20 wherein

X^3 , X^4 , X^5 , X^6 , X^7 , and X^{17} , are independently selected monosaccharyl or oligosaccharyl residues; and

a, b, c, d, e, and x are independently selected from the integers 0, 1 and 2.

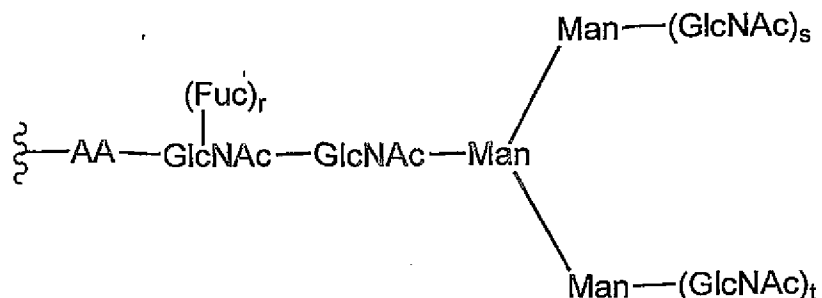- 5 4. The method according to claim 3 wherein said oligosaccharyl residue is a member selected from GlcNAc-Gal-Sia and GlcNAc-Gal.
5. The method according to claim 3 wherein at least one member selected from a, b, c, d, e and x is 1 or 2.
- 10 6. The method of claim 3, wherein said removing of step (a) produces a truncated glycan in which at least one of a, b, c, e and x are 0.
7. The method of claim 6, wherein X^3 , X^5 and X^7 are members
- 15 independently selected from (mannose)_z and (mannose)_z-(X^8)
 wherein
 X^8 is a glycosyl moiety selected from mono- and oligo-saccharides; and
 z is an integer between 1 and 20, wherein
 when z is 3 or greater, each (mannose)_z is independently selected from linear
- 20 and branched structures.
- 8. The method of claim 6 wherein X^4 is selected from the group consisting of GlcNAc and xylose.
- 25 9. The method of claim 6, wherein X^3 , X^5 and X^7 are (mannose)_u
 wherein
 u is selected from the integers between 1 and 20, and when u is 3 or greater,
 each (mannose)_u is independently selected from linear and branched structures.
- 30 10. The method according to claim 3 wherein said peptide has the formula:




wherein

r, s, and t are integers independently selected from 0 and 1.

5

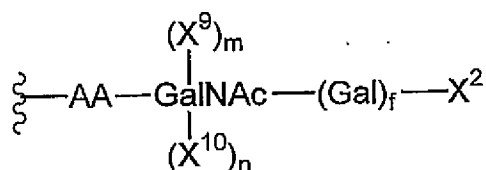11. The method of claim 1, wherein said peptide has the formula:




wherein

X⁹ and X¹⁰ are independently selected monosaccharyl or oligosaccharyl

10 residues; and

m, n and f are integers independently selected from 0 and 1.

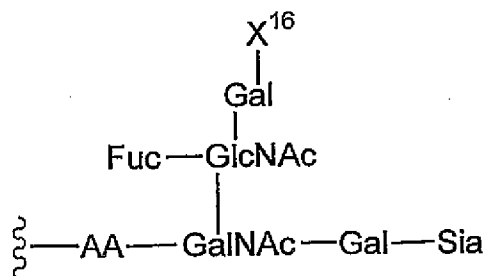
12. The method of claim 11, wherein said peptide has the formula:



15

wherein

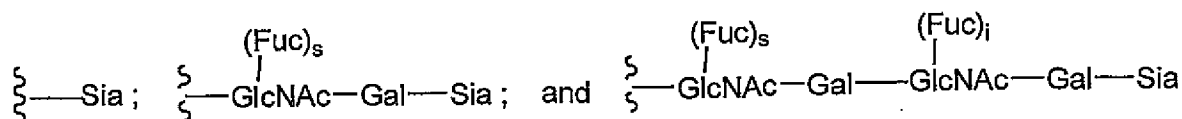
X¹⁶ is a member selected from:



wherein

s and i are integers independently selected from 0 and 1.

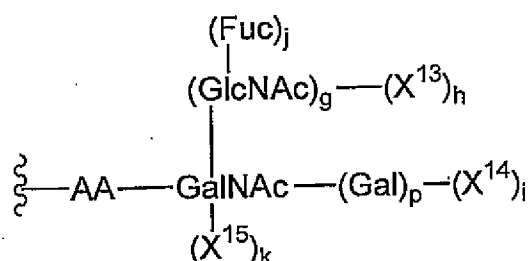- 5 13. The method of claim 12, wherein said peptide has the formula:



wherein

X^{13} , X^{14} , and X^{15} are independently selected glycosyl residues; and
g, h, i, j, k, and p are independently selected from the integers 0 and 1

10

14. The method according to claim 13 wherein at least one of g, h, i, j, k
and p is 1.

15

15. The method of claim 13, wherein
 X^{14} and X^{15} are members independently selected from GlcNAc and Sia; and
i and k are independently selected from the integers 0 and 1.

16. The method according to claim 15 wherein at least one of i and k is 1,
and if k is 1, g, h, and j are 0.

20

17. The method according to claim 1, further comprising:
(b) contacting the truncated glycan with at least one glycosyltransferase
and at least one glycosyl donor under conditions suitable to transfer the at least one glycosyl

donor to the truncated glycan, thereby remodeling said peptide comprising poly(ethylene glycol).

18. The method according to claim 17 wherein said glycosyl donor
5 comprises a modifying group covalently linked thereto.

19. The method of claim 1, further comprising:

(c) removing X^1 , thereby exposing AA.

20. The method according to claim 19, further comprising:

(d) contacting AA with at least one glycosyltransferase and at least one
glycosyl donor under conditions suitable to transfer said at least one glycosyl donor to AA,
thereby remodeling said peptide comprising poly(ethylene glycol).

21. The method according to claim 20 wherein said at least one glycosyl
donor comprises a modifying group covalently linked thereto.

22. The method according to claim 21 wherein said modifying group is
poly(ethylene glycol).

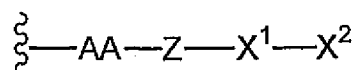
23. The method according to claim 22 wherein said poly(ethylene glycol)
has a molecular weight distribution that is essentially homodisperse.

24. The method of claim 17, further comprising:

(e) prior to step (b), removing a group added to said saccharide during
post-translational modification.

25. The method of claim 24 wherein said group is a member selected from
phosphate, sulfate, carboxylate and esters thereof.


26. The method of claim 1 wherein said peptide has the formula:



wherein

Z is a member selected from O, S, NH and a cross-linker.

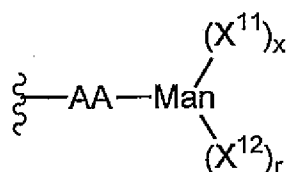- 5 27. The method of claim 1, wherein said peptide has the formula:



wherein

X^{11} and X^{12} are independently selected glycosyl moieties; and

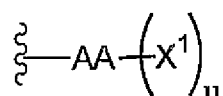
r and x are integers independently selected from 0 and 1.

10

28. The method of claim 27, wherein X^{11} and X^{12} are $(\text{mannose})_q$, wherein q is selected from the integers between 1 and 20, and when q is three or greater, $(\text{mannose})_q$ is selected from linear and branched structures.

- 15 29. A pharmaceutical composition comprising a pharmaceutically acceptable diluent and a remodeled peptide according to claim 1.

30. A cell-free, in vitro method of remodeling a peptide comprising poly(ethylene glycol), said peptide having the formula:



20

wherein

AA is a terminal or internal amino acid residue of said peptide;
X¹ is a glycosyl residue covalently linked to said AA, selected from
monosaccharyl and oligosaccharyl residues; and
u is an integer selected from 0 and 1,

5 said method comprising:

contacting said peptide with at least one glycosyltransferase and at least one
glycosyl donor under conditions suitable to transfer said at least one glycosyl donor to said
truncated glycan, thereby remodeling said peptide.

10 31. The method according to claim 30 wherein said at least one glycosyl
donor comprises a modifying group covalently linked thereto.

32. The method according to claim 30 wherein said modifying group is
poly(ethylene glycol).

15 33. The method according to claim 32 wherein said poly(ethylene glycol)
has a molecular weight distribution that is essentially homodisperse.

20 34. A pharmaceutical composition comprising a pharmaceutically
acceptable diluent and a remodeled peptide according to claim 30.

1/497


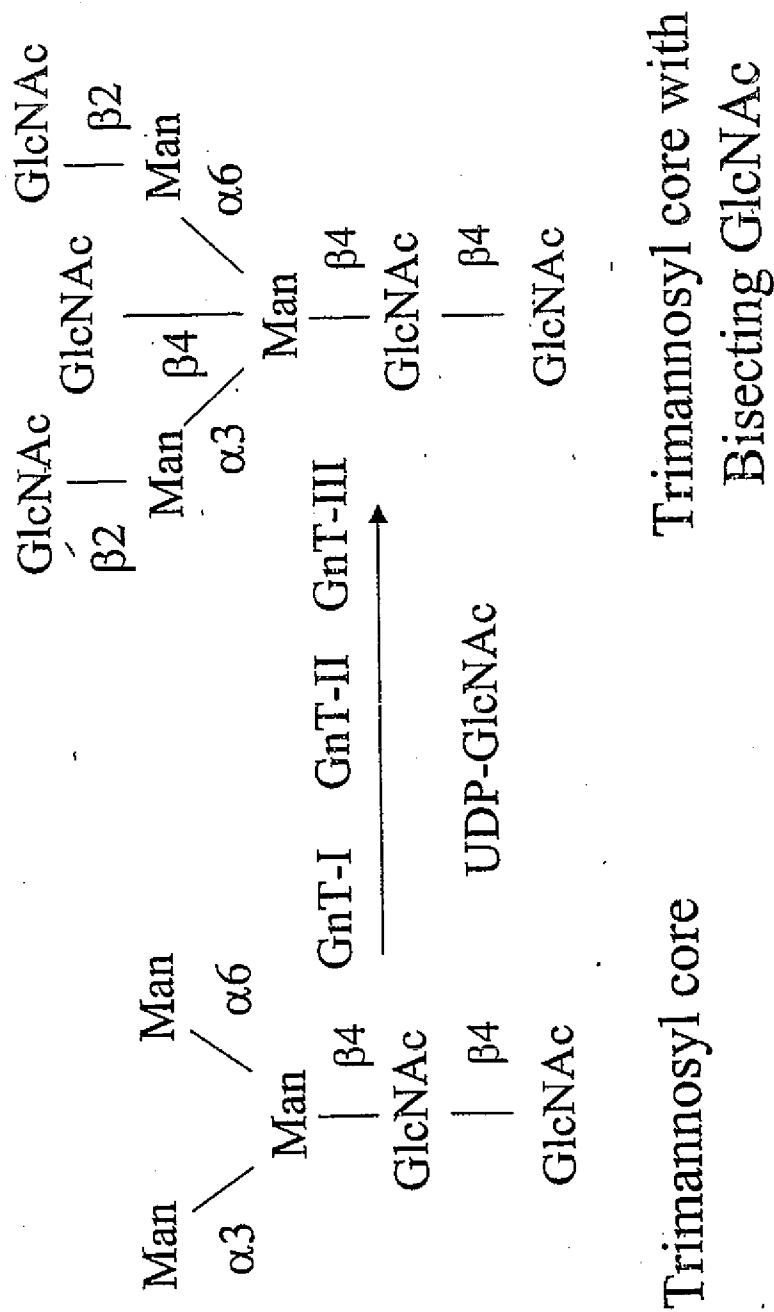


FIG. 1

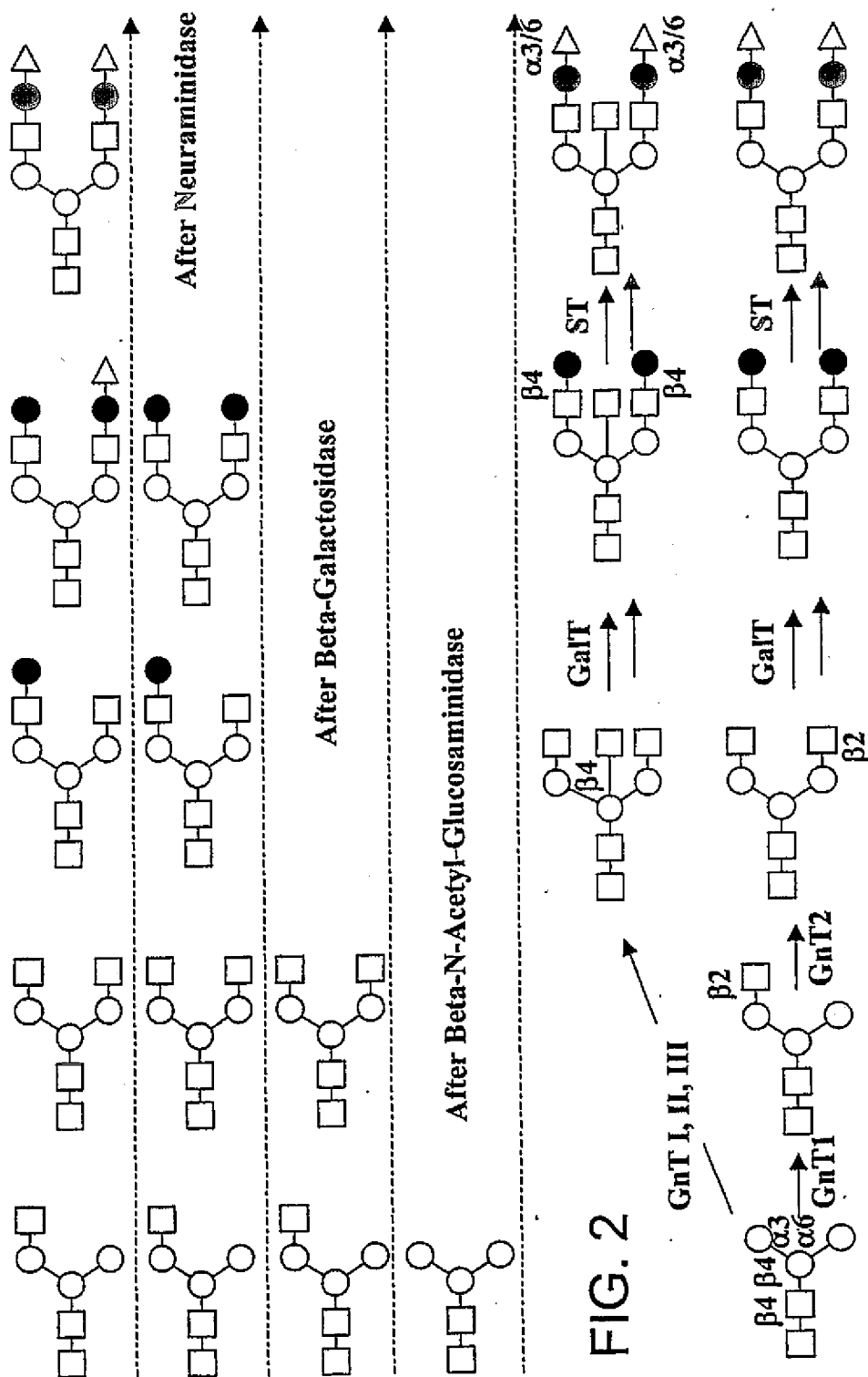2/497



3/497

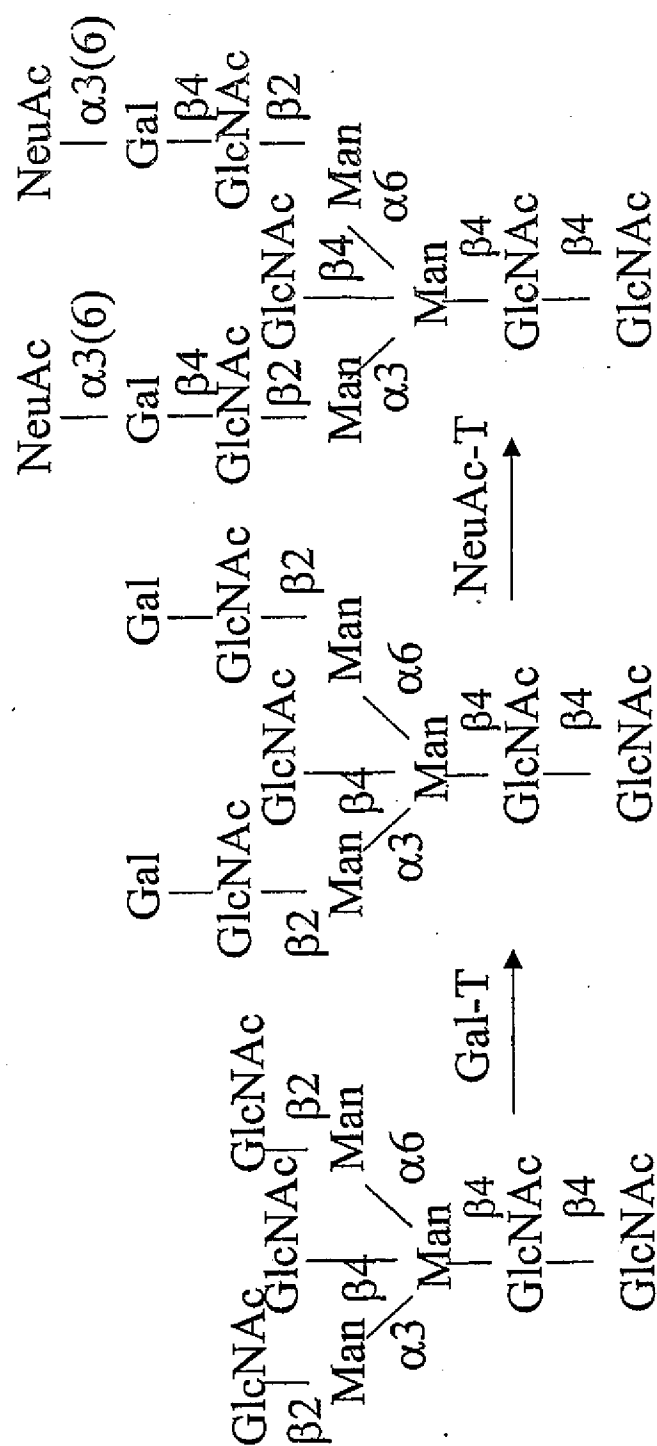


FIG. 3

4/497

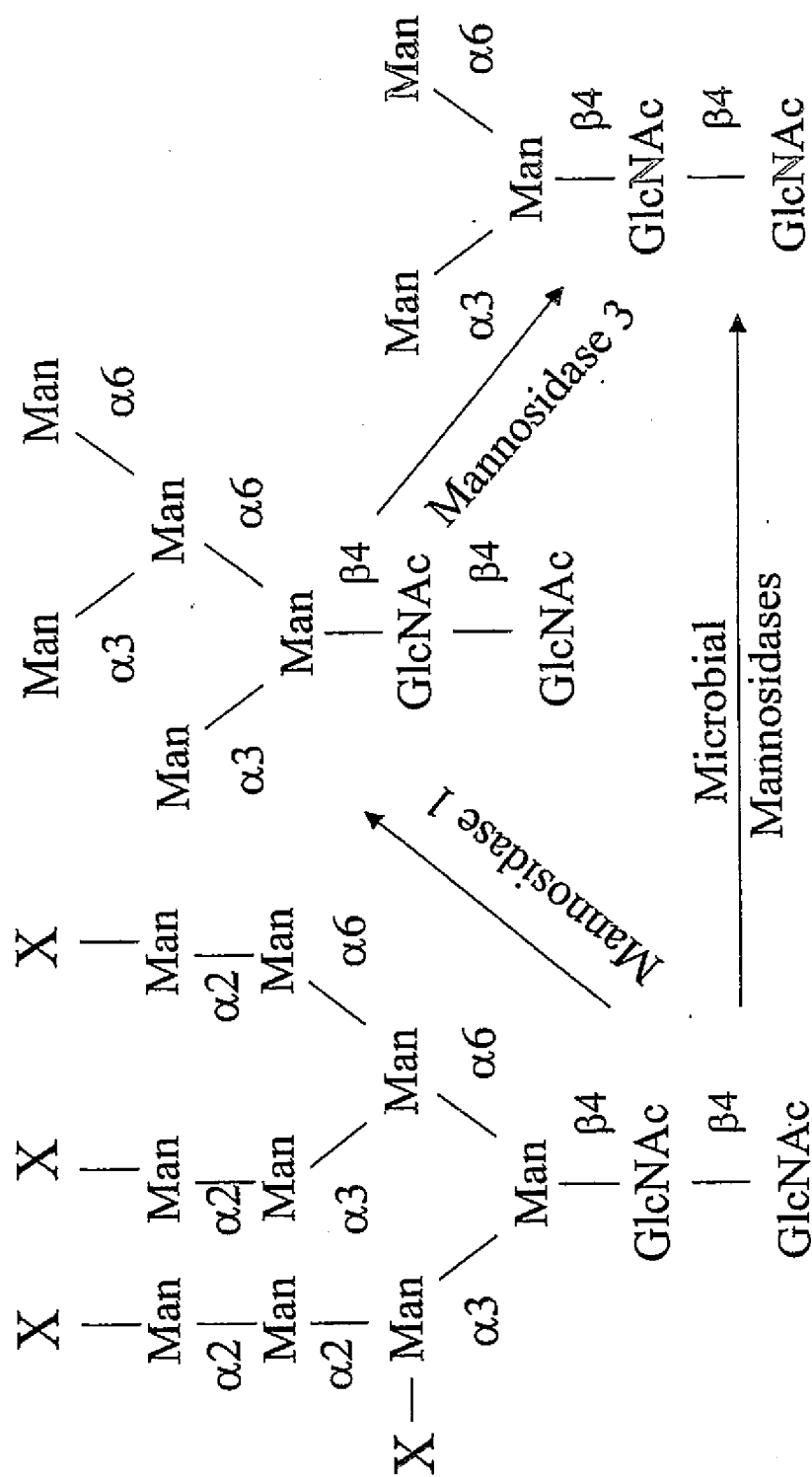


FIG. 4

5/497

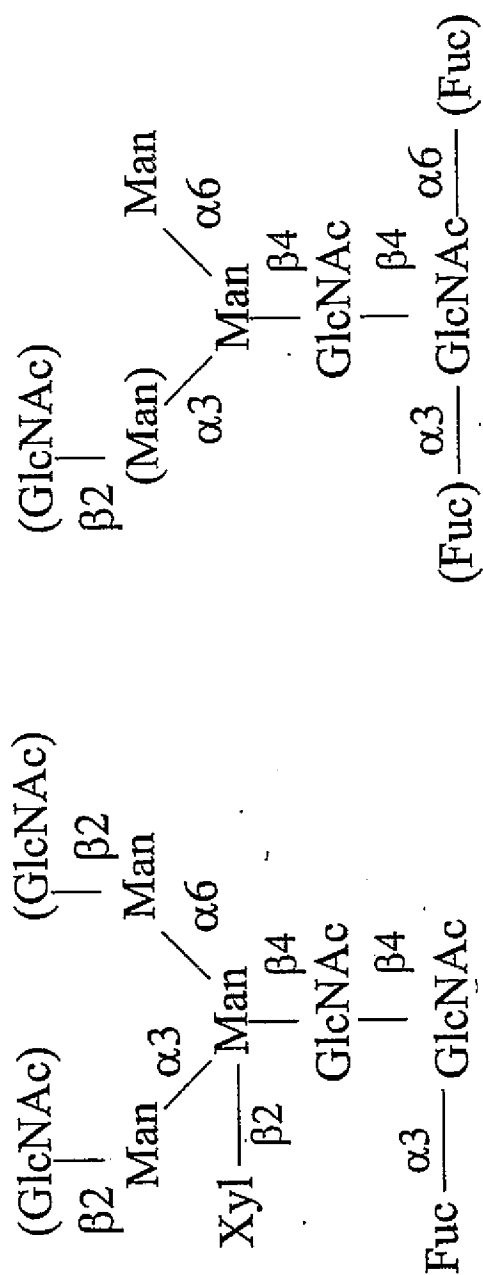


FIG. 6

FIG. 5

6/497

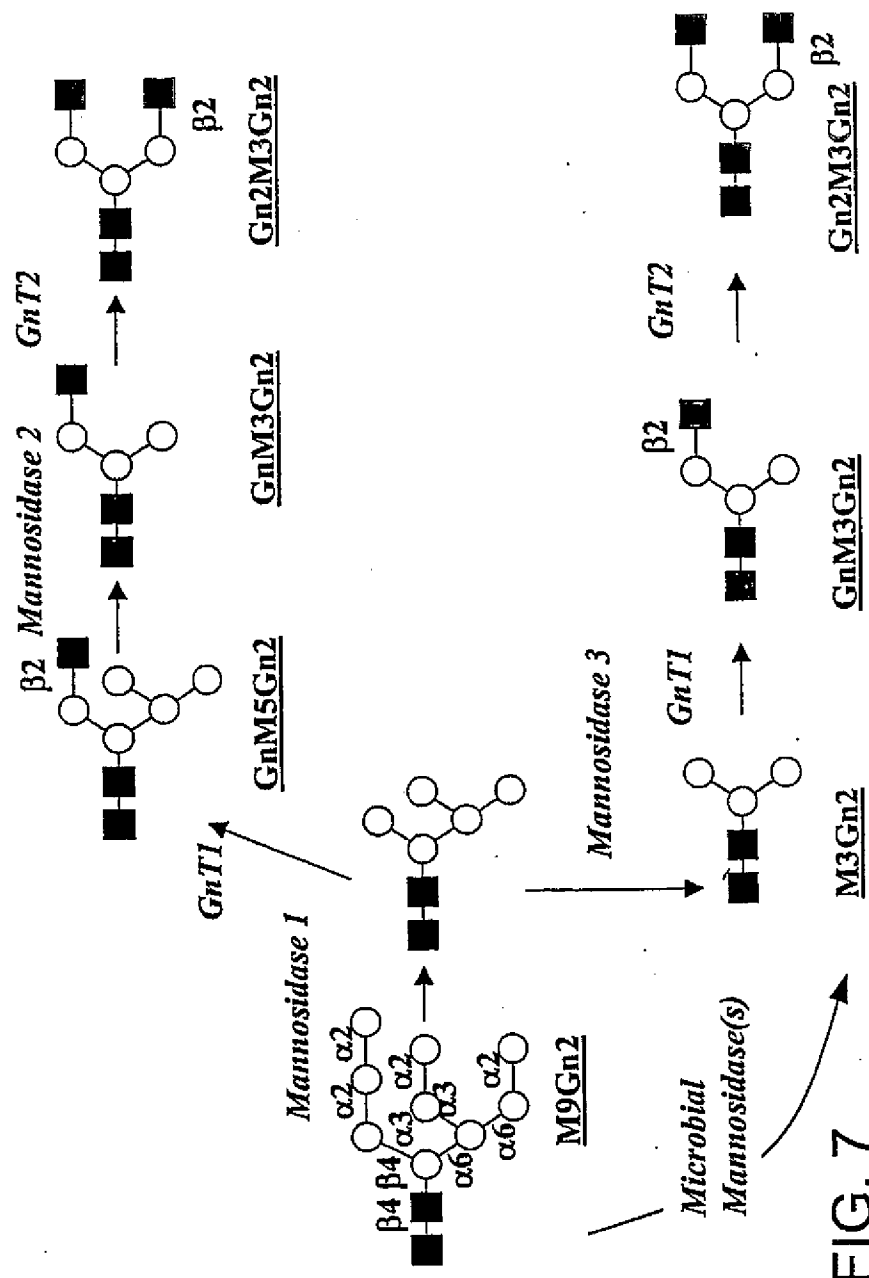


FIG. 7

7/497

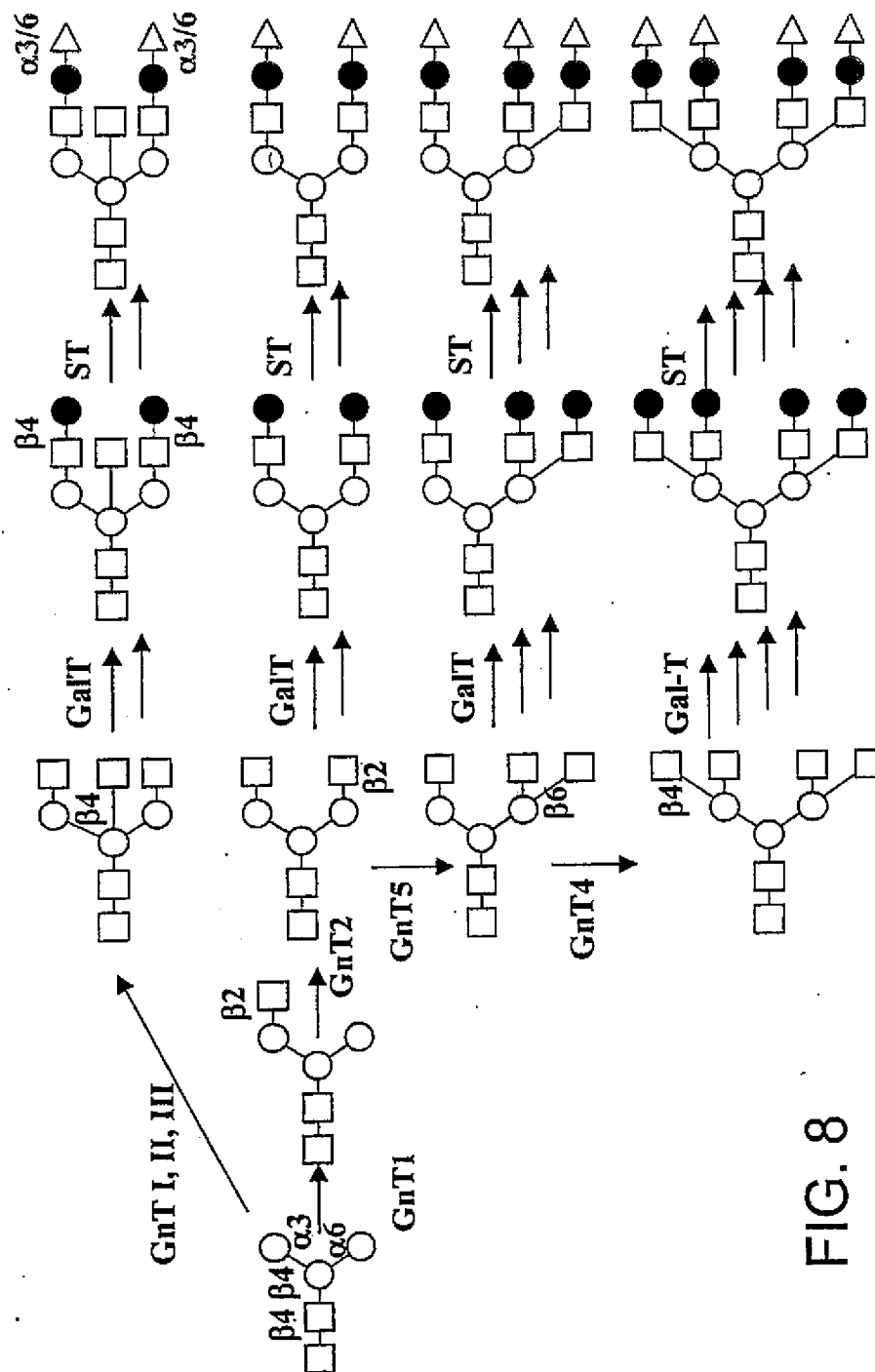


FIG. 8

8/497

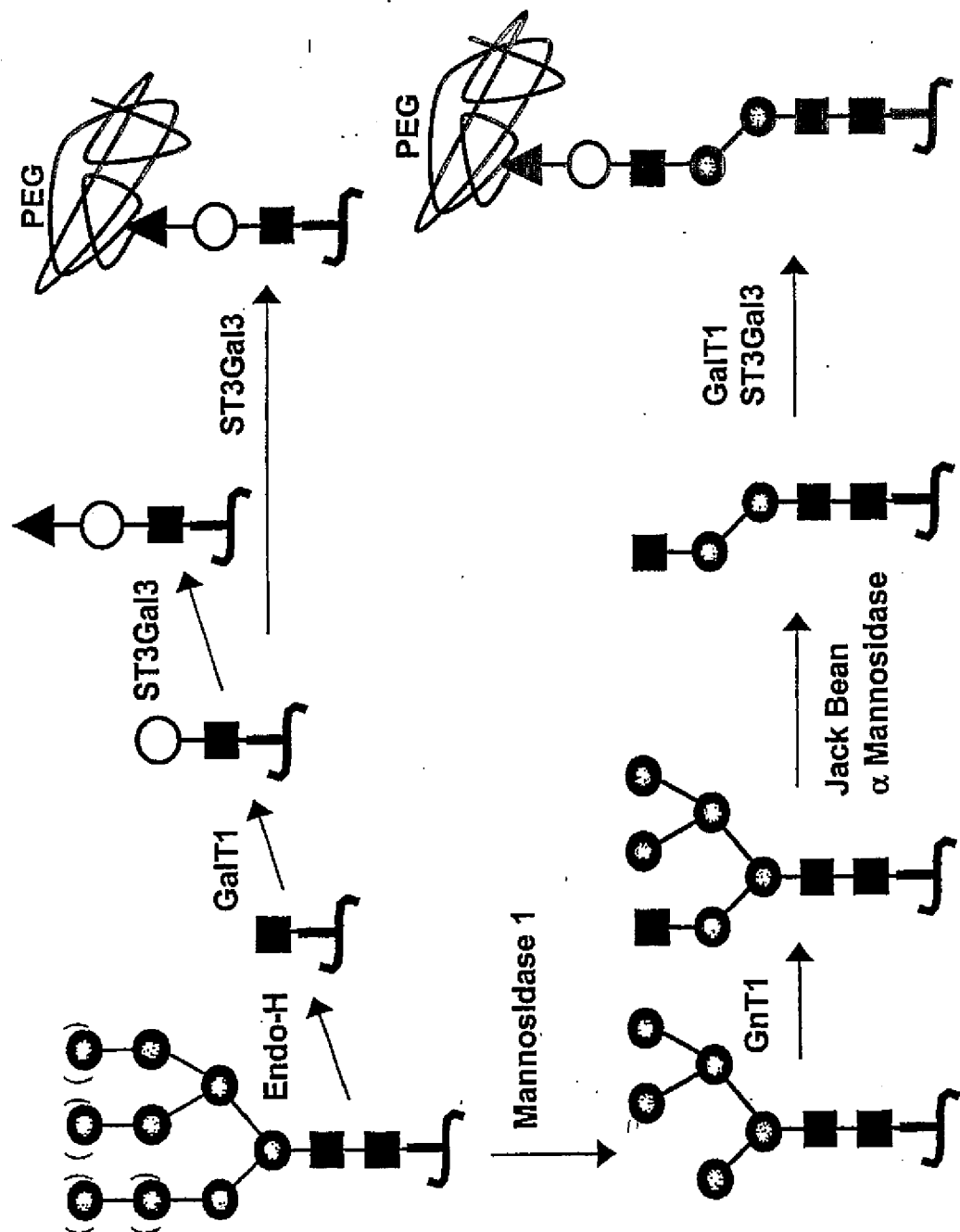


FIG. 9

9/497

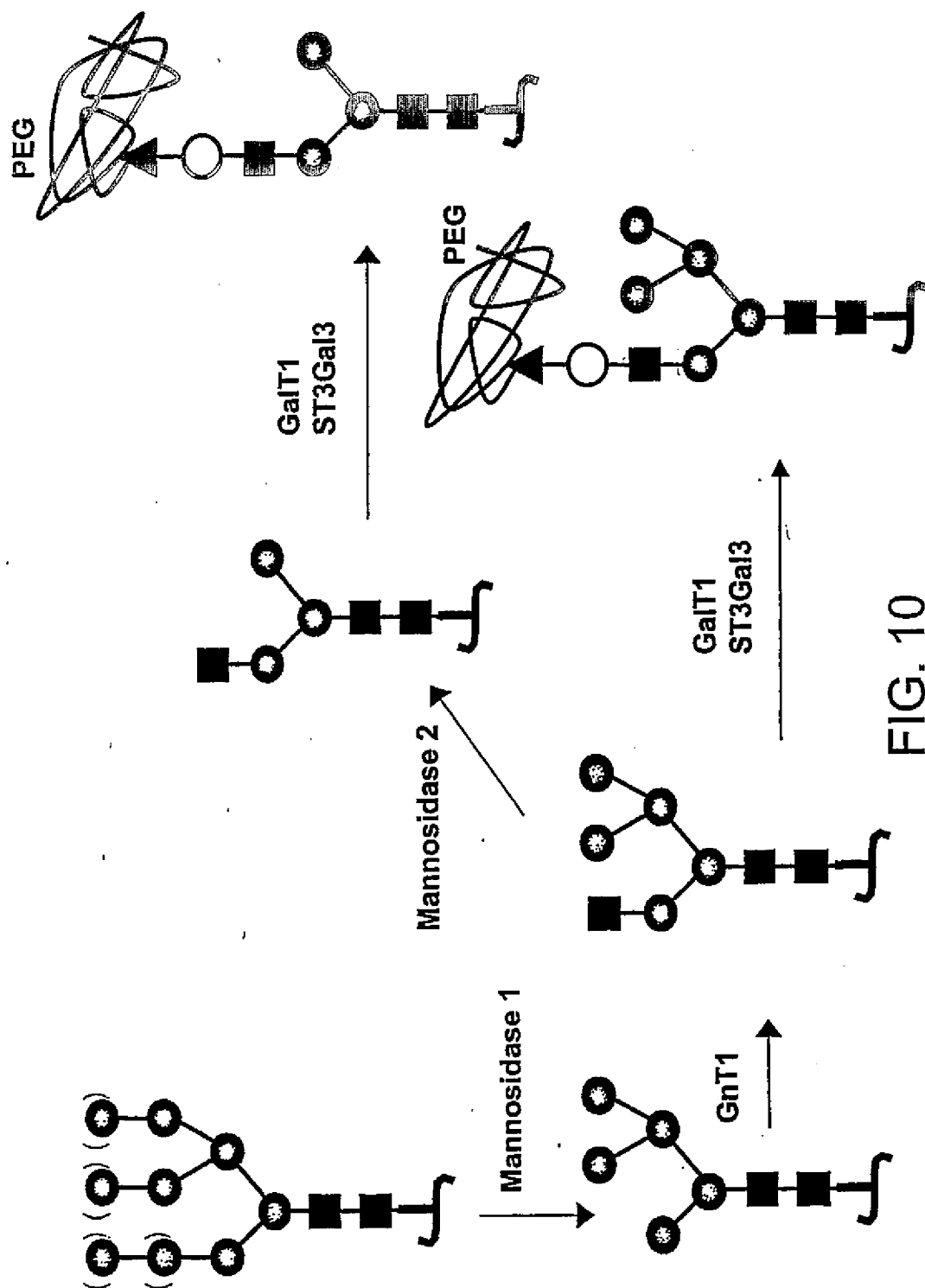


FIG. 10

10/497

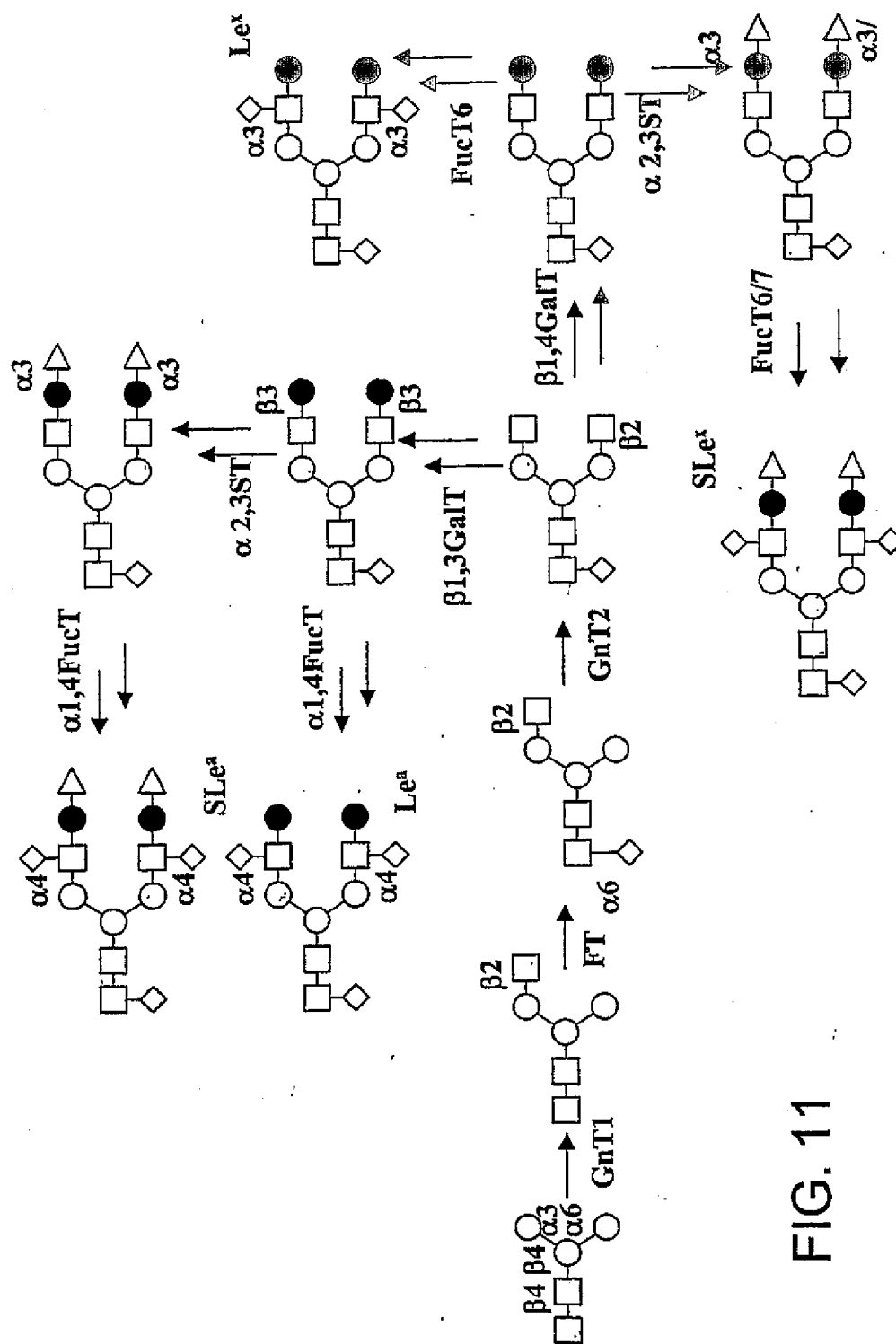


FIG. 11

11/497

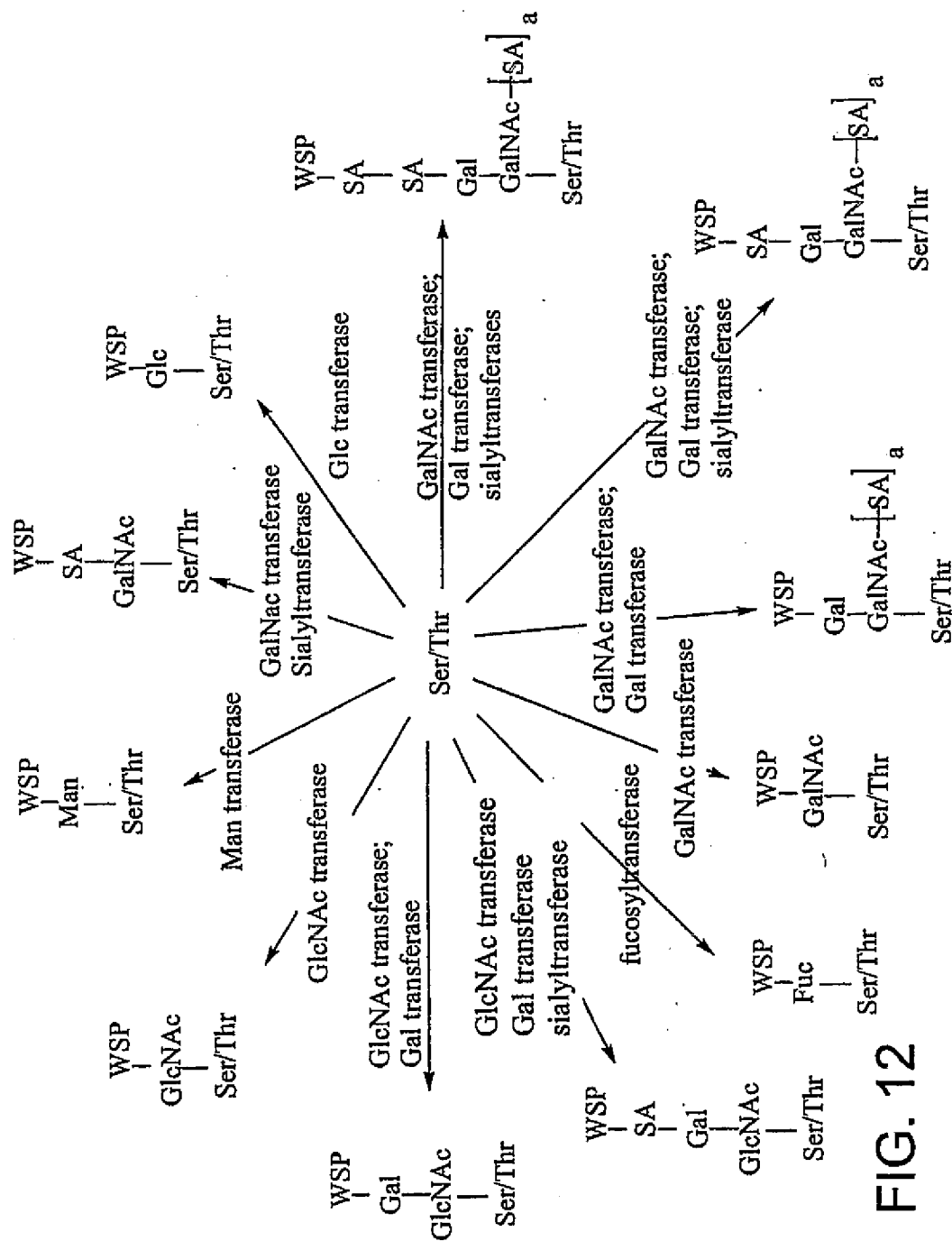


FIG. 12

12/497


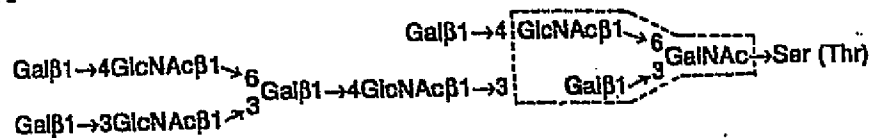
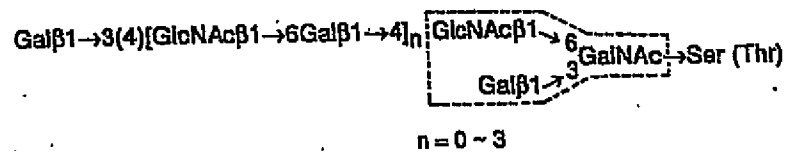
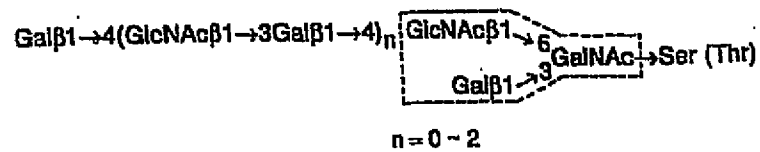
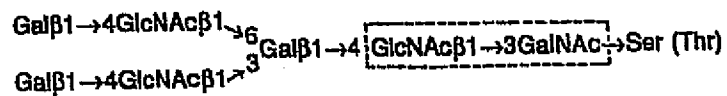
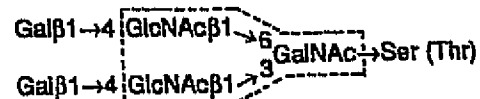
Core 1*Core 2**Core 3**Core 4*

FIG. 13

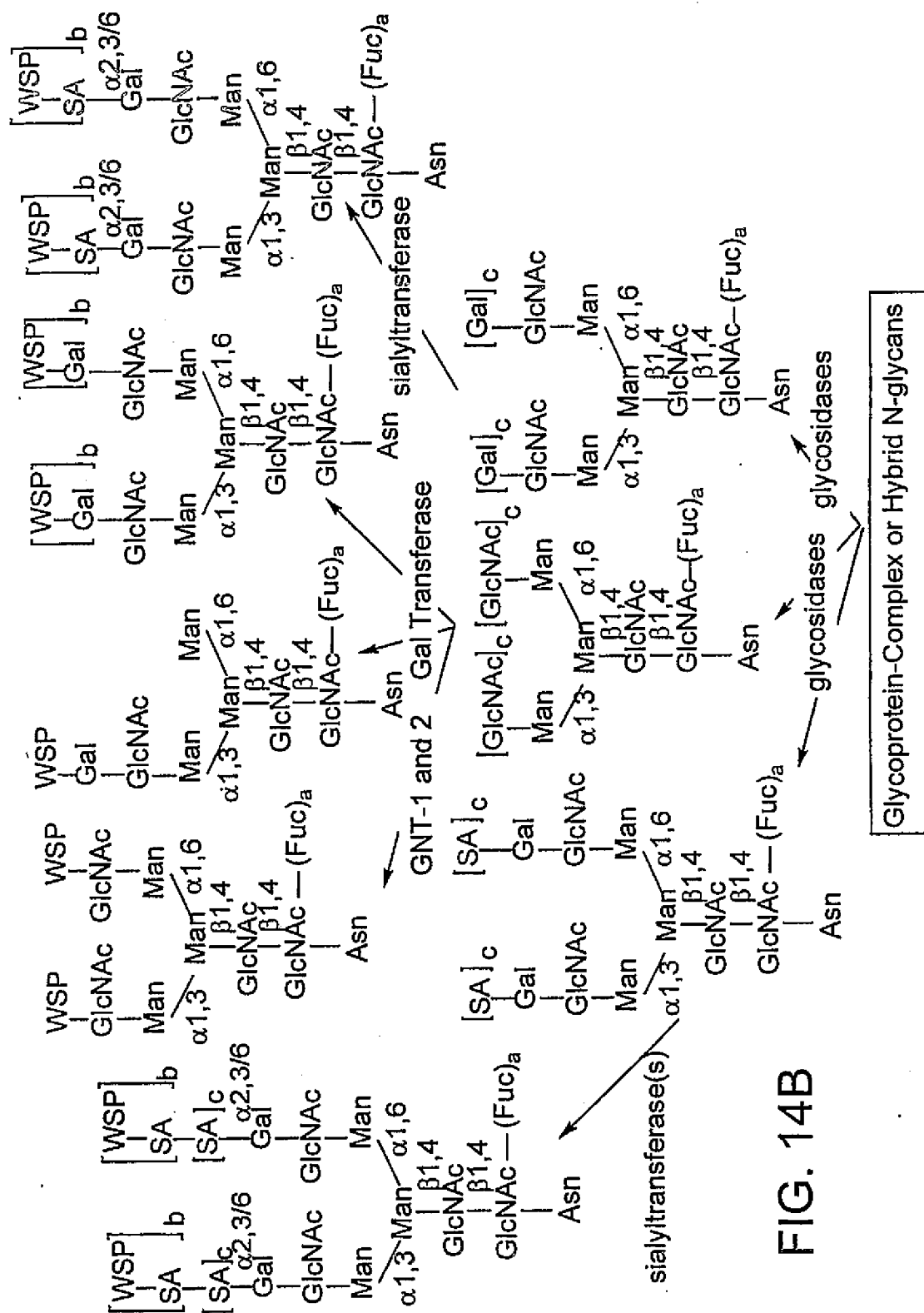14/497



15/497

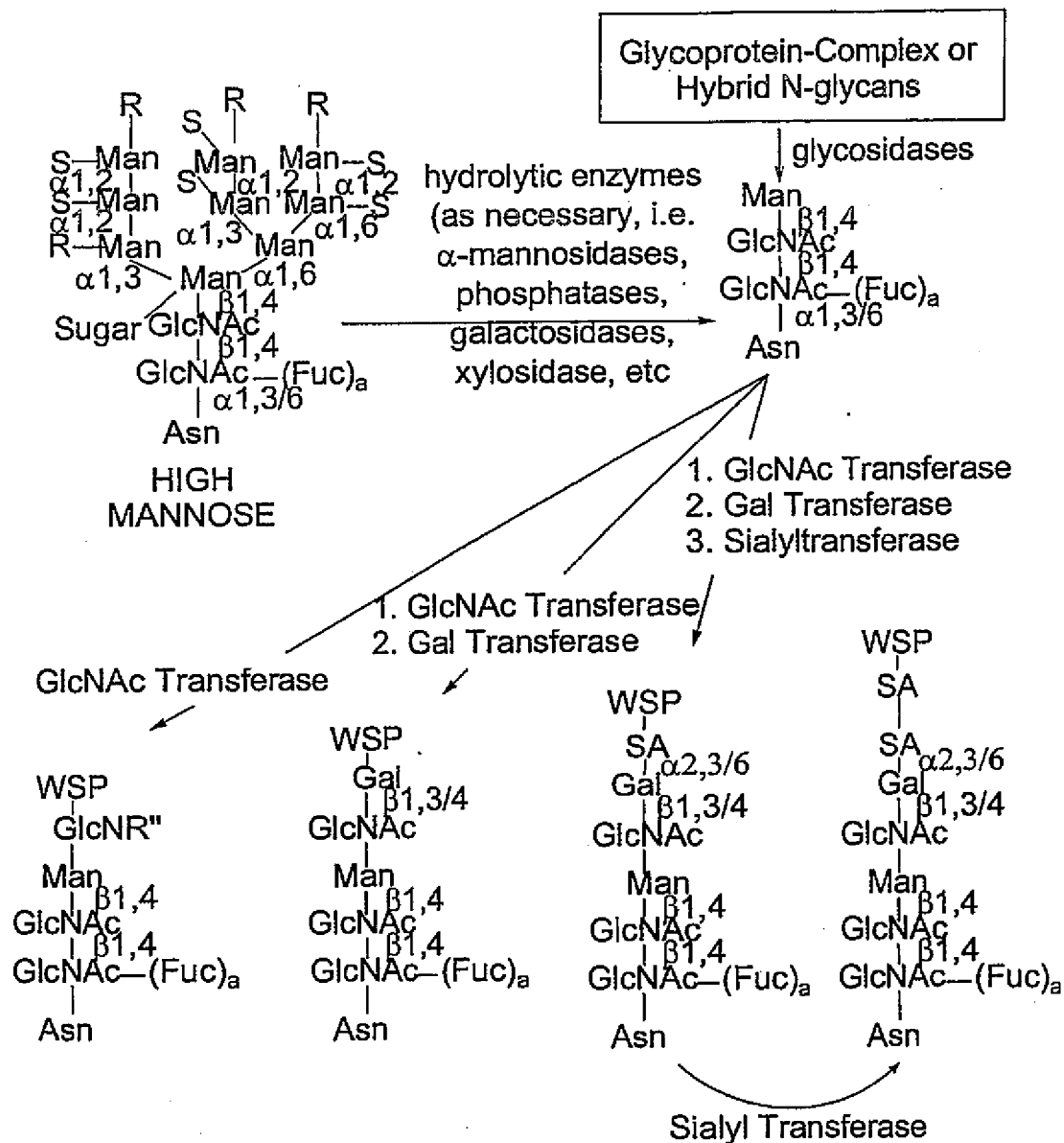


FIG. 15

16/497

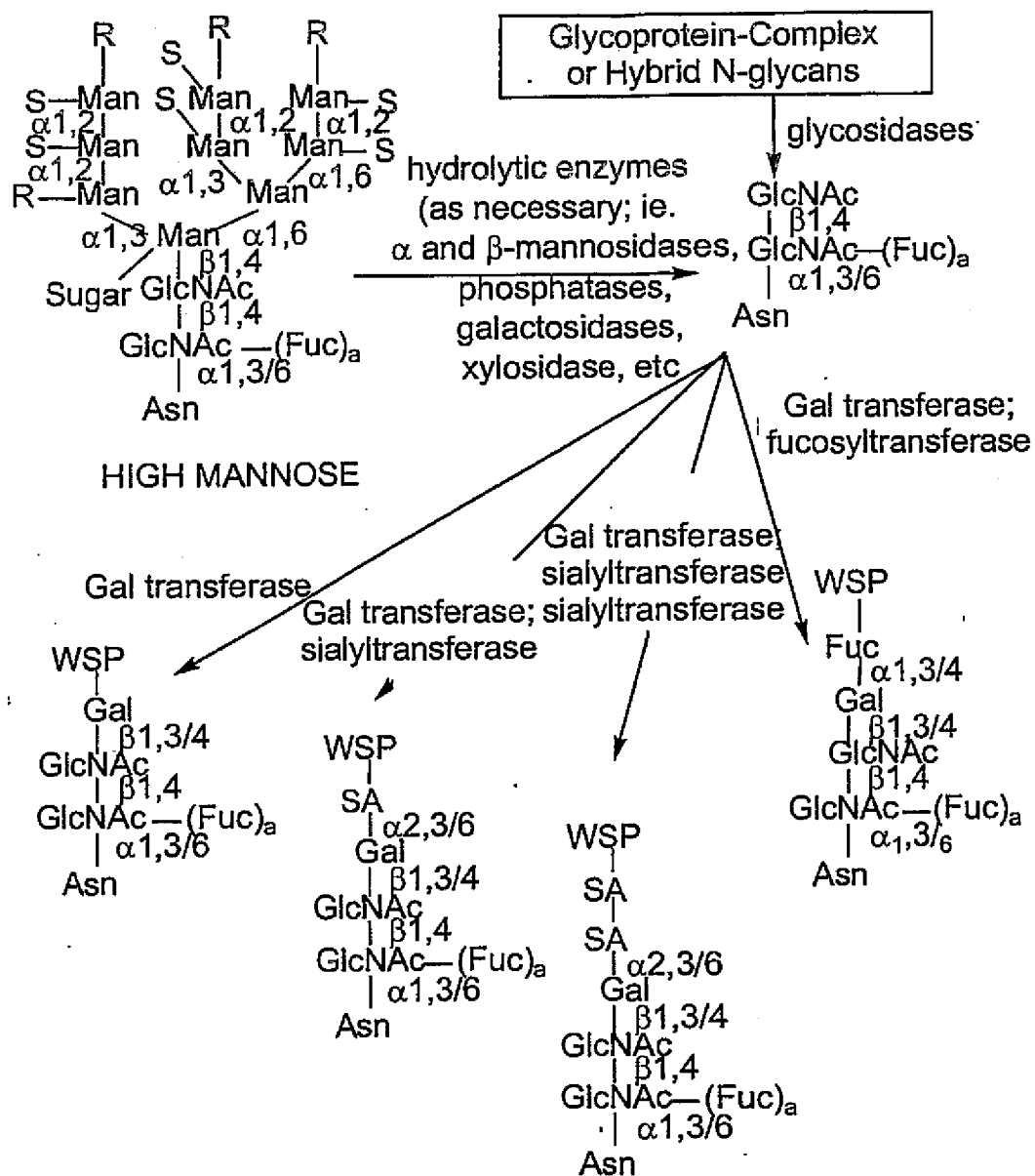


FIG. 16

17/497

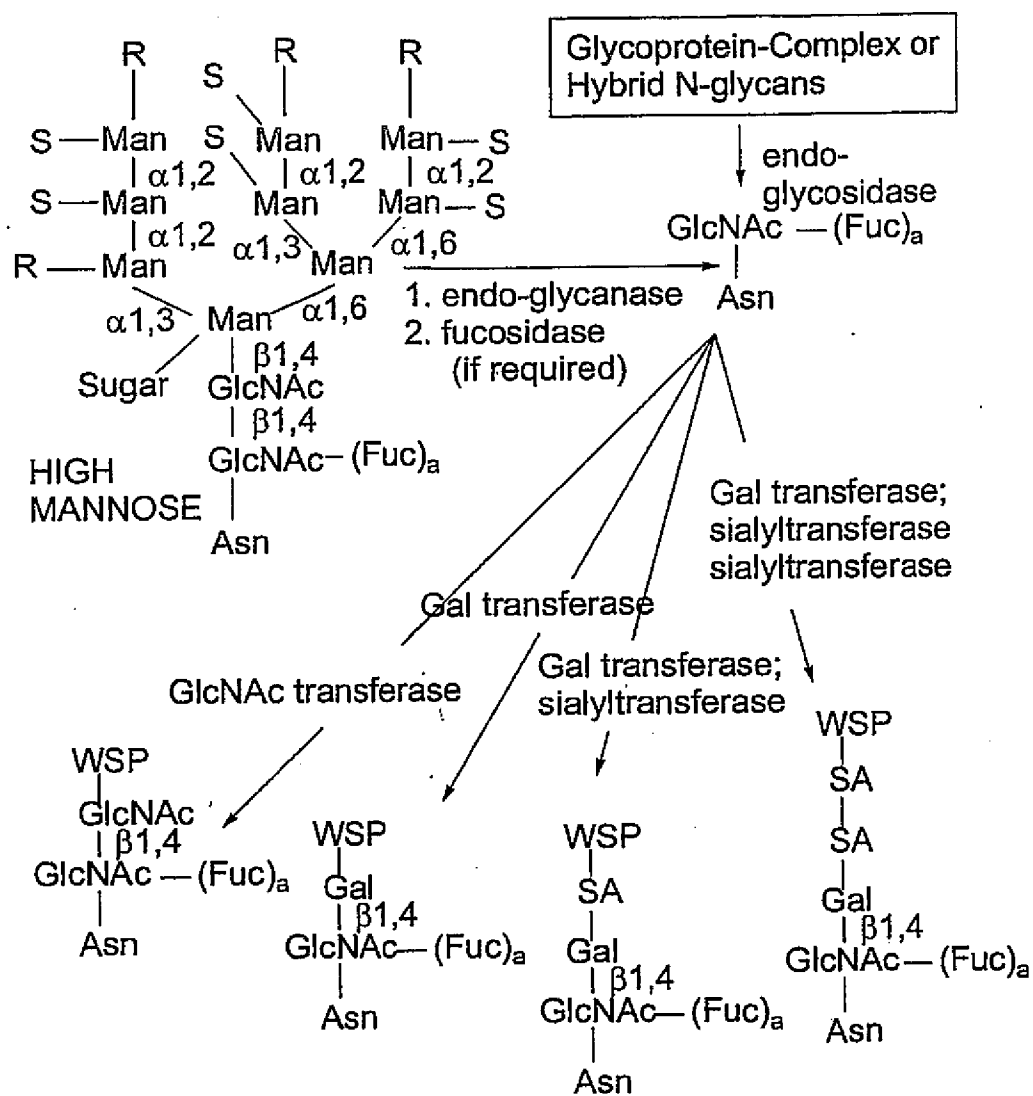


FIG. 17

18/497


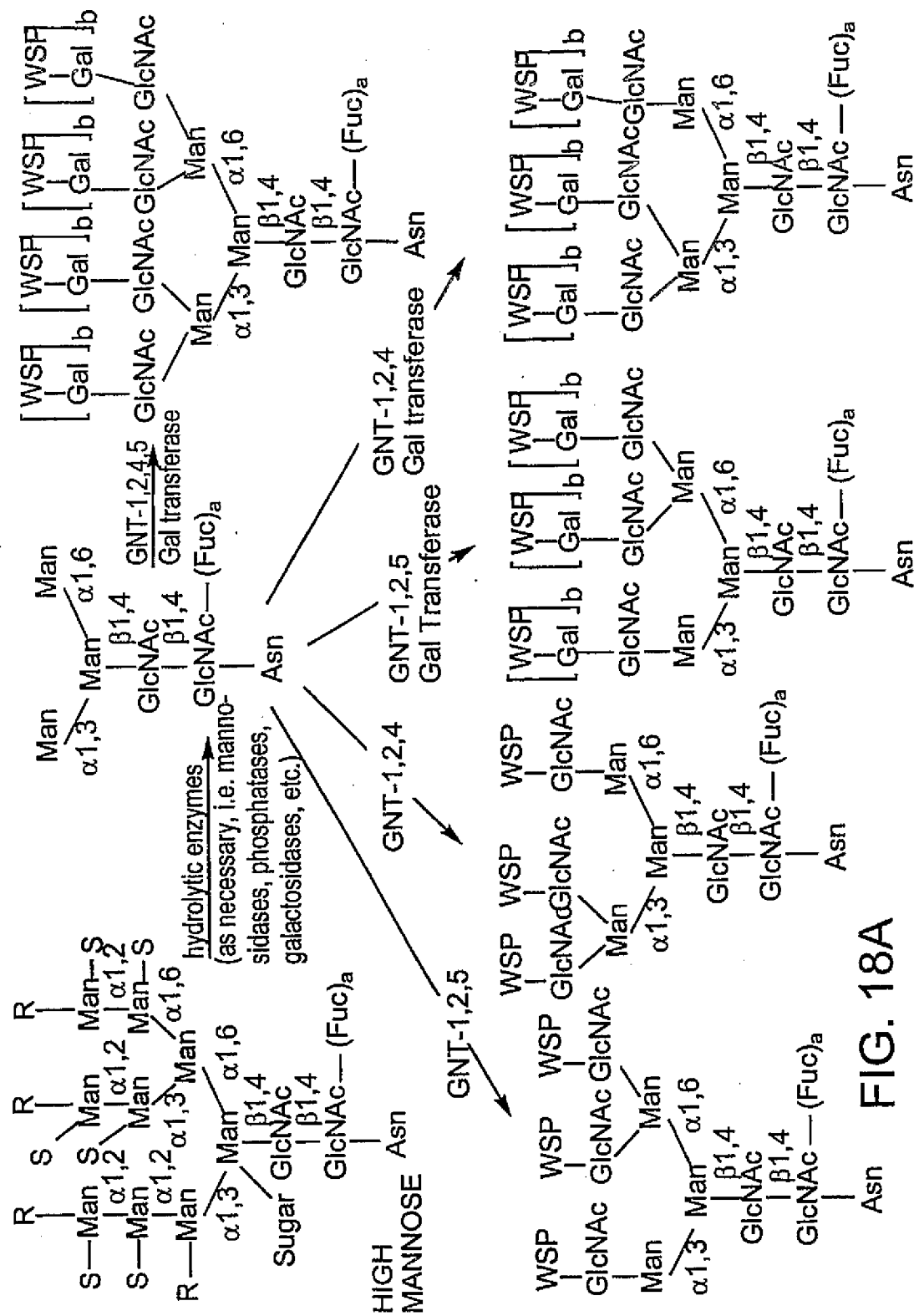



FIG. 18A

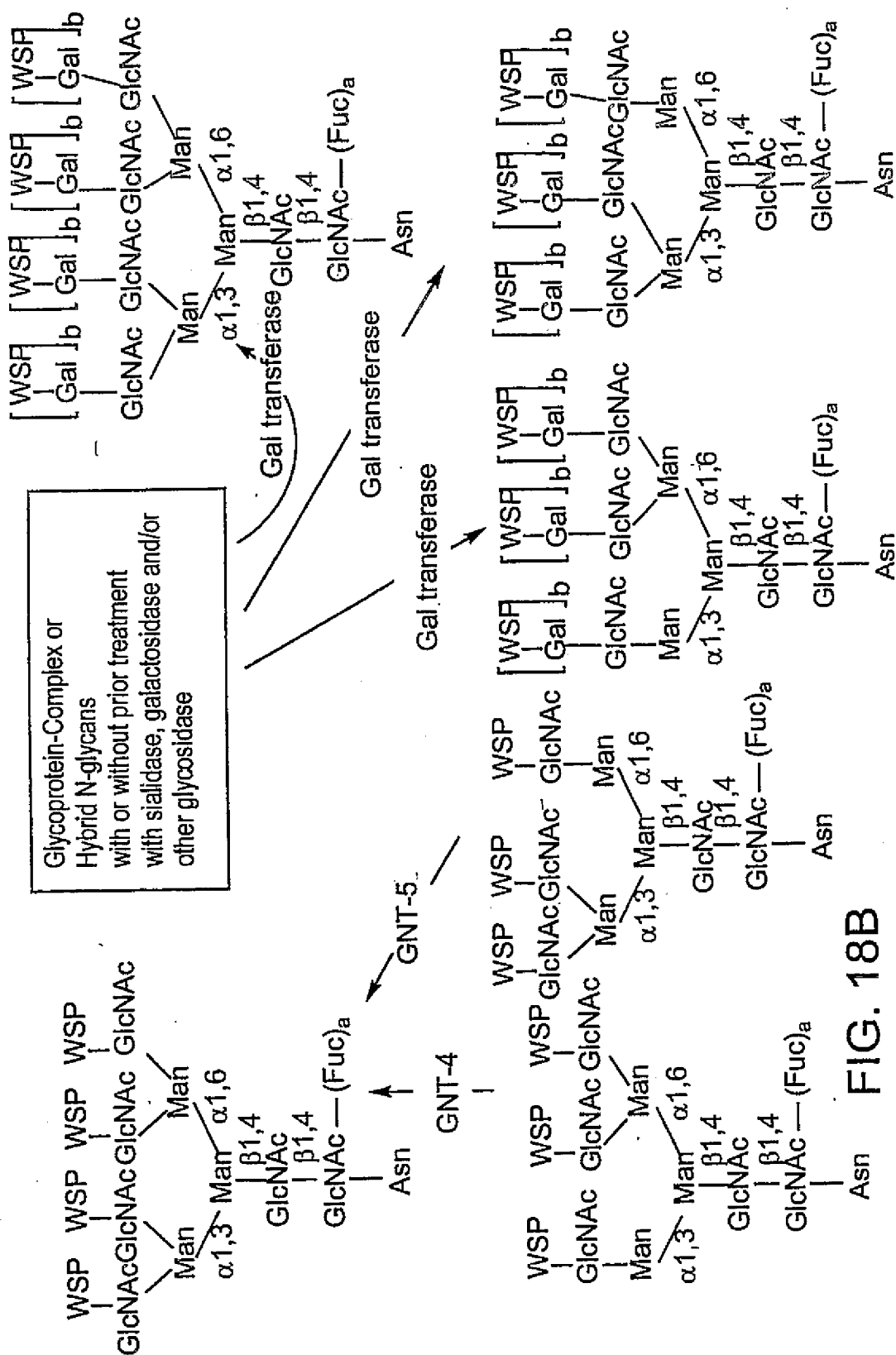19/497



1883

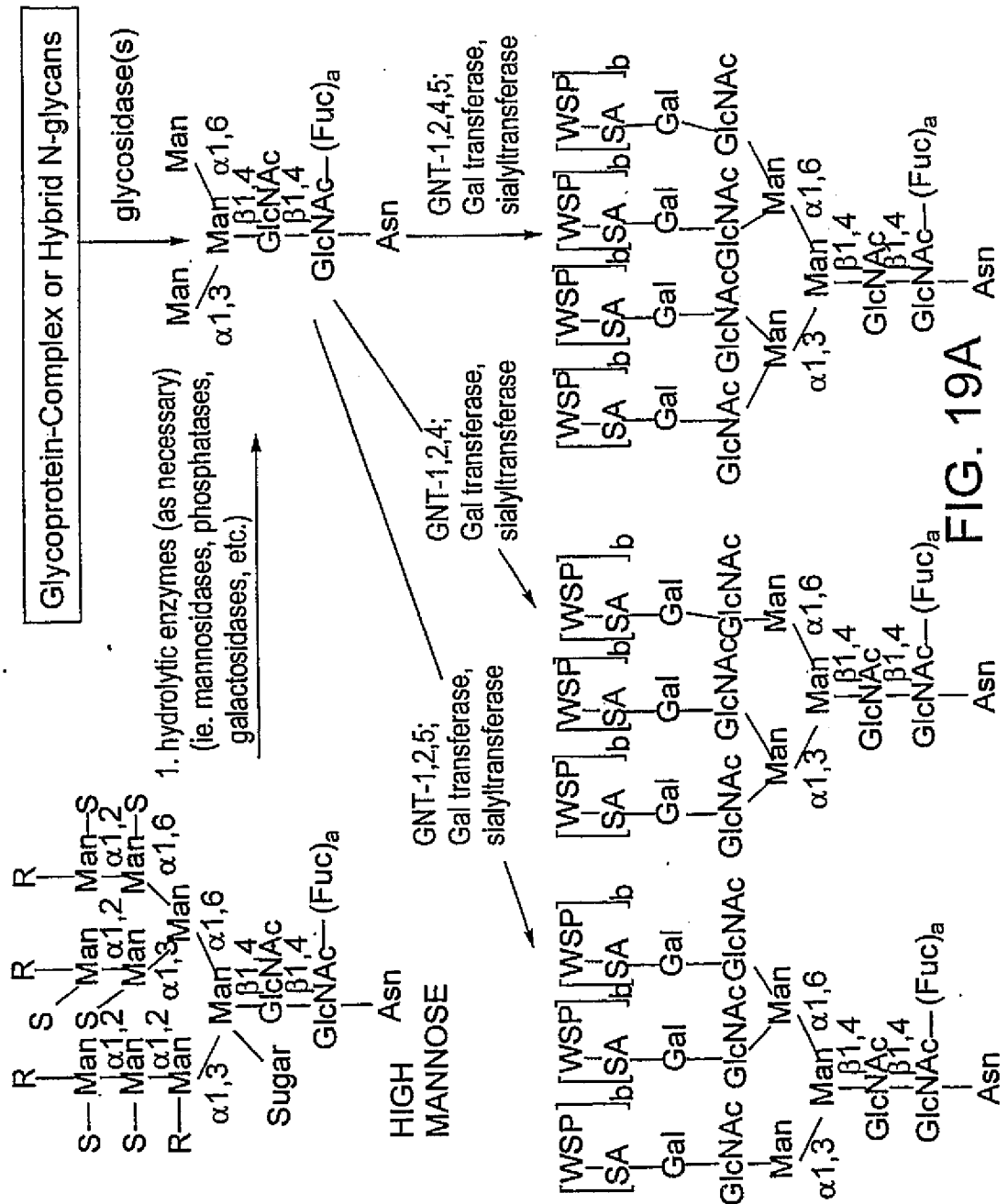20/497



22/497

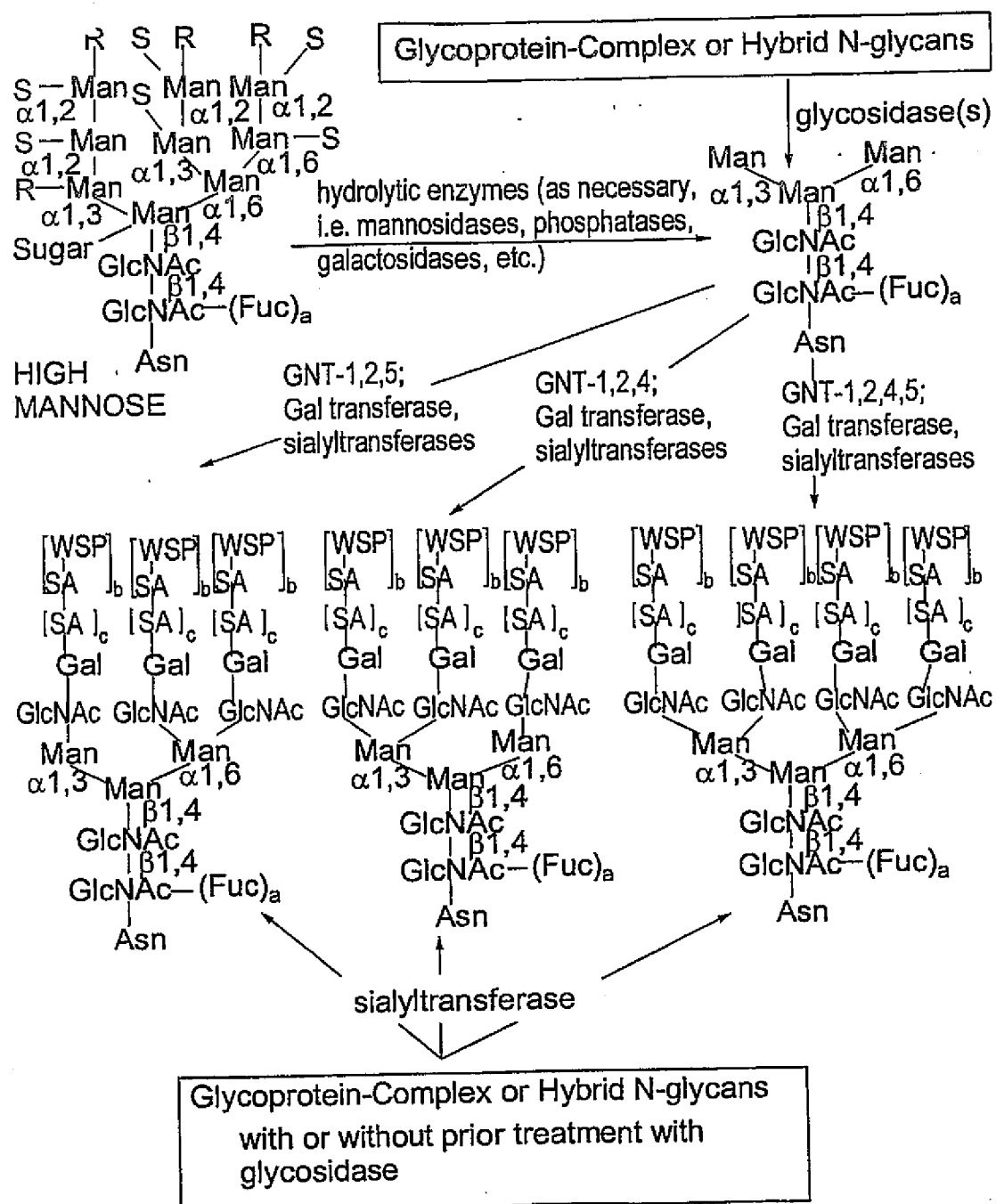


FIG. 20

23/497

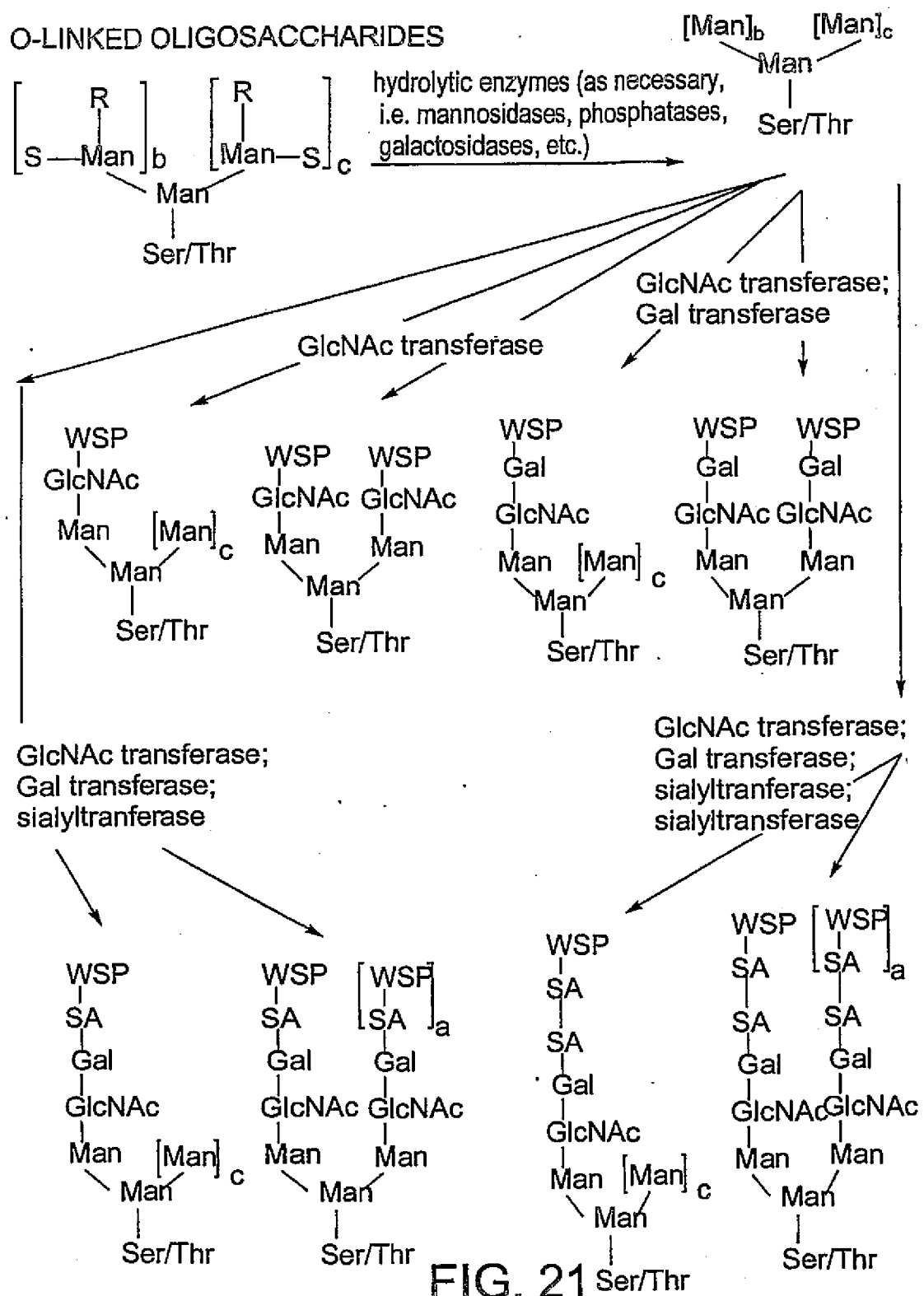


FIG. 21 Ser/Thr

24/497

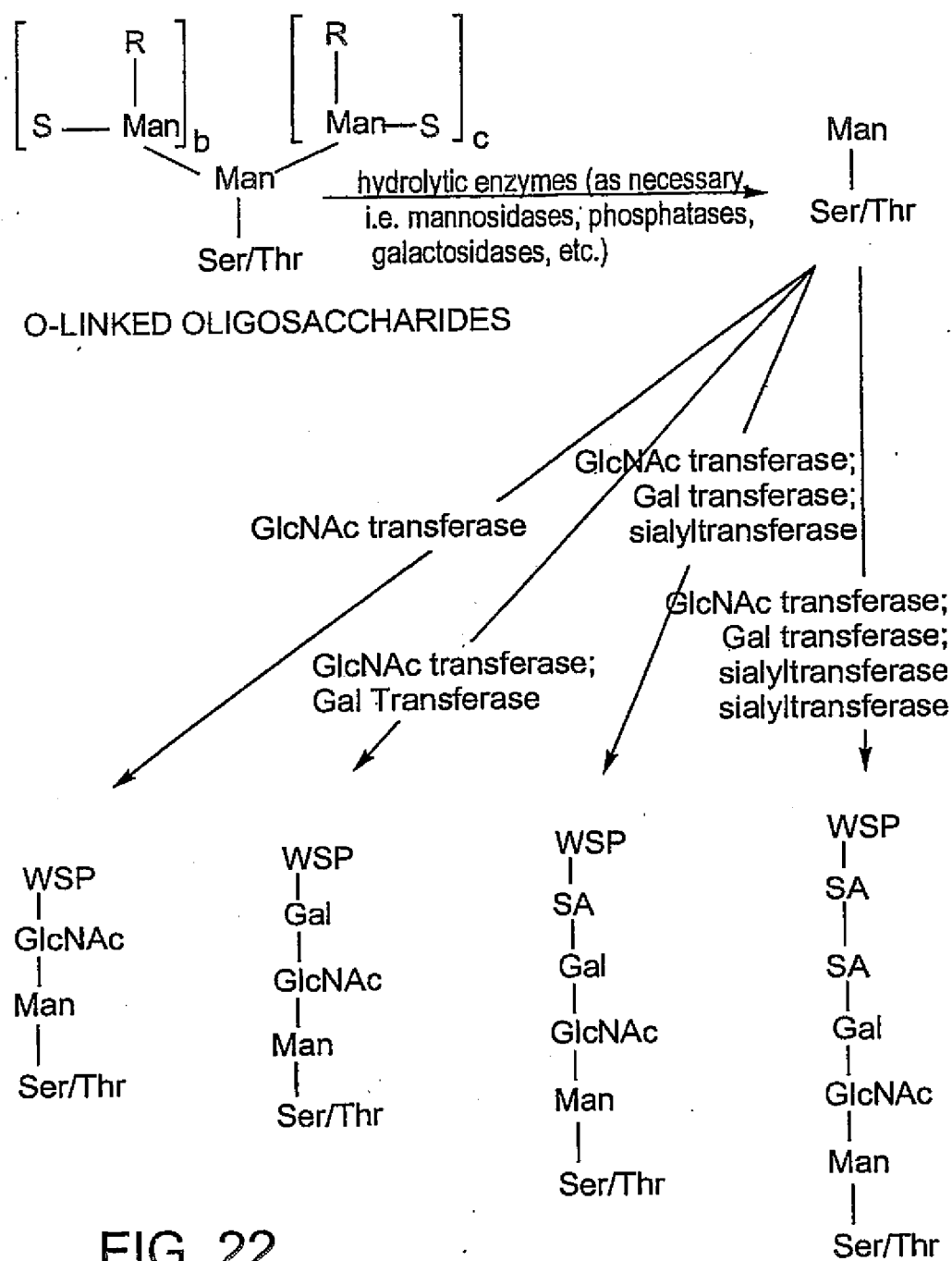


FIG. 22

25/497

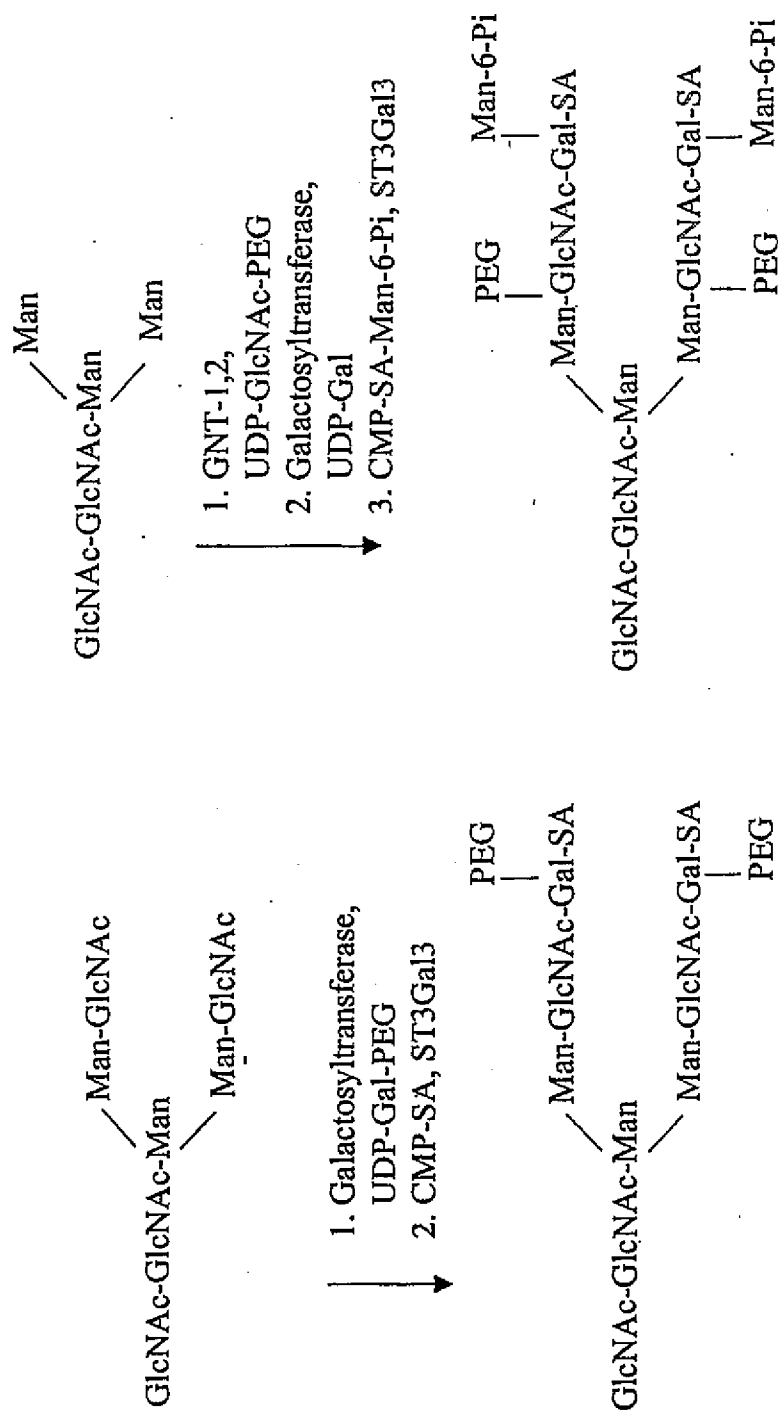


FIG. 23A

FIG. 23B

26/497

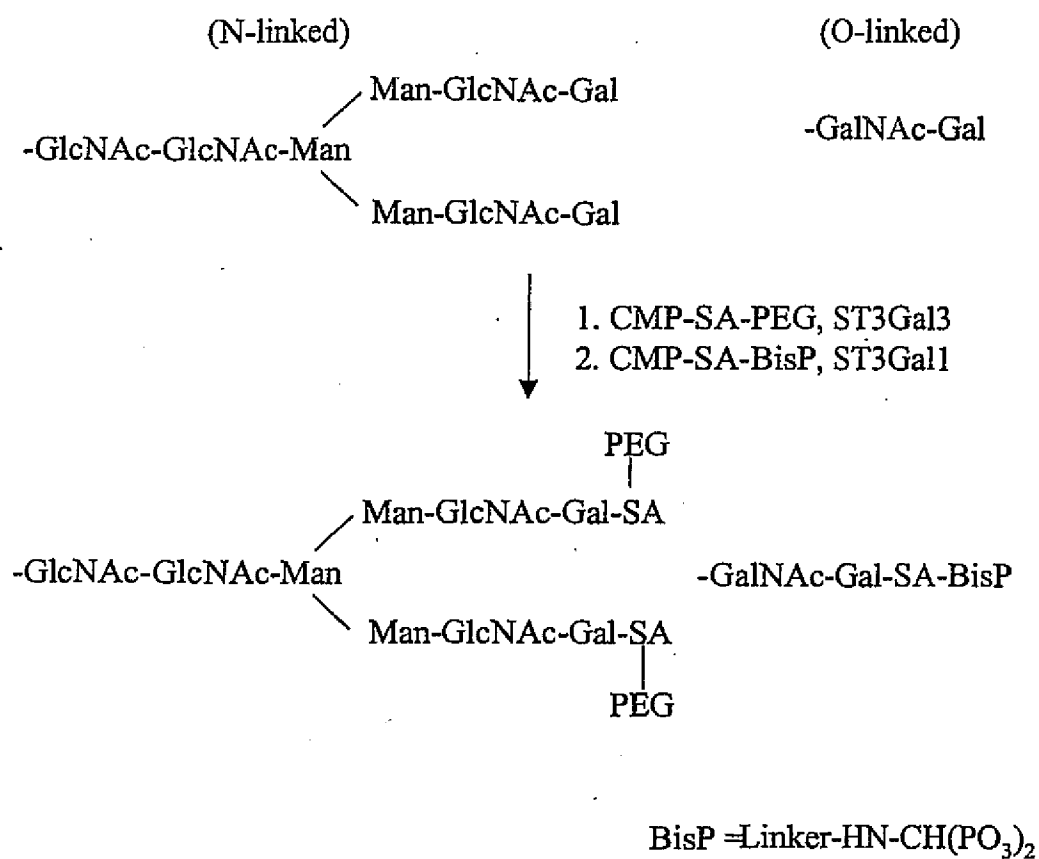


FIG. 23C

27/497

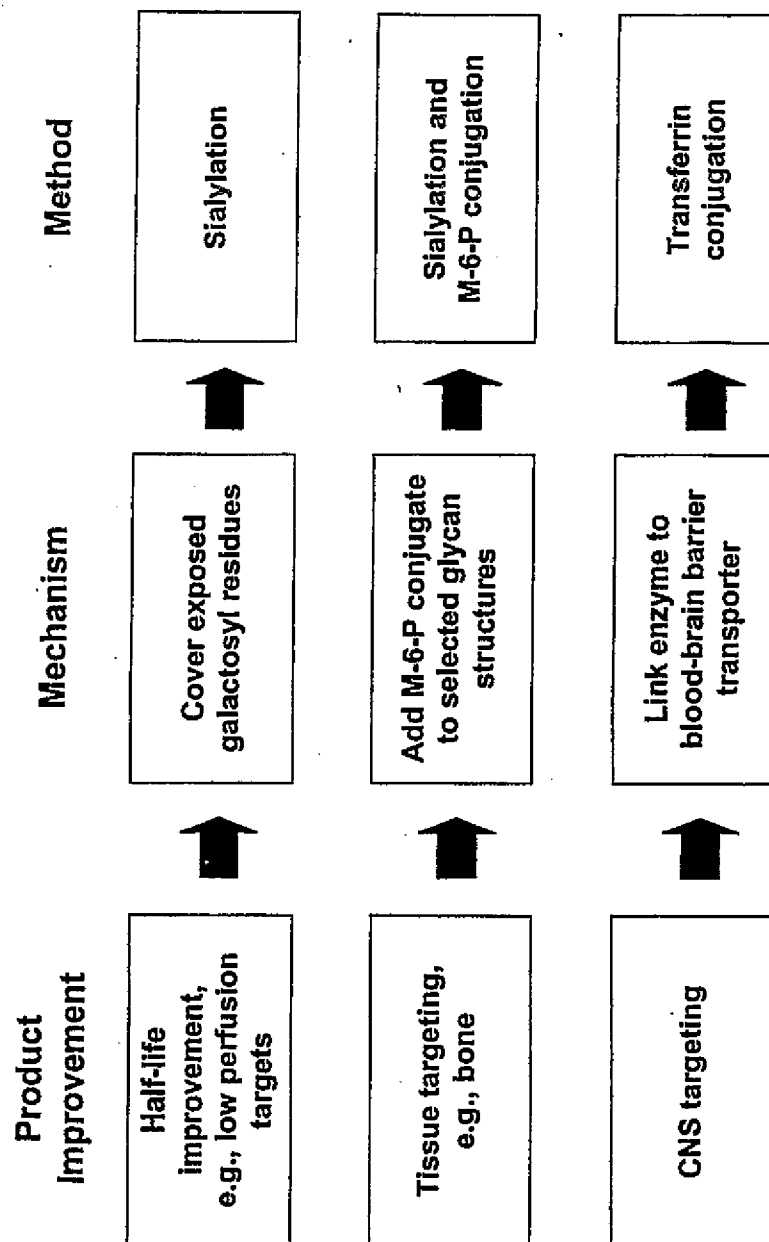


FIG. 24

28/497

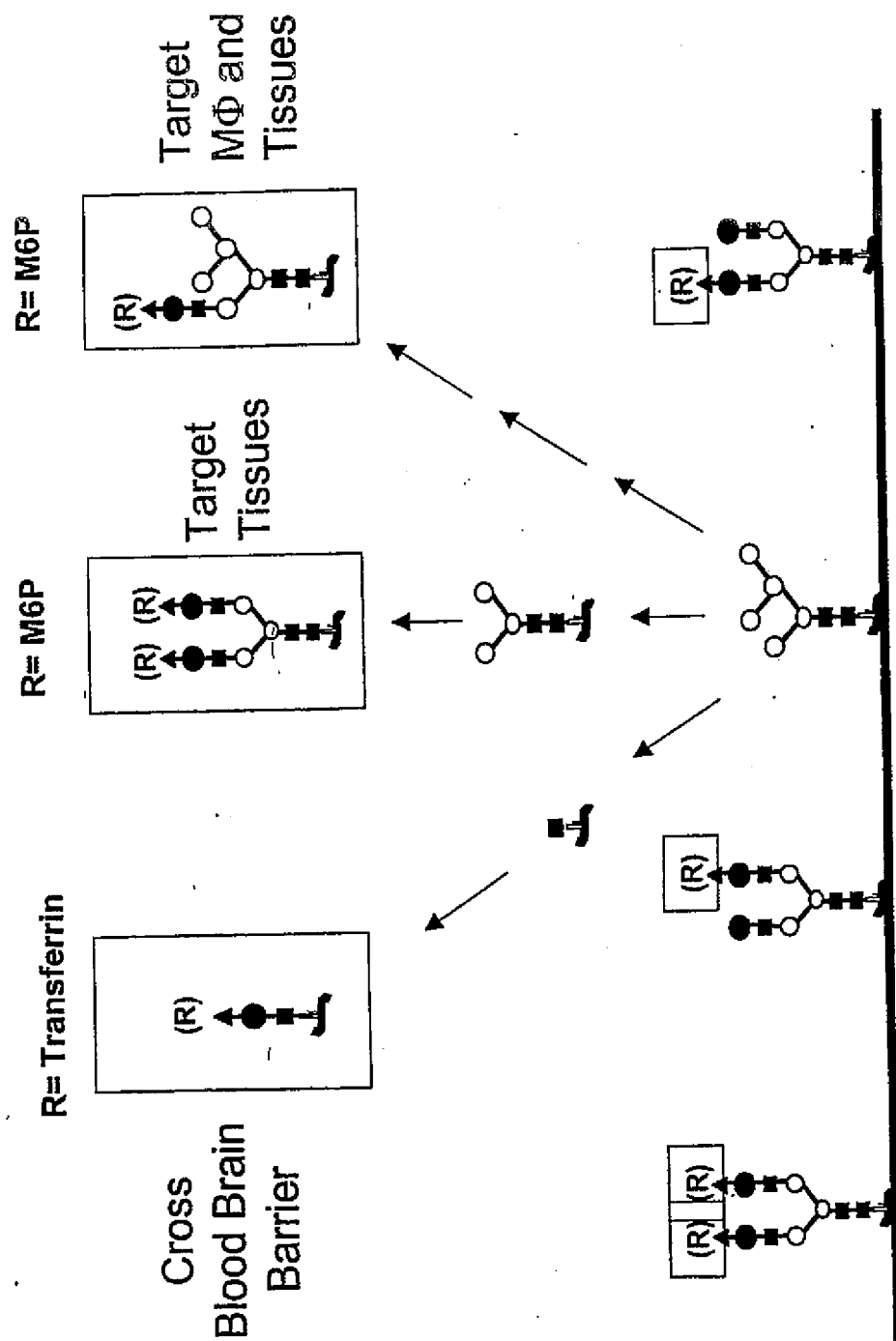


FIG. 25

29/497

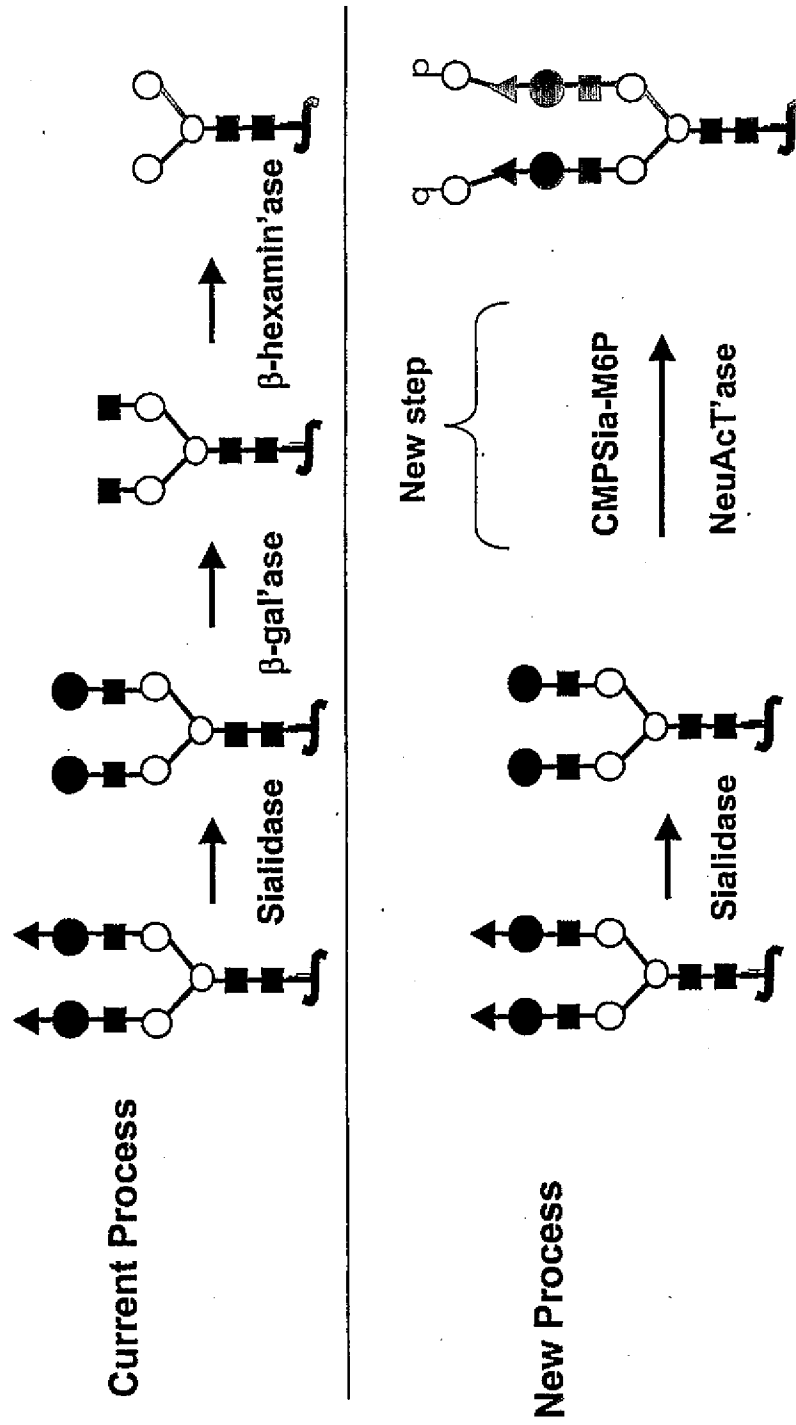


FIG. 26

30/497

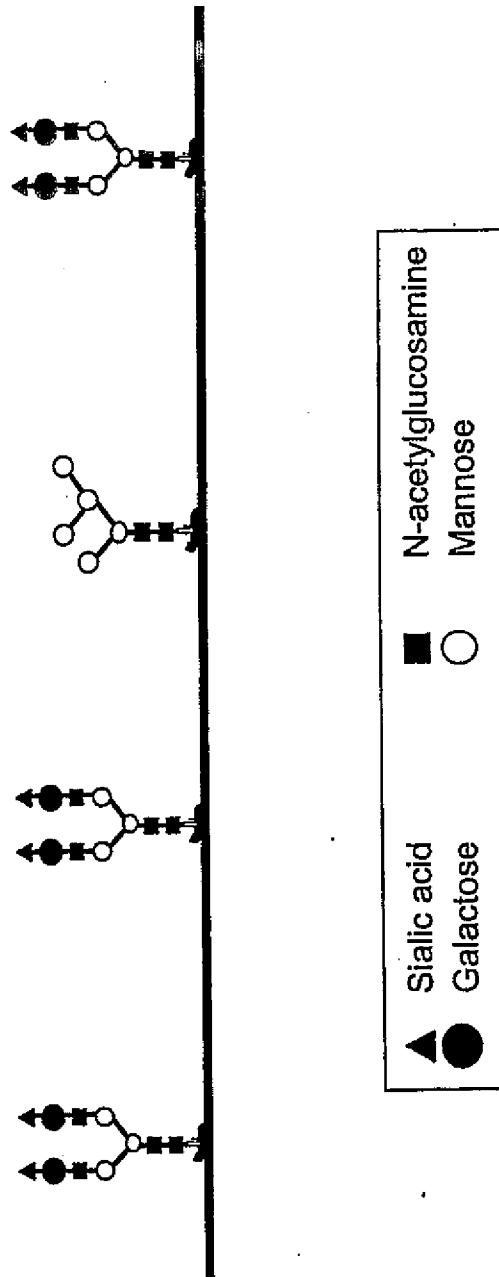


FIG. 27

31/497

12AP1/E5 -- Viventia Biotech	AI-201 -- AutoImmune
1964 -- Aventis	AI-301 -- AutoImmune
20K growth hormone -- AMUR	AIDS vaccine -- ANRS, CIBG, Hesed
28P6/E6 -- Viventia Biotech	Biomed, Hollis-Eden, Rome, United
3-Hydroxyphthaloyl-beta-lactoglobulin --	Biomedical, American Home Products,
4-IBB ligand gene therapy --	Maxygen
64-Cu MAb conjugate TETA-1A3 --	airway receptor ligand -- IC Innovations
Mallinckrodt Institute of Radiology	AJvW 2 -- Ajinomoto
64-Cu MAb conjugate TETA-cT84.66	AK 30 NGF -- Alkermes
64-Cu Trastuzumab TETA conjugate --	Albuferon -- Human Genome Sciences
Genentech	albumin -- Biogen, DSM Anti-Infectives,
A 200 -- Amgen	Genzyme Transgenics, PPL Therapeutics,
A10255 -- Eli Lilly	TranXenoGen, Welfide Corp.
A1PDX -- Hedral Therapeutics	aldesleukin -- Chiron
A6 -- Angstrom	alefacept -- Biogen
aaAT-III -- Genzyme	Alemtuzumab
Abciximab -- Centocor	Allergy therapy -- ALK-Abello/Maxygen,
ABI.001 -- Atlantic BioPharmaceuticals	ALK-Abello/RP Scherer
ABT-828 -- Abbott	allergy vaccines -- Allergy Therapeutics
Accutin	Alnidofibatide -- Aventis Pasteur
Actinohivin	Alnorine -- SRC VB VECTOR
activin -- Biotech Australia, Human	ALP 242 -- Gruenenthal
Therapeutics, Curis	Alpha antitrypsin -- Arriva/Hyland
AD 439 -- Tanox	Immuno/ProMetic/Protease Sciences
AD 519 -- Tanox	Alpha-1 antitrypsin -- Cutter, Bayer, PPL
Adalimumab -- Cambridge Antibody Tech.	Therapeutics, Profile, ZymoGenetics,
Adenocarcinoma vaccine -- Biomira -- NIS	Arriva
Adenosine deaminase -- Enzond	Alpha-1 protease inhibitor -- Genzyme
Adenosine A2B receptor antagonists --	Transgenics, Welfide Corp.
Adenosine Therapeutics	Alpha-galactose fusion protein --
ADP-001 -- Axis Genetics	Immunomedics
AF 13948 -- Affymax	Alpha-galactosidase A -- Research
Afelimomab -- Knoll	Corporation Technologies, Genzyme
AFP-SCAN -- Immunomedics	Alpha-glucosidase -- Genzyme, Novazyme
AG 2195 -- Corixa	Alpha-lactalbumin
agalsidase alfa -- Transkaryotic Therapies	Alpha-L-iduronidase -- Transkaryotic
agalsidase beta -- Genzyme	Therapies, BioMarin
AGENT-- Antisoma	alteplase -- Genentech
AI 300 -- AutoImmune	alvircept sudotox -- NIH
AI-101 -- Teva	ALX1-11 --sNPS Pharmaceuticals
AI-102 -- Teva	Alzheimer's disease gene therapy

FIG. 28A

32/497

AM-133 -- AMRAD	Anti-angiogenesis monoclonal antibodies --
Amb a 1 immunostim conj. -- Dynavax	KS Biomedix/Schering AG
AMD 3100 -- AnorMED -- NIS	Anti-B4 MAb-DC1 conjugate -- ImmunoGen
AMD 3465 -- AnorMED -- NIS	Anti-B7 antibody PRIMATIZED -- IDEC
AMD 3465 -- AnorMED -- NIS	Anti-B7-1 MAb 16-10A1
AMD Fab -- Genentech	Anti-B7-1 MAb 1G10
Amediplase -- Menarini, Novartis	Anti-B7-2 MAb GL-1
AM-F9	Anti-B7-2-gelonin immunotoxin --
Amoebiasis vaccine	Antibacterials/antifungals --
Amphiregulin -- Octagene	Diversa/IntraBiotics
anakinra -- Amgen	Anti-beta-amyloid monoclonal antibodies --
analgesic -- Nobex	Cambridge Antibody Tech., Wyeth-Ayerst
ancestim -- Amgen	Anti-BLyS antibodies -- Cambridge
AnergiX.RA -- Corixa, Organon	Antibody Tech. /Human Genome Sciences
Angiocidin -- InKine	Antibody-drug conjugates -- Seattle
angiogenesis inhibitors -- ILEX	Genetics/Eos
AngioMab -- Antisoma	Anti-C5 MAb BB5-1 -- Alexion
Angiopoietins -- Regeneron/Procter &	Anti-C5 MAb N19-8 -- Alexion
Gamble	Anti-C8 MAb
angiostatin -- EntreMed	anticancer cytokines -- BioPulse
Angiostatin/endostatin gene therapy --	anticancer matrix -- Telios Integra
Genetix Pharmaceuticals	Anticancer monoclonal antibodies -- ARIUS,
angiotensin-II, topical -- Maret	Immunex
Anthrax -- EluSys Therapeutics/US Army	anticancer peptides -- Maxygen, Micrologix
Medical Research Institute	Anticancer prodrug Tech. -- Alexion
Anthrax vaccine	Antibody Technologies
Anti platelet-derived growth factor D human	anticancer Troy-Bodies -- Affite -- Affitech
monoclonal antibodies -- CuraGen	anticancer vaccine -- NIH
Anti-17-1A MAb 3622W94 --	anticancers -- Epimmune
GlaxoSmithKline	Anti-CCR5/CXCR4 sheep MAb -- KS
Anti-2C4 MAb -- Genentech	Biomedix Holdings
anti-4-1BB monoclonal antibodies -- Bristol-	Anti-CD11a MAb KBA --
Myers Squibb	Anti-CD11a MAb M17
Anti-Adhesion Platform Tech. -- Cytovax	Anti-CD11a MAb TA-3 --
Anti-adipocyte MAb -- Cambridge Antibody	Anti-CD11a MAb WT.1 --
Tech./ObeSys	Anti-CD11b MAb -- Pharmacia
antiallergics -- Maxygen	Anti-CD11b MAb LM2
antiallergy vaccine -- Acambis	Anti-CD154 MAb -- Biogen
Anti-alpha-4-integrin MAb	Anti-CD16-anti-CD30 MAb -- Biotest
Anti-alpha-vb3 integrin MAb -- Applied	Anti-CD18 MAb -- Pharmacia
Molecular Evolution	Anti-CD19 MAb B43 --

FIG. 28B

33/497

Anti-CD19 MAb -liposomal sodium butyrate conjugate –	Anti-CD4 MAb 4162W94 – GlaxoSmithKline
Anti-CD147	Anti-CD4 MAb B-F5 – Diaclone
Anti-CD19 MAb-saporin conjugate –	Anti-CD4 MAb GK1-5
Anti-CD19-dsFv-PE38-immunotoxin –	Anti-CD4 MAb KT6
Anti-CD2 MAb 12-15 –	Anti-CD4 MAb OX38
Anti-CD2 MAb B-E2 – Diaclone	Anti-CD4 MAb PAP conjugate -- Bristol-Myers Squibb
Anti-CD2 MAb OX34 –	Anti-CD4 MAb RIB 5-2
Anti-CD2 MAb OX54 –	Anti-CD4 MAb W3/25
Anti-CD2 MAb OX55 –	Anti-CD4 MAb YTA 3.1.2
Anti-CD2 MAb RM2-1	Anti-CD4 MAb YTS 177-9
Anti-CD2 MAb RM2-2	Anti-CD40 ligand MAb 5c8 -- Biogen
Anti-CD2 MAb RM2-4	Anti-CD40 MAb
Anti-CD20 MAb BCA B20	Anti-CD40 MAb 5D12 – Tanox
Anti-CD20-anti-Fc alpha RI bispecific MAb – Medarex, Tenovus	Anti-CD44 MAb A3D8
Anti-CD22 MAb-saporin-6 complex –	Anti-CD44 MAb GKWA3
Anti-CD3 immunotoxin –	Anti-CD44 MAb IM7
Anti-CD3 MAb 145-2C11 – Pharming	Anti-CD44 MAb KM81
Anti-CD3 MAb CD4IgG conjugate -- Genentech	Anti-CD44 variant monoclonal antibodies -- Corixa/Hebrew University
Anti-CD3 MAb humanised – Protein Design, RW Johnson	Anti-CD45 MAb BC8-I-131
Anti-CD3 MAb WT32	Anti-CD45RB MAb
Anti-CD3 MAb-ricin-chain-A conjugate –	Anti-CD48 MAb HuLy-m3
Anti-CD3 MAb-xanthine-oxidase conjugate –	Anti-CD48 MAb WM-63
Anti-CD30 MAb BerH2 -- Medac	Anti-CD5 MAb – Becton Dickinson
Anti-CD30 MAb-saporin conjugate	Anti-CD5 MAb OX19
Anti-CD30-scFv-ETA'-immunotoxin	Anti-CD6 MAb
Anti-CD38 MAb AT13/5	Anti-CD7 MAb-PAP conjugate
Anti-CD38 MAb-saporin conjugate	Anti-CD7 MAb-ricin-chain-A conjugate
Anti-CD3-anti-CD19 bispecific MAb	Anti-CD8 MAb – Amerimmune, Cytodyn, Becton Dickinson
Anti-CD3-anti-EGFR MAb	Anti-CD8 MAb 2-43
Anti-CD3-anti-interleukin-2-receptor MAb	Anti-CD8 MAb OX8
Anti-CD3-anti-MOV18 MAb – Centocor	Anti-CD80 MAb P16C10 – IDEC
Anti-CD3-anti-SCLC bispecific MAb	Anti-CD80 MAb P7C10 – ID Vaccine
Anti-CD4 idiotype vaccine	Anti-CD8-idarubicin conjugate
Anti-CD4 MAb – Centocor, IDEC Pharmaceuticals, Xenova Group	Anti-CEA MAb CE-25
Anti-CD4 MAb 16H5	Anti-CEA MAb MN 14 – Immunomedics
	Anti-CEA MAb MN14-PE40 conjugate – Immunomedics

FIG. 28C

34/497

Anti-CEA MAb T84.66-interleukin-2 conjugate	Anti-heparanase human monoclonal antibodies -- Oxford
Anti-CEA sheep MAb -- KS Biomedix Holdings	Glycosciences/Medarex
Anti-cell surface monoclonal antibodies -- Cambridge Antibody Tech. /Pharmacia	Anti-hepatitis C virus human monoclonal antibodies -- XTL Biopharmaceuticals
Anti-c-erbB2-anti-CD3 bifunctional MAb -- Otsuka	Anti-HER-2 antibody gene therapy
Anti-CMV MAb -- Scotgen	Anti-herpes antibody -- Epicyte
Anti-complement	Anti-HIV antibody -- Epicyte
Anti-CTLA-4 MAb	anti-HIV catalytic antibody -- Hersed Biomed
Anti-EGFR catalytic antibody -- Hersed Biomed	anti-HIV fusion protein -- Idun
anti-EGFR immunotoxin -- IVAX	anti-HIV proteins -- Cangene
Anti-EGFR MAb -- Abgenix	Anti-HM1-24 MAb -- Chugai
Anti-EGFR MAb 528	Anti-hR3 MAb
Anti-EGFR MAb KSB 107 -- KS Biomedix	Anti-Human-Carcinoma-Antigen MAb -- Epicyte
Anti-EGFR MAb-DM1 conjugate -- ImmunoGen	Anti-ICAM-1 MAb -- Boehringer Ingelheim
Anti-EGFR MAb-LA1 --	Anti-ICAM-1 MAb 1A-29 -- Pharmacia
Anti-EGFR sheep MAb -- KS Biomedix	Anti-ICAM-1 MAb HA58
Anti-FAP MAb F19-I-131	Anti-ICAM-1 MAb YN1/1.7.4
Anti-Fas IgM MAb CH11	Anti-ICAM-3 MAb ICM3 -- ICOS
Anti-Fas MAb Jo2	Anti-idiotypic breast cancer vaccine 11D10
Anti-Fas MAb RK-8	Anti-idiotypic breast cancer vaccine ACA14C5 --
Anti-Fit-1 monoclonal antibodies -- ImClone	Anti-idiotypic cancer vaccine -- ImClone Systems/Merck KGaA ImClone, Viventia Biotech
Anti-fungal peptides -- State University of New York	Anti-idiotypic cancer vaccine 1A7 -- Titan
antifungal tripeptides -- BTG	Anti-idiotypic cancer vaccine 3H1 -- Titan
Anti-ganglioside GD2 antibody-interleukin-2 fusion protein -- Lexigen	Anti-idiotypic cancer vaccine TriAb -- Titan
Anti-GM2 MAb -- Kyowa	Anti-idiotypic Chlamydia trachomatis vaccine
Anti-GM-CSF receptor monoclonal antibodies -- AMRAD	Anti-idiotypic colorectal cancer vaccine -- Novartis
Anti-gp130 MAb -- Tosoh	Anti-idiotypic colorectal cancer vaccine -- Onyvax
Anti-HCA monoclonal antibodies -- AltaRex/Epigen	Anti-idiotypic melanoma vaccine -- IDEC Pharmaceuticals
Anti-hCG antibodies -- Abgenix/AVI BioPharma	Anti-idiotypic ovarian cancer vaccine ACA 125
	Anti-idiotypic ovarian cancer vaccine AR54 - AltaRex

FIG. 28D

35/497

Anti-idiotypic ovarian cancer vaccine CA-125 -- AltaRex, Biomira	Anti-L-selectin monoclonal antibodies -- Protein Design Labs, Abgenix, Stanford University
Anti-IgE catalytic antibody -- Hersed Biomed	Anti-MBL monoclonal antibodies -- Alexion/Brigham and Women's Hospital
Anti-IgE MAb E26 -- Genentech	Anti-MHC monoclonal antibodies
Anti-IGF-1 MAb	Anti-MIF antibody humanised -- IDEC, Cytokine PharmaSciences
anti-inflammatory -- GeneMax	Anti-MRSA/VRSA sheep MAb -- KS Biomedix Holdings
anti-inflammatory peptide -- BTG	Anti-mu MAb -- Novartis
anti-integrin peptides -- Burnha	Anti-MUC-1 MAb
Anti-interferon-alpha-receptor MAb 64G12 -- Pharma Pacific Management	Anti-MUC 18
Anti-interferon-gamma MAb -- Protein Design Labs	Anti-Nogo-A MAb IN1
Anti-interferon-gamma polyclonal antibody - Advanced Biotherapy	Anti-nuclear autoantibodies -- Procyon
Anti-interleukin-10 MAb --	Anti-ovarian cancer monoclonal antibodies - Dompe
Anti-interleukin-12 MAb --	Anti-p185 monoclonal antibodies
Anti-interleukin-1-beta polyclonal antibody -- R&D Systems	Anti-p43 MAb
Anti-interleukin-2 receptor MAb 2A3	Antiparasitic vaccines
Anti-interleukin-2 receptor MAb 33B3-1 -- Immunotech	Anti-PDGF/bFGF sheep MAb -- KS Biomedix
Anti-interleukin-2 receptor MAb ART-18	Anti-properdin monoclonal antibodies -- Abgenix/Giatech
Anti-interleukin-2 receptor MAb LO-Tact-1	Anti-PSMA (prostate specific membrane antigen)
Anti-interleukin-2 receptor MAb Mikbeta1	Anti-PSMA MAb J591 -- BZL Biologics
Anti-interleukin-2 receptor MAb NDS61	Anti-Rev MAb gene therapy --
Anti-interleukin-4 MAb 11B11	Anti-RSV antibodies -- Epicyte, Intracell
Anti-interleukin-5 MAb -- Wallace Laboratories	Anti-RSV monoclonal antibodies -- Medarex/MedImmune, Applied Molecular Evolution/MedImmune
Anti-interleukin-6 MAb -- Centocor, Diaclone, Pharmadigm	Anti-RSV MAb, inhalation -- Alkermes/MedImmune
Anti-interleukin-8 MAb -- Abgenix	Anti-RT gene therapy
Anti-interleukin-8 MAb -- Xenotech	Antisense K-ras RNA gene therapy
Anti-JL1 MAb	Anti-SF-25 MAb
Anti-Klebsiella sheep MAb -- KS Biomedix Holdings	Anti-sperm antibody -- Epicyte
Anti-Laminin receptor MAb-liposomal doxorubicin conjugate	Anti-Tac(Fv)-PE38 conjugate
Anti-LCG MAb -- Cytoclonal	Anti-TAPA/CD81 MAb AMP1
Anti-lipopolysaccharide MAb -- VitaResc	Anti-tat gene therapy

FIG. 28E

36/497

Anti-TCR-alphabeta MAb H57-597	AOP-RANTES -- Senetek
Anti-TCR-alphabeta MAb R73	Apan-CH -- Praecis Pharmaceuticals
Anti-tenascin MAb BC-4-I-131	APC-8024 -- Demegen
Anti-TGF-beta human monoclonal antibodies -- Cambridge Antibody Tech., Genzyme	ApoA-1 -- Milano, Pharmacia
Anti-TGF-beta MAb 2G7 -- Genentech	Apogen -- Alexion
Antithrombin III -- Genzyme Transgenics, Aventis, Bayer, Behringwerke, CSL, Myriad	apolipoprotein A1 -- Avanir
Anti-Thy1 MAb	Apolipoprotein E -- Bio-Tech. General
Anti-Thy1.1 MAb	Applaggin -- Biogen
Anti-tissue factor/factor VIIa sheep MAb -- KS Biomedix	aprotinin -- ProdiGene
Anti-TNF monoclonal antibodies -- Centocor, Chiron, Peptech, Pharacia, Serono	APT-070C -- AdProTech
Anti-TNF sheep MAb -- KS Biomedix Holdings	AR 177 -- Aronex Pharmaceuticals
Anti-TNFalpha MAb -- Genzyme	AR 209 -- Aronex Pharmaceuticals, Antigenics
Anti-TNFalpha MAb B-C7 -- Diacione	AR545C
Anti-tooth decay MAb -- Planet BioTech.	ARGENT gene delivery systems -- ARIAD
Anti-TRAIL receptor-1 MAb -- Takeda	Arresten
Antitumour RNases -- NIH	ART-123 -- Asahi Kasei
Anti-VCAM MAb 2A2 -- Alexion	arylsulfatase B -- BioMarin
Anti-VCAM MAb 3F4 -- Alexion	Arylsulfatase B, Recombinant human -- BioMarin
Anti-VCAM-1 MAb	AS 1051 -- Ajinomoto
Anti-VEC MAb -- ImClone	ASI-BCL -- Intracell
Anti-VEGF MAb -- Genentech	Asparaginase - Merck
Anti-VEGF MAb 2C3	ATL-101 -- Alizyme
Anti-VEGF sheep MAb -- KS Biomedix Holdings	Atrial natriuretic peptide -- Pharis
Anti-VLA-4 MAb HP1/2 -- Biogen	Aurintricarboxylic acid-high molecular weight
Anti-VLA-4 MAb PS/2	Autoimmune disorders -- GPC
Anti-VLA-4 MAb R1-2	Biotech/MorphoSys
Anti-VLA-4 MAb TA-2	Autoimmune disorders and transplant rejection -- Bristol-Myers Squibb/Genzyme
Anti-VAP-1 human MAb	Tra
Anti-VRE sheep MAb -- KS Biomedix Holdings	Autoimmune disorders/cancer -- Abgenix/Chiron, CuraGen
ANUP -- TranXenoGen	Autotaxin
ANUP-1 -- Pharis	Avicidin -- NeoRx
	axogenesis factor-1 -- Boston Life Sciences
	Axokine -- Regeneron
	B cell lymphoma vaccine -- Biomira
	B7-1 gene therapy --
	BABS proteins -- Chiron

FIG. 28F

37/497

BAM-002 -- Novelos Therapeutics	BMP 2 -- Genetics Institute/Medtronic-
Basiliximab (anti CD25 MAb) -- Novartis	Sofamor Danek, Genetics Institute/
Bay-16-9996 -- Bayer	Collagenesis, Genetics
Bay-39-9437 -- Bayer	Institute/Yamanouch
Bay-50-4798 -- Bayer	BMP 2 gene therapy
BB-10153 -- British Biotech	BMP 52 -- Aventis Pasteur, Biopharm
BBT-001 -- Bolder BioTech.	BMP-2 -- Genetics Institute
BBT-002 -- Bolder BioTech.	BMS 182248 -- Bristol-Myers Squibb
BBT-003 -- Bolder BioTech.	BMS 202448 -- Bristol-Myers Squibb
BBT-004 -- Bolder BioTech.	bone growth factors -- IsoTis
BBT-005 -- Bolder BioTech.	BPC-15 -- Pfizer
BBT-006 -- Bolder BioTech.	brain natriuretic peptide --
BBT-007 -- Bolder BioTech.	Breast cancer -- Oxford
BCH-2763 -- Shire	GlycoSciences/Medarex
BCSF -- Millenium Biologix	Breast cancer vaccine -- Therion Biologics,
BDNF -- Regeneron -- Amgen	Oregon
Becaplermin -- Johnson & Johnson, Chiron	BSSL -- PPL Therapeutics
Bectumomab -- Immunomedics	BST-2001 -- BioStratum
Beriplast -- Aventis	BST-3002 -- BioStratum
Beta-adrenergic receptor gene therapy --	BTI 322 --
University of Arkansas	butyrylcholinesterase -- Shire
bFGF -- Scios	C 6822 -- COR Therapeutics
BI 51013 -- Behringwerke AG	C1 esterase inhibitor -- Pharming
BIBH 1 -- Boehringer Ingelheim	C3d adjuvant -- AdProTech
BIM-23190 -- Beaufour-Ipsen	CAB-2.1 -- Millennium
birch pollen immunotherapy -- Pharmacia	calcitonin -- Inhale Therapeutics Systems,
bispecific fusion proteins -- NIH	Aventis, Genetronics, TranXenoGen,
Bispecific MAb 2B1 -- Chiron	Unigene, Rhone Poulenc Rohrer
Bitistatin	calcitonin -- oral -- Nobex, Emisphere,
BIWA 4 -- Boehringer Ingelheim	Pharmaceutical Discovery
blood substitute -- Northfield, Baxter Intl.	Calcitonin gene-related peptide -- Asahi
BLP-25 -- Biomira	Kasei -- Unigene
BLS-0597 -- Boston Life Sciences	calcitonin, human -- Suntory
BLyS -- Human Genome Sciences	calcitonin, nasal -- Novartis, Unigene
BLyS radiolabelled -- Human Genome	calcitonin, Panoderm -- Elan
Sciences	calcitonin, Peptitrol -- Shire
BM 06021 -- Boehringer Mannheim	calcitonin, salmon -- Therapicon
BM-202 -- BioMarin	calin -- Biopharm
BM-301 -- BioMarin	Calphobindin I
BM-301 -- BioMarin	calphobindin I -- Kowa
BM-302 -- BioMarin	calreticulin -- NYU

FIG. 28G

38/497

Campath-1G	CD4 fusion toxin -- Senetek
Campath-1M	CD4 IgG -- Genentech
cancer therapy -- Cangene	CD4 receptor antagonists --
cancer vaccine -- Aixlie, Aventis Pasteur,	PharmacoPeia/Progenics
Center of Molecular Immunology, YM	CD4 soluble -- Progenics
BioSciences, Cytos, Genzyme,	CD4, soluble -- Genzyme Transgenics
Transgenics, GlobelImmune, Igeneon,	CD40 ligand -- Immunex
ImClone, Virogenetics, InterCell, Iomai,	CD4-ricin chain A -- Genentech
Jenner Biotherapies, Memorial Sloan-	CD59 gene therapy -- Alexion
Kettering Cancer Center, Sydney Kimmel	CD8 TIL cell therapy -- Aventis Pasteur
Cancer Center, Novavax, Protein	CD8, soluble -- Avidex
Sciences, Argonex, SIGA	CD95 ligand -- Roche
Cancer vaccine ALVAC-CEA B7.1 --	CDP 571 -- Celltech
Aventis Pasteur/Therion Biologics	CDP 850 -- Celltech
Cancer vaccine CEA-TRICOM -- Aventis	CDP-860 (PEG-PDGF MAb) -- Celltech
Pasteur/Therion Biologics	CDP 870 -- Celltech
Cancer vaccine gene therapy -- Cantab	CDS-1 -- Ernest Orlando
Pharmaceuticals	Cedelizumab -- Ortho-McNeil
Cancer vaccine HER-2/neu -- Corixa	Cetermin -- Insmed
Cancer vaccine THERATOPE -- Biomira	CETP vaccine -- Avant
cancer vaccine, PolyMASC -- Valentis	Cetrorelix
Candida vaccine -- Corixa, Inhibitex	Cetuximab
Canstatin -- ILEX	CGH 400 -- Novartis
CAP-18 -- Panorama	CGP 42934 -- Novartis
Cardiovascular gene therapy -- Collateral	CGP 51901 -- Tanox
Therapeutics	CGRP -- Unigene
carperitide -- Suntory	CGS 27913 -- Novartis
Casocidin-1 -- Pharis	CGS 32359 -- Novartis
CAT 152 -- Cambridge Antibody Tech.	Chagas disease vaccine -- Corixa
CAT 192 -- Cambridge Antibody Tech.	chemokines -- Immune Response
CAT 213 -- Cambridge Antibody Tech.	CHH 380 -- Novartis
Catalase -- Enzon	chitinase -- Genzyme, ICOS
Cat-PAD -- Circassia	Chlamydia pneumoniae vaccine -- Antex
CB 0006 -- Celltech	Biologics
CCK(27-32) -- Akzo Nobel	Chlamydia trachomatis vaccine -- Antex
CCR2-64I -- NIH	Biologics
CD, Procept -- Paligent	Chlamydia vaccine -- GlaxoSmithKline
CD154 gene therapy	Cholera vaccine CVD 103-HgR -- Swiss
CD39 -- Immunex	Serum and Vaccine Institute Berne
CD39-L2 -- Hyseq	Cholera vaccine CVD 112 -- Swiss Serum
CD39-L4 -- Hyseq	and Vaccine Institute Berne

FIG. 28H

39/497

Cholera vaccine inactivated oral – SBL	CRL 1605 – CytRx
Vaccin	CS-560 – Sankyo
Chrysalin – Chrysalis BioTech.	CSF – ZymoGenetics
CI-782 – Hitachi Kase	CSF-G – Hangzhou, Dong-A, Hanmi
Ciliary neurotrophic factor – Fidia, Roche	CSF-GM – Cangene, Hunan, LG Chem
CIM project – Active Biotech	CSF-M – Zarix
CL 329753 – Wyeth-Ayerst	CT 1579 – Merck Frosst
CL22, Cobra – ML Laboratories	CT 1786 – Merck Frosst
Clenoliximab – IDEC	CT-112 [^] – BTG
Clostridium difficile antibodies – Epicyte	CTB-134L – Xenova
clotting factors – Octagene	CTC-111 – Kaketsuken
CMB 401 – Celltech	CTGF – FibroGen
CNTF – Sigma-Tau	CTLA4-Ig – Bristol-Myers Squibb
Cocaine abuse vaccine – Cantab,	CTLA4-Ig gene therapy –
ImmuLogic, Scripps	CTP-37 – AVI BioPharma
coccidiomycosis vaccine – Arizo	C-type natriuretic peptide – Suntory
collagen – Type I – Pharming	CVS 995 – Corvas Intl.
Collagen formation inhibitors – FibroGen	CX 397 – Nikko Kyodo
Collagen/hydroxyapatite/bone growth factor	CY 1747 – Epimmune
– Aventis Pasteur, Biopharm, Orquest	CY 1748 – Epimmune
collagenase – BioSpecifics	Cyanovirin-N
Colorectal cancer vaccine – Wistar Institute	Cystic fibrosis therapy – CBR/IVAX
Component B, Recombinant – Serono	CYT 351
Connective tissue growth factor inhibitors –	cytokine Traps – Regeneron
FibroGen/Taisho	cytokines – Enzon, Cytoclonal
Contortrostatin	Cytomegalovirus glycoprotein vaccine –
contraceptive vaccine – Zonagen	Chiron, Aquila Biopharmaceuticals,
Contraceptive vaccine hCG	Aventis Pasteur, Virogenetics
Contraceptive vaccine male reversible –	Cytomegalovirus vaccine live – Aventis
IMMUCON	Pasteur
Contraceptive vaccine zona pellucida –	Cytosine deaminase gene therapy –
Zonagen	GlaxoSmithKline
Copper-64 labelled MAb TETA-1A3 – NCI	DA-3003 – Dong-A
Coralyne	DAB389interleukin-6 – Senetek
Corsevin M	DAB389interleukin-7
C-peptide analogues – Schwarz	Daclizumab (anti-IL2R MAb) – Protein
CPI-1500 – Consensus	Design Labs
CRF – Neurobiological Tech.	DAMP [^] – Incyte Genomics
cRGDFV pentapeptide –	Daniplestim – Pharmacia
CRL 1095 – CytRx	darbepoetin alfa – Amgen
CRL 1336 – CytRx	DBI-3019 – Diabetogen

FIG. 28I

40/497

DCC -- Genzyme	Duteplase -- Baxter Intl.
DDF -- Hyseq	DWP-401 -- Daewoong
decorin -- Integra, Telios	DWP-404 -- Daewoong
defensins -- Large Scale Biology	DWP-408 -- Daewoong
DEGR-VIIa	Dx 88 (Epi-KAL2) -- Dyax
DeImmunised antibody 3B6/22 AGEN	Dx 890 (elastin inhibitors) -- Dyax
Deimmunised anti-cancer antibodies -- Biovation/Viragen	E coli O157 vaccine -- NIH
Dendroamide A	E21-R -- BresaGen
Dengue vaccine -- Bavarian Nordic, Merck	Eastern equine encephalitis virus vaccine --
denileukin diftitox -- Ligand	Echicetin --
DES-1101 -- Desmos	Echinhibin 1 --
desirudin -- Novartis	Echistatin -- Merck
desmopressin -- Unigene	Echitamine --
Desmoteplase -- Merck, Schering AG	Ecromeximab -- Kyowa Hakko
Destabilase	EC-SOD -- PPL Therapeutics
Diabetes gene therapy -- DeveloGen, Pfizer	Eculizumab (5G1.1) -- Alexion
Diabetes therapy -- Crucell	EDF -- Ajinomoto
Diabetes type 1 vaccine -- Diamyd Therapeutics	EDN derivative -- NIH
DiaCIM -- YM BioSciences	EDNA -- NIH
dialytic oligopeptides -- Research Corp	Edobacomab -- XOMA
Diamyd -- Diamyd Therapeutics	Edrecolomab -- Centocor
DiaPep227 -- Pepgen	EF 5077
DiavaX -- Corixa	Efalizumab -- Genentech
Digoxin MAb -- Glaxo	EGF fusion toxin -- Seragen, Ligand
Diphtheria tetanus pertussis-hepatitis B vaccine -- GlaxoSmithKline	EGF-P64k vaccine -- Center of Molecular Immunology
DIR therapy -- Solis Therapeutics --	EL 246 -- LigoCyte
DNase -- Genentech	elastase inhibitor -- Synergen
Dornase alfa -- Genentech	elcatonin -- Therapicon
Dornase alfa, inhalation -- Genentech	EMD 72000 -- Merck KGaA
Doxorubicin-anti-CEA MAb conjugate -- Immunomedics	Emdogain -- BIORA
DP-107 -- Trimeris	emfilermin -- AMRAD
drotrecogin alfa -- Eli Lilly	Emoctakin -- Novartis
DTctGMCSF	enamel matrix protein -- BIORA
DTP-polio vaccine -- Aventis Pasteur	Endo III -- NYU
DU 257-KM231 antibody conjugate -- Kyowa	endostatin -- EntreMed, Pharis
dural graft matrix -- Integra	Enhancins -- Micrologix
	Enlimomab -- Isis Pharm.
	Enoxaparin sodium -- Pharmuka
	enzyme linked antibody nutrient depletion therapy -- KS Biomedix Holdings

FIG. 28J

41/497

Eosinophil-derived neutralizing agent –	Factor VII – Novo Nordisk, Bayer, Baxter
EP-51216 – Asta Medica	Intl.
EP-51389 – Asta Medica	Factor VIIa – PPL Therapeutics,
EPH family ligands – Regeneron	ZymoGenetics
Epidermal growth factor – Hitachi Kasei,	Factor VIII – Bayer Genentech, Beaufour-
Johnson & Johnson	Ipsen, CLB, Inex, Octagen, Pharmacia,
Epidermal growth factor fusion toxin –	Pharming
Senetek	Factor VIII – PEGylated – Bayer
Epidermal growth factor-genistein –	Factor VIII fragments – Pharmacia
EPI-HNE-4 – Dyax	Factor VIII gene therapy – Targeted
EPI-KAL2 – Dyax	Genetics
Epoetin-alfa – Amgen, Dragon	Factor VIII sucrose formulation – Bayer,
Pharmaceuticals, Nanjing Huaxin	Genentech
Epratuzumab – Immunomedics	Factor VIII-2 – Bayer
Epstein-Barr virus vaccine –	Factor VIII-3 – Bayer
Aviron/SmithKline Beecham, Bioresearch	Factor Xa inhibitors – Merck, Novo Nordisk,
Eptacog alfa – Novo Nordisk	Mochida
Eptifibatide – COR Therapeutics	Factor XIII – ZymoGenetics
erb-38 –	Factors VIII and IX gene therapy – Genetics
Erlizumab – Genentech	Institute/Targeted Genetics
erythropoietin – Alkermes, ProLease, Dong-	Famoxin – Genset
A, Elanex, Genetics Institute, LG Chem,	Fas (delta) TM protein – LXR BioTech.
Protein Sciences, Serono, Snow Brand,	Fas TR – Human Genome Sciences
SRC VB VECTOR, Transkaryotic	Felvizumab – Scotgen
Therapies	FFR-VIIa – Novo Nordisk
Erythropoietin Beta – Hoffman La Roche	FG-001 – F-Gene
Erythropoietin/Epoetin alfa – Chugai	FG-002 – F-Gene
Escherichia coli vaccine – North American	FG-004 – F-Gene
Vaccine, SBL Vaccin, Swiss Serum and	FG-005 – F-Gene
Vaccine Institute Berne	FGF + fibrin – Repair
etanercept – Immunex	Fibrimage – Bio-Tech. General
examorelin – Mediolanum	fibrin-binding peptides – ISIS Innovation
Exendin 4 – Amylin	fibrinogen – PPL Therapeutics, Pharming
exonuclease VII	fibroblast growth factor – Chiron, NYU,
F 105 – Centocor	Ramot, ZymoGenetics
F-992 – Fornix	fibrolase conjugate – Schering AG
Factor IX – Alpha Therapeutics, Welfide	Filgrastim – Amgen
Corp., CSL, enetics Institute/AHP,	filgrastim – PDA modified – Xencor
Pharmacia, PPL Therapeutics	FLT-3 ligand – Immunex
Factor IX gene therapy – Cell Genesys	FN18 CRM9 –

FIG. 28K

42/497

follistatin -- Biotech Australia, Human Therapeutics	glutamate decarboxylase -- Genzyme Transgenics
follitropin alfa -- Alkermes, ProLease, PowderJect, Serono, Akzo Nobel	Glycoprotein S3 -- Kureha
Follitropin Beta -- Bayer, Organon	GM-CSF -- Immunex
FP 59	GM-CSF tumour vaccine -- PowderJect
FSH -- Ferring	GnRH immunotherapeutic -- Protherics
FSH + LH -- Ferring	Goserelin (LhRH antagonist) -- AstraZeneca
F-spondin -- CeNeS	gp75 antigen -- ImClone
fusion protein delivery system -- UAB Research Foundation	gp96 -- Antigenics
fusion toxins -- Boston Life Sciences	GPI 0100 -- Galenica
G 5598 -- Genentech	GR 4991W93 -- GlaxoSmithKline
GA-II -- Transkaryotic Therapies	Granulocyte colony-stimulating factor -- Dong-A
Gamma-interferon analogues -- SRC VB VECTOR	Granulocyte colony-stimulating factor conjugate
Ganirelix -- Roche	grass allergy therapy -- Dynavax
gastric lipase -- Meristem	GRF1-44 -- ICN
Gavilimomab --	Growth Factor -- Chiron, Atrigel, Atrix, Innogenetics, ZymoGenetics, Novo
G-CSF -- Amgen, SRC VB VECTOR	growth factor peptides -- Biotherapeutics
GDF-1 -- CeNeS	growth hormone -- LG Chem
GDF-5 -- Biopharm	growth hormone, Recombinant human -- Serono
GDNF (glial derived neurotrophic factor) -- Amgen	GT 4086 -- Gliatech
gelsolin -- Biogen	GW 353430 -- GlaxoSmithKline
Gemtuzumab ozogamicin -- Celltech	GW-278884 -- GlaxoSmithKline
Gene-activated epoetin-alfa -- Aventis Pharma -- Transkaryotic Therapies	H 11 -- Viventia Biotech
Glanzmann thrombasthenia gene therapy --	H5N1 influenza A virus vaccine -- Protein Sciences
Glatiramer acetate -- Yeda	haemoglobin -- Biopure
glial growth factor 2 -- CeNeS	haemoglobin 3011, Recombinant -- Baxter Healthcare
GLP-1 -- Amylin, Suntory, TheraTech, Watson	haemoglobin crosfumaril -- Baxter Intl.
GLP-1 peptide analogues -- Zealand Pharmaceuticals	haemoglobin stabilized -- Ajinomoto
glucagon -- Eli Lilly, ZymoGenetics	haemoglobin, recombinant -- Apex
Glucagon-like peptide-1 7-36 amide -- Suntory	HAF -- Immune Response
Glucogen-like peptide -- Amylin	Hantavirus vaccine
Glucocerebrosidase -- Genzyme	HB 19
	HBNF -- Regeneron
	HCC-1 -- Pharis
	hCG -- Milkhaus

FIG. 28L

43/497

hCG vaccine -- Zonagen	Herpes simplex glycoprotein DNA vaccine --
HE-317 -- Hollis-Eden Pharmaceuticals	Merck, Wyeth-Lederle Vaccines-Malvern,
Heat shock protein cancer and influenza	Genentech, GlaxoSmithKline, Chiron,
vaccines -- StressGen	Takeda
Helicobacter pylori vaccine -- Acambis,	Herpes simplex vaccine -- Cantab
AstraZeneca/CSL, Chiron, Provalis	Pharmaceuticals, CEL-SCI, Henderson
Helistat-G -- GalaGen	Morley
Hemolink -- Hemosol	Herpes simplex vaccine live -- ImClone
hepapoietin -- Snow Brand	Systems/Wyeth-Lederle, Aventis Pasteur
heparanase -- InSight	HGF derivatives -- Dompe
heparinase I -- Ibex	hIAPP vaccine -- Crucell
heparinase III -- Ibex	Hib-hepatitis B vaccine -- Aventis Pasteur
Hepatitis A vaccine -- American Biogenetic	HIC 1
Sciences	HIP -- Altachem
Hepatitis A vaccine inactivated	Hirudins -- Biopharma, Cangene, Dongkook,
Hepatitis A vaccine Nothav -- Chiron	Japan Energy Corporation, Pharmacia
Hepatitis A-hepatitis B vaccine --	Corporation, SIR International, Sanofi-
GlaxoSmithKline	Synthelabo, Sotragene, Rhein Biotech
hepatitis B therapy -- Tripep	HIV edible vaccine -- ProdiGene
Hepatitis B vaccine -- Amgen, Chiron SpA,	HIV gp120 vaccine -- Chiron, Ajinomoto,
Meiji Milk, NIS, Prodeva, PowderJect,	GlaxoSmithKline, ID Vaccine, Progenics,
Rhein Biotech	VaxGen
Hepatitis B vaccine recombinant -- Evans	HIV gp120 vaccine gene therapy --
Vaccines, Epitex Combiotech, Genentech,	HIV gp160 DNA vaccine -- PowderJect,
MedImmune, Merck Sharp & Dohme,	Aventis Pasteur, Oncogen, Hyland
Rhein Biotech, Shantha Biotechnics,	Immuno, Protein Sciences
Vector, Yeda	HIV gp41 vaccine -- Panacos
Hepatitis B vaccine recombinant TGP 943 --	HIV HGP-30W vaccine -- CEL-SCI
Takeda	HIV immune globulin -- Abbott, Chiron
Hepatitis C vaccine -- Bavarian Nordic,	HIV peptides -- American Home Products
Chiron, Innogenetics Acambis,	HIV vaccine -- Applied bioTech., Axis
Hepatitis D vaccine -- Chiron Vaccines	Genetics, Biogen, Bristol-Myers Squibb,
Hepatitis E vaccine recombinant --	Genentech, Korea Green Cross, NIS,
Genelabs/GlaxoSmithKline, Novavax	Oncogen, Protein Sciences Corporation,
hepatocyte growth factor -- Panorama,	Terumo, Tonen Corporation, Wyeth-
Sosei	Ayerst, Wyeth-Lederle Vaccines-Malvern,
hepatocyte growth factor kringle fragments --	Advanced BioScience Laboratories,
- EntreMed	Bavarian Nordic, Bavarian Nordic/Statens
Her-2/Neu peptides -- Corixa	Serum Institute, GeneCure, Immune
	Response, Progenics, Therion Biologics,
	United Biomedical, Chiron

FIG. 28M

44/497

HIV vaccine vCP1433 -- Aventis Pasteur	Human monoclonal antibodies --
HIV vaccine vCP1452 -- Aventis Pasteur	Medarex/Northwest Biotherapeutics,
HIV vaccine vCP205 -- Aventis Pasteur	Medarex/Seattle Genetics
HL-9 -- American BioScience	human netrin-1 -- Exelixis
HM-9239 -- Cytran	human papillomavirus antibodies -- Epicyte
HML-103 -- Hemosol	Human papillomavirus vaccine -- Biotech
HML-104 -- Hemosol	Australia, IDEC, StressGen
HML-105 -- Hemosol	Human papillomavirus vaccine MEDI 501 --
HML-109 -- Hemosol	MedImmune/GlaxoSmithKline
HML-110 -- Hemosol	Human papillomavirus vaccine MEDI
HML-121 -- Hemosol	503/MEDI 504 --
hNLP -- Pharis	MedImmune/GlaxoSmithKline
Hookworm vaccine	Human papillomavirus vaccine TA-CIN --
host-vector vaccines -- Henogen	Cantab Pharmaceuticals
HPM 1 -- Chugai	Human papillomavirus vaccine TA-HPV --
HPV vaccine -- MediGene	Cantab Pharmaceuticals
HSA -- Meristem	Human papillomavirus vaccine TH-GW --
HSF -- StressGen	Cantab/GlaxoSmithKline
HSP carriers -- Weizmann, Yeda, Peptor	human polyclonal antibodies -- Biosite/Eos
HSPPC-70 -- Antigenics	BioTech./ Medarex
HSPPC-96, pathogen-derived -- Antigenics	human type II anti factor VIII monoclonal
HSV 863 -- Novartis	antibodies -- ThromboGenics
HTLV-I DNA vaccine	humanised anti glycoprotein Ib murine
HTLV-I vaccine	monoclonal antibodies -- ThromboGenics
HTLV-II vaccine -- Access	HumaRAD -- Intracell
HU 901 -- Tanox	HuMax EGFR -- Genmab
Hu23F2G -- ICOS	HuMax-CD4 -- Medarex
HuHMFG1	HuMax-IL15 -- Genmab
HumaLYM -- Intracell	HYB 190 -- Hybridon
Human krebs statika -- Yamanouchi	HYB 676 -- Hybridon
human monoclonal antibodies --	I-125 MAb A33 -- Celltech
Abgenix/Biogen, Abgenix/ Corixa,	Ibritumomab tiuxetan -- IDEC
Abgenix/Immunex, Abgenix/Lexicon,	IBT-9401 -- Ibex
Abgenix/ Pfizer, Athersys/Medarex,	IBT-9402 -- Ibex
Biogen/MorphoSys, CAT/Searle,	IC 14 -- ICOS
Centocor/Medarex, Corixa/Kirin Brewery,	Idarubicin anti-Ly-2.1 --
Corixa/Medarex, Eos BioTech./Medarex,	IDEC 114 -- IDEC
Eos/Xenerex, Exelixis/Protein Design	IDEC 131 -- IDEC
Labs, ImmunoGen/ Raven, Medarex/	IDEC 152 -- IDEC
B.Twelve, MorphoSys/ImmunoGen, XTL	IDM 1 -- IDM
Biopharmaceuticals/Dyax,	IDPS -- Hollis-Eden Pharmaceuticals

FIG. 28N

45/497

iduronate-2-sulfatase -- Transkaryotic Therapies	insulin -- AutoImmune, Altea, Biobras, BioSante, Bio-Tech. General, Chong Kun Dang, Emisphere, Flamel, Provalis, Rhein Biotech, TranXenoGen
IGF/IBP-2-13 -- Pharis	insulin (bovine) -- Novartis
IGN-101 -- Igeneon	insulin analogue -- Eli Lilly
IK HIR02 -- Iketon	Insulin Aspart -- Novo Nordisk
IL-11 -- Genetics Institute/AHP	insulin detemir -- Novo Nordisk
IL-13-PE38 -- NeoPharm	insulin glargine -- Aventis
IL-17 receptor -- Immunex	insulin inhaled -- Inhale Therapeutics Systems, Alkermes
IL-18BP -- Yeda	insulin oral -- Inovax
IL-1Hy1 -- Hyseq	insulin, AeroDose -- AeroGen
IL-1 β -- Celltech	insulin, AERx -- Aradigm
IL-1 β adjuvant -- Celltech	insulin, BEODAS -- Elan
IL-2 -- Chiron	insulin, Biphasix -- Helix
IL-2 + IL-12 -- Hoffman La-Roche	insulin, buccal -- Generex
IL-6/sIL-6R fusion -- Hadasit	insulin, I2R -- Flemington
IL-6R derivative -- Tosoh	insulin, intranasal -- Bentley
IL-7-Dap 389 fusion toxin -- Ligand	insulin, oral -- Nobex, Unigene
IM-862 -- Cytran	insulin, Orasome -- Endorex
IMC-1C11 -- ImClone	insulin, ProMaxx -- Epic
imiglucerase -- Genzyme	insulin, Quadrant -- Elan
Immune globulin intravenous (human) -- Hoffman La Roche	insulin, recombinant -- Aventis
immune privilege factor -- Proneuron	insulin, Spiros -- Elan
Immunocal -- Immunotec	insulin, Transfersome -- IDEA
Immunogene therapy -- Briana Bio-Tech	insulin, Zymo, recombinant -- Novo Nordisk
Immunoliposomal 5-fluorodeoxyuridine-dipalmitate --	insulinotropin -- Scios
immunosuppressant vaccine -- Aixdie	Insulysin gene therapy --
immunotoxin -- Antisoma, NIH	integrin antagonists -- Merck
ImmuRAIT-Re-188 -- Immunomedics	interferon (Alpha2) -- SRC VB VECTOR, Viragen, Dong-A, Hoffman La-Roche, Genentech
imreg-1 -- Imreg	interferon -- BioMedicines, Human Genome Sciences
infertility -- Johnson & Johnson, E-TRANS	interferon (Alfa-n3) -- Interferon Sciences Intl.
Infliximab -- Centocor	interferon (Alpha), Biphasix -- Helix
Influenza virus vaccine -- Aventis Pasteur, Protein Sciences	
inhibin -- Biotech Australia, Human Therapeutics	
Inhibitory G protein gene therapy	
INKP-2001 -- InKine	
Inolimomab -- Diaclone	

FIG. 280

46/497

interferon (Alpha)—Amgen, BioNative,	IL-2/ diphtheria toxin — Ligand
Novartis, Genzyme Transgenics,	Interleukin-3 — Cangene
Hayashibara, Inhale Therapeutics	Interleukin-4 — Immunology Ventures,
Systems, Medusa, Flamel, Dong-A,	Sanofi Winthrop, Schering-Plough,
GeneTrol, Nastech, Shantha,	Immunex/ Sanofi Winthrop, Bayer, Ono
Wassermann, LG Chem, Sumitomo,	interleukin-4 + TNF-Alpha — NIH
Aventis, Behring EGIS, Pepgen, Servier,	interleukin-4 agonist — Bayer
Rhein Biotech,	interleukin-4 fusion toxin — Ligand
interferon (Alpha2A)	Interleukin-4 receptor — Immunex, Immun
interferon (Alpha2B) — Enzon, Schering-	Interleukin-6 — Ajinomoto, Cangene, Yeda,
Plough, Biogen, IDEA	Genetics Institute, Novartis
interferon (Alpha-N1) — GlaxoSmithKline	interleukin-6 fusion protein
interferon (beta) — Rentschler, GeneTrol,	Interleukin-6 fusion toxin — Ligand, Serono
Meristem, Rhein Biotech, Toray, Yeda,	interleukin-7 — IC Innovations
Daiichi, Mochida	interleukin-7 receptor — Immunex
interferon (Beta1A) — Serono, Biogen	interleukin-8 antagonists — Kyowa
interferon (beta1A), inhale — Biogen	Hakko/Millennium/Pfizer
interferon (B1b)— Chiron	interleukin-9 antagonists — Genaera
interferon (tau)— Pepgen	Interleukin-10 — DNAX, Schering-Plough
Interferon alfacon-1 — Amgen	Interleukin-10 gene therapy —
Interferon alpha-2a vaccine	interleukin-12 — Genetics Institute, Hoffman
Interferon Beta 1b — Schering/Chiron,	La-Roche
InterMune	interleukin-13 — Sanofi
Interferon Gamma — Boehringer Ingelheim,	interleukin-13 antagonists — AMRAD
Sheffield, Rentschler, Hayashibara	Interleukin-13-PE38QQR
interferon receptor, Type I — Serono	interleukin-15 — Immunex
interferon(Gamma1B) — Genentech	interleukin-16 — Research Corp
Interferon-alpha-2b + ribavirin — Biogen,	interleukin-18 — GlaxoSmithKline
ICN	Interleukin-18 binding protein — Serono
Interferon-alpha-2b gene therapy —	lor-P3 — Center of Molecular Immunology
Schering-Plough	IP-10 — NIH
Interferon-con1 gene therapy —	IPF — Metabolex
interleukin-1 antagonists — Dompe	IR-501 — Immune Response
Interleukin-1 receptor antagonist — Abbott	ISIS 9125 — Isis Pharmaceuticals
Bioresearch, Pharmacia	ISURF No. 1554 — Millennium
Interleukin-1 receptor type I — Immunex	ISURF No. 1866 — Iowa State Univer.
interleukin-1 receptor Type II — Immunex	ITF-1697 — Italfarmaco
Interleukin-1 trap — Regeneron	IxC 162 — Ixion
Interleukin-1-alpha — Immunex/Roche	J 695 — Cambridge Antibody Tech.,
interleukin-2 — SRC VB VECTOR,	Genetics Inst., Knoll
Ajinomoto, Biomira, Chiron	Jagged + FGF — Repair

FIG. 28P

47/497

JKC-362 -- Phoenix Pharmaceuticals	leptin, 2nd-generation -- Amgen
JTP-2942 -- Japan Tobacco	leridistim -- Pharmacia
Juman monoclonal antibodies -- Medarex/Raven	leuprolide, ProMaxx -- Epic
K02 -- Axy's Pharmaceuticals	leuprorelin, oral -- Unigene
Keliximab -- IDEC	LeuTech -- Papatin
Keyhole limpet haemocyanin	LEX 032 -- SuperGen
KGF -- Amgen	LiDEPT -- Novartis
KM 871 -- Kyowa	Lintuzumab (anti-CD33 MAb) -- Protein Design Labs
KPI 135 -- Scios	lipase -- Altus Biologics
KPI-022 -- Scios	lipid A vaccine -- EntreMed
Kringle 5	lipid-linked anchor Tech. -- ICRT, ID Biomedical
KSB 304	liposome-CD4 Tech. -- Sheffield
KSB-201 -- KS Biomedix	Listeria monocytogenes vaccine
L 696418 -- Merck	LMB 1
L 703801 -- Merck	LMB 7
L1 -- Acorda	LMB 9 -- Battelle Memorial Institute, NIH
L-761191 -- Merck	LM-CD45 -- Cantab Pharmaceuticals
lactoferrin -- Meristem, Pharming, Agennix	lovastatin -- Merck
lactoferrin cardio -- Pharming	LSA-3
LAG-3 -- Serono	LT- β receptor -- Biogen
LAIT -- GEMMA	lung cancer vaccine -- Corixa
LAK cell cytotoxin -- Arizona	lusupultide -- Scios
lamellarins -- PharmaMar/University of Malaga	L-Vax -- AVAX
laminin A peptides -- NIH	LY 355455 -- Eli Lilly
lanotepase -- Genetics Institute	LY 366405 -- Eli Lilly
laronidase -- BioMarin	LY-355101 -- Eli Lilly
Lassa fever vaccine	Lyme disease DNA vaccine -- Vical/Aventis Pasteur
LCAT -- NIH	Lyme disease vaccine -- Aquila
LDP 01 -- Millennium	Biopharmaceuticals, Aventis, Pasteur, Symbicom, GlaxoSmithKline, Hyland
LDP 02 -- Millennium	Immuno, MedImmune
Lecithinized superoxide dismutase -- Seikagaku	Lymphocytic choriomeningitis virus vaccine
LeIF adjuvant -- Corixa	lymphoma vaccine -- Biomira, Genitope
leishmaniasis vaccine -- Corixa	LYP18
lenercept -- Hoffman La-Roche	lys plasminogen, recombinant
Lenograstim -- Aventis, Chugai	Lysosomal storage disease gene therapy -- Avigen
lepirudin -- Aventis	lysostaphin -- Nutrition 21
leptin -- Amgen, IC Innovations	
Leptin gene therapy -- Chiron Corporation	

FIG. 28Q

48/497

M 23 -- Gruenenthal	MEDI 507 -- BioTransplant
M1 monoclonal antibodies -- Acorda	melanin concentrating hormone --
Therapeutics	Neurocrine Biosciences
MA 16N7C2 -- Corvas Intl.	melanocortins -- OMRF
malaria vaccine -- GlaxoSmithKline,	Melanoma monoclonal antibodies -- Viragen
AdProTech, Antigenics, Apovia, Aventis	melanoma vaccine -- GlaxoSmithKline,
Pasteur, Axis Genetics, Behringwerke,	Akzo Nobel, Avant, Aventis Pasteur,
CDCP, Chiron Vaccines, Genzyme	Bavarian Nordic, Biovector, CancerVax,
Transgenics, Hawaii, MedImmune, NIH,	Genzyme Molecular Oncology, Humbolt,
NYU, Oxon, Roche/Saramane, Biotech	ImClone Systems, Memorial, NYU, Oxon
Australia, Rx Tech	Melanoma vaccine Magevac -- Therion
Malaria vaccine CDC/NIIMALVAC-1	memory enhancers -- Scios
malaria vaccine, multicomponent	meningococcal B vaccine -- Chiron
mammaglobin -- Corixa	meningococcal vaccine -- CAMR
mammastatin -- Biotherapeutics	Meningococcal vaccine group B conjugate -
mannan-binding lectin -- NatImmu	- North American Vaccine
mannan-MUC1 -- Psiron	Meningococcal vaccine group B
MAP 30	recombinant -- BioChem Vaccines,
Marinovir -- Phytera	Microscience
MARstem -- Maret	Meningococcal vaccine group Y conjugate -
MB-015 -- Mochida	- North American Vaccine
MBP -- ImmuLogic	Meningococcal vaccine groups A B and C
MCI-028 -- Mitsubishi-Tokyo	conjugate -- North American Vaccine
MCIF -- Human Genome Sciences	Mepolizumab -- GlaxoSmithKline
MDC -- Advanced BioScience -- Akzo	Metastatin -- EntreMed, Takeda
Nobel, ICOS	Met-CkB7 -- Human Genome Sciences
MDX 11 -- Medarex	met-enkephalin -- TNI
MDX 210 -- Medarex	METH-1 -- Human Genome Sciences
MDX 22 -- Medarex	methioninase -- AntiCancer
MDX 22	Methionine lyase gene therapy --
MDX 240 -- Medarex	AntiCancer
MDX 33	Met-RANTES -- Genexa Biomedical,
MDX 44 -- Medarex	Serono
MDX 447 -- Medarex	Metreleptin
MDX H210 -- Medarex	Microtubule inhibitor MAb
MDX RA -- Houston BioTech., Medarex	Immunogen/Abgenix
ME-104 -- Pharmexa	MGDF -- Kirin
Measles vaccine	MGV -- Progenics
Mecasemin -- Cephalon/Chiron, Chiron	micrin -- Endocrine
MEDI 488 -- MedImmune	microplasmin -- ThromboGenics
MEDI 500	MIF -- Genetics Institute

FIG. 28R

49/497

migration inhibitory factor -- NIH	MAB 45-2D9- -- haematoporphyrin conjugate
Mim CD4.1 -- Xycte Therapies	MAB 4B4
mirostipen -- Human Genome Sciences	MAB 4E3-CPA conjugate -- BCM Oncologia
Mitumomab (BEC-2) -- ImClone Systems, Merck KGaA	MAB 4E3-daunorubicin conjugate
MK 852 -- Merck	MAB 50-6
MLN 1202 (Anti-CCR2 monoclonal antibody) -- Millenium Pharmaceuticals	MAB 50-61A -- Institut Pasteur
Mobenakin -- NIS	MAB 5A8 -- Biogen
molgramostim -- Genetics Institute, Novartis	MAB 791T/36-methotrexate conjugate
monoclonal antibodies -- Abgenix/Celltech, Immusol/ Medarex, Viragen/ Roslin Institute, Cambridge Antibody Tech./Elan	MAB 7c11.e8
MAB 108 --	MAB 7E11 C5-selenocystamine conjugate
MAB 10D5 --	MAB 93KA9 -- Novartis
MAB 14.18-interleukin-2 immunocytokine -- Lexigen	MAB A5B7-cisplatin conjugate -- Biodynamics Research, Pharmacia
MAB 14G2a --	MAB A5B7-I-131
MAB 15A10 --	MAB A7
MAB 170 -- Biomira	MAB A717 -- Exocell
MAB 177Lu CC49 --	MAB A7-zinostatin conjugate
MAB 17F9	MAB ABX-RB2 -- Abgenix
MAB 1D7	MAB ACA 11
MAB 1F7 -- Immune Network	MAB AFP-I-131 -- Immunomedics
MAB 1H10-doxorubicin conjugate	MAB AP1
MAB 26-2F	MAB AZ1
MAB 2A11	MAB B3-LysPE40 conjugate
MAB 2E1 -- RW Johnson	MAB B4 -- United Biomedical
MAB 2F5	MAB B43 Genistein-conjugate
MAB 31.1 -- International BioImmune Systems	MAB B43.13-Tc-99m -- Biomira
MAB 32 -- Cambridge Antibody Tech., Peptech	MAB B43-PAP conjugate
MAB 323A3 -- Centocor	MAB B4G7-gelonin conjugate
MAB 3C5	MAB BCM 43-daunorubicin conjugate -- BCM Oncologia
MAB 3F12	MAB BIS-1
MAB 3F8	MAB BMS 181170 -- Bristol-Myers Squibb
MAB 42/6	MAB BR55-2
MAB 425 -- Merck KGaA	MAB BW494
MAB 447-52D -- Merck Sharp & Dohme	MAB C 242-DM1 conjugate -- ImmunoGen
	MAB C242-PE conjugate
	MAB c30-6
	MAB CA208-cytorhodin-S conjugate -- Hoechst Japan
	MAB CC49 -- Enzon

FIG. 28S

50/497

MAb ch14.18 --	MAb LL2-I-131 -- Immunomedics
MAb CH14.18-GM-CSF fusion protein --	MAb LL2-Y-90
Lexigen	MAb LS2D617 -- Hybritech
MAb chCE7	MAb LYM-1-gelonin conjugate
MAb CI-137 -- AMRAD	MAb LYM-1-I-131
MAb cisplatin conjugate	MAb LYM-1-Y-90
MAb CLB-CD19	MAb LYM-2 -- Peregrine
MAb CLB-CD19v	MAb M195
MAb CLL-1 -- Peregrine	MAb M195-bismuth 213 conjugate --
MAb CLL-1-GM-CSF conjugate	Protein Design Labs
MAb CLL-1-IL-2 conjugate -- Peregrine	MAb M195-gelonin conjugate
MAb CLN IgG -- doxorubicin conjugates	MAb M195-I-131
MAb conjugates -- Tanox	MAb M195-Y-90
MAb D612	MAb MA 33H1 -- Sanofi
MAb Dal B02	MAb MAD11
MAb DC101 -- ImClone	MAb MGB2
MAb EA 1 --	MAb MINT5
MAb EC708 -- Biovation	MAb MK2-23
MAb EP-5C7 -- Protein Design Labs	MAb MOC31 ETA(252-613) conjugate
MAb ERIC-1 -- ICRT	MAb MOC-31-In-111
MAb F105 gene therapy	MAb MOC-31-PE conjugate
MAb FC 2.15	MAb MR6 --
MAb G250 -- Centocor	MAb MRK-16 -- Aventis Pasteur
MAb GA6	MAb MS11G6
MAb GA733	MAb MX-DTPA BrE-3
MAb Gliomab-H -- Viventia Biotech	MAb MY9
MAb HB2-saporin conjugate	MAb Nd2 -- Tosoh
MAb HD 37 --	MAb NG-1 -- Hygeia
MAb HD37-ricin chain-A conjugate	MAb NM01 -- Nissin Food
MAb HNK20 -- Acambis	MAb OC 125
MAb huN901-DM1 conjugate --	MAb OC 125-CMA conjugate
ImmunoGen	MAb OKI-1 -- Ortho-McNeil
MAb I-131 CC49 -- Corixa	MAb OX52 -- Bioproducts for Science
MAb ICO25	MAb PMA5
MAb ICR12-CPG2 conjugate	MAb PR1
MAb ICR-62	MAb prost 30
MAb IRac-ricin A conjugate	MAb R-24
MAb K1	MAb R-24 α Human GD3 -- Celltech
MAb KS1-4-methotrexate conjugate	MAb RFB4-ricin chain A conjugate
MAb L6 -- Bristol-Myers Squibb, Oncogen	MAb RFT5-ricin chain A conjugate
MAb LiCO 16-88	MAb SC 1

FIG. 28T

51/497

MAB SM-3 -- ICRT	Muc-1 vaccine -- Corixa
MAB SMART 1D10 -- Protein Design Labs	mucosal tolerance -- Aberdeen
MAB SMART ABL 364 -- Novartis	mullerian inhibiting subst
MAB SN6f	muplestim -- Genetics Institute, Novartis,
MAB SN6f-deglycosylated ricin A chain	DSM Anti-Infectives
conjugate --	murine MAB -- KS Biomedix
MAB SN6j	Mutant somatropin -- JCR Pharmaceutical
MAB SN7-ricin chain A conjugate	MV 833 -- Toagosei
MAB T101-Y-90 conjugate -- Hybritech	Mycoplasma pulmonis vaccine
MAB T-88 -- Chiron	Mycoprex -- XOMA
MAB TB94 -- Cancer ImmunoBiology	myeloperoxidase -- Henogen
MAB TEC 11	myostatin -- Genetics Institute
MAB TES-23 -- Chugai	Nacolomab tafenatox -- Pharmacia
MAB TM31 -- Avant	Nagrecor -- Scios
MAB TNT-1 -- Cambridge Antibody Tech.,	nagrestipen -- British Biotech
Peregrine	NAP-5 -- Corvas Intl.
MAB TNT-3	NAPc2 -- Corvas Intl.
MAB TNT-3 -- IL2 fusion protein --	nartograstim -- Kyowa
MAB TP3-At-211	Natalizumab -- Protein Design Labs
MAB TP3-PAP conjugate --	Nateplase -- NIH, Nihon Schering
MAB UJ13A -- ICRT	nateplase -- Schering AG
MAB UN3	NBI-3001 -- Neurocrine Biosci.
MAB ZME-018-gelonin conjugate	NBI-5788 -- Neurocrine Biosci.
MAB-BC2 -- GlaxoSmithKline	NBI-6024 -- Neurocrine Biosci.
MAB-DM1 conjugate -- ImmunoGen	Nef inhibitors -- BRI
MAB-ricin-chain-A conjugate -- XOMA	Neisseria gonorrhoea vaccine -- Antex
MAB-temoporfin conjugates	Biologics
Monopharm C -- Viventia Biotech	Neomycin B-arginine conjugate
monteplase -- Eisai	Nerelimomab -- Chiron
montirelin hydrate -- Gruenenthal	Nerve growth factor -- Amgen -- Chiron,
moroctocog alfa -- Genetics Institute	Genentech
Moroctocog-alfa -- Pharmacia	Nerve growth factor gene therapy
MP 4	nesiritide citrate -- Scios
MP-121 -- Biopharm	neuregulin-2 -- CeNeS
MP-52 -- Biopharm	neurocan -- NYU
MRA -- Chugai	neuronal delivery system -- CAMR
MS 28168 -- Mitsui Chemicals, Nihon	Neurophil inhibitory Factor -- Corvas
Schering	Neuroprotective vaccine -- University of
MSH fusion toxin -- Ligand	Auckland
MSI-99 -- Genaera	neurotrophic chimaeras -- Regeneron
MT 201 -- Micromet	neurotrophic factor -- NsGene, CereMedix

FIG. 28U

52/497

NeuroVax -- Immune Response	Oncophage -- Antigenics
neurturin -- Genentech	Oncostatin M -- Bristol-Myers Squibb
neutral endopeptidase -- Genentech	OncoVax-CL -- Jenner Biotherapies
NGF enhancers -- NeuroSearch	OncoVax-P -- Jenner Biotherapies
NHL vaccine -- Large Scale Biology	oncept -- Yeda
NIP45 -- Boston Life Sciences	onychomycosis vaccine -- Boehringer
NKI-B20	Ingelheim
NM 01 -- Nissin Food	opebecan -- XOMA
NMI-139 -- NitroMed	opioids -- Arizona
NMMP -- Genetics Institute	Oprelvekin -- Genetics Institute
NN-2211 -- Novo Nordisk	Oregovomab -- AltaRex
Noggin -- Regeneron	Org-33408 b -- Akzo Nobel
Nonacog alfa	Orolip DP -- EpiCept
Norelin -- Biostar	oryzacystatin
Norwalk virus vaccine	OSA peptides -- GenSci Regeneration
NRLU 10 -- NeoRx	osteoblast-cadherin GF -- Pharis
NRLU 10 PE -- NeoRx	Osteocalcin-thymidine kinase gene therapy
NT-3 -- Regeneron	osteogenic protein -- Curis
NT-4/5 -- Genentech	osteopontin -- OraPharma
NU 3056	osteoporosis peptides -- Integra, Telios
NU 3076	osteoprotegerin -- Amgen, SnowBrand
NX 1838 -- Gilead Sciences	otitis media vaccines -- Antex Biologics
NY ESO-1/CAG-3 antigen -- NIH	ovarian cancer -- University of Alabama
NYVAC-7 -- Aventis Pasteur	OX40-IgG fusion protein -- Cantab, Xenova
NZ-1002 -- Novazyme	P 246 -- Diatide
obesity therapy -- Nobex	P 30 -- Alfacell
OC 10426 -- Ontogen	p1025 -- Active Biotech
OC 144093 -- Ontogen	P-113 ^A -- Demegen
OCIF -- Sankyo	P-16 peptide -- Transition Therapeutics
Oct-43 -- Otsuka	p43 -- Ramot
Odulimomab -- Immunotech	P-50 peptide -- Transition Therapeutics
OK PSA - liposomal	p53 + RAS vaccine -- NIH, NCI
OKT3-gamma-1-ala-ala	PACAP(1-27) analogue
OM 991	paediatric vaccines -- Chiron
OM 992	Pafase -- ICOS
Omalizumab -- Genentech	PAGE-4 plasmid DNA -- IDEC
oncoimmunin-L -- NIH	PAI-2 -- Biotech Australia, Human
Oncolysin B -- ImmunoGen	Therapeutics
Oncolysin CD6 -- ImmunoGen	Palifermin (keratinocyte growth factor) --
Oncolysin M -- ImmunoGen	Amgen
Oncolysin S -- ImmunoGen	Palivizumab -- MedImmune

FIG. 28V

53/497

PAM 4 -- Merck	PEG-uricase -- Mountain View
pamiteplase -- Yamanouchi	Pegvisomant -- Genentech
pancreatin, Minitabs -- Eurand	PEGylated proteins, PolyMASC -- Valentis
Pangen -- Fournier	PEGylated recombinant native human leptin
Pantarin -- Selective Genetics	-- Roche
Parainfluenza virus vaccine -- Pharmacia,	Pemtumomab
Pierre Fabre	Penetratin -- Cyclacel
paraoxanase -- Esperion	Pepscan -- Antisoma
parathyroid hormone -- Abiogen, Korea	peptide G -- Peptech, ICRT
Green Cross	peptide vaccine -- NIH, NCI
Parathyroid hormone (1-34) --	Pexelizumab
Chugai/Suntory	pexiganan acetate -- Genaera
Parkinson's disease gene therapy -- Cell	Pharmaprojects No. 3179 -- NYU
Genesys/ Ceregene	Pharmaprojects No. 3390 -- Ernest Orlando
Parvovirus vaccine -- MedImmune	Pharmaprojects No. 3417 -- Sumitomo
PCP-Scan -- Immunomedics	Pharmaprojects No. 3777 -- Acambis
PDGF -- Chiron	Pharmaprojects No. 4209 -- XOMA
PDGF cocktail -- Theratechnologies	Pharmaprojects No. 4349 -- Baxter Intl.
peanut allergy therapy -- Dynavax	Pharmaprojects No. 4651
PEG anti-ICAM MAb -- Boehringer	Pharmaprojects No. 4915 -- Avanir
Ingelheim	Pharmaprojects No. 5156 -- Rhizogenics
PEG asparaginase -- Enzon	Pharmaprojects No. 5200 -- Pfizer
PEG glucocerebrosidase	Pharmaprojects No. 5215 -- Origene
PEG hirudin -- Knoll	Pharmaprojects No. 5216 -- Origene
PEG interferon-alpha-2a -- Roche	Pharmaprojects No. 5218 -- Origene
PEG interferon-alpha-2b + ribavirin --	Pharmaprojects No. 5267 -- ML
Biogen, Enzon, ICN Pharmaceuticals,	Laboratories
Schering-Plough	Pharmaprojects No. 5373 -- MorphoSys
PEG MAb A5B7 --	Pharmaprojects No. 5493 -- Metabolex
Pegacaristim -- Amgen -- Kirin Brewery --	Pharmaprojects No. 5707 -- Genentech
ZymoGenetics	Pharmaprojects No. 5728 -- Autogen
Pegaldesleukin -- Research Corp	Pharmaprojects No. 5733 -- BioMarin
pegaspargase -- Enzon	Pharmaprojects No. 5757 -- NIH
pegfilgrastim -- Amgen	Pharmaprojects No. 5765 -- Gryphon
PEG-interferon Alpha -- Viragen	Pharmaprojects No. 5830 -- AntiCancer
PEG-interferon Alpha 2A -- Hoffman La-	Pharmaprojects No. 5839 -- Dyax
Roche	Pharmaprojects No. 5849 -- Johnson &
PEG-interferon Alpha 2B -- Schering-	Johnson
Plough	Pharmaprojects No. 5860 -- Mitsubishi-
PEG-r-hirudin -- Abbott	Tokyo
PEG-rHuMGDF -- Amgen	

FIG. 28W

54/497

Pharmaprojects No. 5869 -- Oxford GlycoSciences	Plasminogen activators -- Abbott Laboratories, American Home Products, Boehringer Mannheim, Chiron
Pharmaprojects No. 5883 -- Asahi Brewery	Corporation, DuPont Pharmaceuticals, Eli Lilly, Shionogi, Genentech, Genetics Institute, GlaxoSmithKline, Hemispherx
Pharmaprojects No. 5947 -- StressGen	Biopharma, Merck & Co, Novartis, Pharmacia Corporation, Wakamoto, Yeda
Pharmaprojects No. 5961 -- Theratechnologies	plasminogen-related peptides -- Bio-Tech. General/MGH
Pharmaprojects No. 5962 -- NIH	platelet factor 4 -- RepliGen
Pharmaprojects No. 5966 -- NIH	Platelet-derived growth factor -- Amgen -- ZymoGenetics
Pharmaprojects No. 5994 -- Pharming	plusonemin-- Hayashibara
Pharmaprojects No. 5995 -- Pharming	PMD-2850 -- Protherics
Pharmaprojects No. 6023 -- IMMUCON	Pneumococcal vaccine -- Antex Biologics, Aventis Pasteur
Pharmaprojects No. 6063 -- Cytoclonal	Pneumococcal vaccine intranasal -- BioChem Vaccines/Biovector
Pharmaprojects No. 6073 -- SIDDCO	PR1A3
Pharmaprojects No. 6115 -- Genzyme	PR-39
Pharmaprojects No. 6227 -- NIH	pralmorelin -- Kaken
Pharmaprojects No. 6230 -- NIH	Pretarget-Lymphoma -- NeoRx
Pharmaprojects No. 6236 -- NIH	Priliximab -- Centocor
Pharmaprojects No. 6243 -- NIH	PRO 140 -- Progenics
Pharmaprojects No. 6244 -- NIH	PRO 2000 -- Procept
Pharmaprojects No. 6281 -- Senetek	PRO 367 -- Progenics
Pharmaprojects No. 6365 -- NIH	PRO 542 -- Progenics
Pharmaprojects No. 6368 -- NIH	pro-Apo A-I -- Esperion
Pharmaprojects No. 6373 -- NIH	prolactin -- Genzyme
Pharmaprojects No. 6408 -- Pan Pacific	Prosaptide TX14(A) -- Bio-Tech. General
Pharmaprojects No. 6410 -- Athersys	prostate cancer antibodies -- Immunex, UroCor
Pharmaprojects No. 6421 -- Oxford GlycoSciences	prostate cancer antibody therapy -- Genentech/UroGenesys, Genotherapeutics
Pharmaprojects No. 6522 -- Maxygen	prostate cancer immunotherapeutics -- The PSMA Development Company
Pharmaprojects No. 6523 -- Pharis	prostate cancer vaccine -- Aventis Pasteur, Zonagen, Corixa, Dendreon, Jenner
Pharmaprojects No. 6538 -- Maxygen	Biotherapies, Therion Biologics
Pharmaprojects No. 6554 -- APALEXO	
Pharmaprojects No. 6560 -- Ardana	
Pharmaprojects No. 6562 -- Bayer	
Pharmaprojects No. 6569 -- Eos	
Phenoxazine	
Phenylase -- Ibex	
Pigment epithelium derived factor -- plasminogen activator inhibitor-1, recombinant -- DuPont Pharmaceuticals	

FIG. 28X

55/497

prostate-specific antigen -- Entremed	RD 62198
protein A -- Repligen	rDnase -- Genentech
protein adhesives -- Enzon	RDP-58 -- SangStat
protein C -- Baxter Intl., PPL Therapeutics,	RecepTox-Fce -- Keryx
ZymoGenetics	RecepTox-GnRH -- Keryx, MTR
protein C activator -- Gilead Sciences	Technologies
protein kinase R antags -- NIH	RecepTox-MBP -- Keryx, MTR
protirelin -- Takeda	Technologies
protocadherin 2 -- Caprion	recFSH -- Akzo Nobel, Organon
Pro-urokinase -- Abbott, Bristol-Myers	REGA 3G12
Squibb, Dainippon, Tosoh -- Welfide	Regavirumab -- Teijin
P-selectin glycoprotein ligand-1 -- Genetics	relaxin -- Connetics Corp
Institute	Renal cancer vaccine -- Macropharm
pseudomonal infections -- InterMune	repifermin -- Human Genome Sciences
Pseudomonas vaccine -- Cytovax	Respiratory syncytial virus PFP-2 vaccine --
PSGL-Ig -- American Home Products	Wyeth-Lederle
PSP-94 -- Procyon	Respiratory syncytial virus vaccine --
PTH 1-34 -- Nobex	GlaxoSmithKline, Pharmacia, Pierre Fabre
Quilimmune-M -- Antigenics	Respiratory syncytial virus vaccine
R 744 -- Roche	inactivated
R 101933	Respiratory syncytial virus-parainfluenza
R 125224 -- Sankyo	virus vaccine -- Aventis Pasteur,
RA therapy -- Cardion	Pharmacia
Rabies vaccine recombinant -- Aventis	Reteplase -- Boehringer Mannheim,
Pasteur, BioChem Vaccines, Kaketsuken	Hoffman La-Roche
Pharmaceuticals	Retropep -- Retroscreen
RadioTheraCIM -- YM BioSciences	RFB4 (dsFv) PE38
Ramot project No. 1315 -- Ramot	RFI 641 -- American Home Products
Ramot project No. K-734A -- Ramot	RFTS -- UAB Research Foundation
Ramot project No. K-734B -- Ramot	RG 12986 -- Aventis Pasteur
Ranibizumab (Anti-VEGF fragment) --	RG 83852 -- Aventis Pasteur
Genentech	RG-1059 -- Repligen
RANK -- Immunex	rGCR -- NIH
ranpirnase -- Alfacell	rGLP-1 -- Restoragen
ranpirnase-anti-CD22 MAb -- Alfacell	rGRF -- Restoragen
RANTES inhibitor -- Milan	rh Insulin -- Eli Lilly
RAPID drug delivery systems -- ARIAD	RHAMM targeting peptides -- Cangene
rasburicase -- Sanofi	rHb1.1 -- Baxter Intl.
rBPI-21, topical -- XOMA	rhCC10 -- Claragen
RC 529 -- Corixa	rhCG -- Serono
rCFTR -- Genzyme Transgenics	Rheumatoid arthritis gene therapy

FIG. 28Y

56/497

Rheumatoid arthritis vaccine -- Veterans Affairs Medical Center	SB RA 31012 --
rhLH -- Serono	SC 56929 -- Pharmacia
Ribozyme gene therapy -- Genset	SCA binding proteins -- Curis, Enzon
Rickettsial vaccine recombinant	scFv(14E1)-ETA Berlex Laboratories, Schering AG
RIGScan CR -- Neoprobe	ScFv(FRP5)-ETA --
RIP-3 -- Rigel	ScFv6C6-PE40 --
Rituximab -- Genentech	SCH 55700 -- Celltech
RK-0202 -- RxKinetix	Schistosomiasis vaccine -- Glaxo Wellcome/Medeva, Brazil
RLT peptide -- Esperion	SCPF -- Advanced Tissue Sciences
rM/NEI -- IVAX	scuPA-suPAR complex -- Hadasit
rmCRP -- Immtech	SD-9427 -- Pharmacia
RN-1001 -- Renovo	SDF-1 -- Ono
RN-3 -- Renovo	SDZ 215918 -- Novartis
RNase conjugate -- Immunomedics	SDZ 280125 -- Novartis
RO 631908 -- Roche	SDZ 89104 -- Novartis
Rotavirus vaccine -- Merck	SDZ ABL 364 -- Novartis
RP 431 -- DuPont Pharmaceuticals	SDZ MMA 383 -- Novartis
RP-128 -- Resolution	Secretin -- Ferring, Repligen
RPE65 gene therapy --	serine protease inhbs -- Pharis
RPR 110173 -- Aventis Pasteur	sermorelin acetate -- Serono
RPR 115135 -- Aventis Pasteur	SERP-1 -- Viron
RPR 116258A -- Aventis Pasteur	sertenef -- Daiippon
rPSGL-Ig -- American Home Products	serum albumin, Recombinant human -- Aventis Behring
r-SPC surfactant -- Byk Gulden	serum-derived factor -- Hadasit
RSV antibody -- Medimmune	Sevirumab -- Novartis
Ruplizumab -- Biogen	SGN 14 -- Seattle Genetics
rV-HER-2/neu -- Therion Biologics	SGN 15 -- Seattle Genetics
SA 1042 -- Sankyo	SGN 17/19 -- Seattle Genetics
sacrosidase -- Orphan Medical	SGN 30 -- Seattle Genetics
Sant 7	SGN-10 -- Seattle Genetics
Sargramostim -- Immunex	SGN-11 -- Seattle Genetics
saruplase -- Gruenenthal	SH 306 -- DuPont Pharmaceuticals
Satumomab -- Cytogen	Shanvac-B -- Shantha
SB 1 -- COR Therapeutics	Shigella flexneri vaccine -- Avant, Acambis, Novavax
SB 207448 -- GlaxoSmithKline	Shigella sonnei vaccine --
SB 208651 -- GlaxoSmithKline	siCAM-1 -- Boehringer Ingelheim
SB 240683 -- GlaxoSmithKline	Silteplase -- Genzyme
SB 249415 -- GlaxoSmithKline	
SB 249417 -- GlaxoSmithKline	
SB 6 -- COR Therapeutics	

FIG. 28Z

57/497

SIV vaccine -- Endocon, Institut Pasteur	Staphylococcus aureus vaccine conjugate --
SK 896 -- Sanwa Kagaku Kenkyusho	Nabi
SK-827 -- Sanwa Kagaku Kenkyusho	Staphylococcus therapy -- Tripep
Skeletex -- CellFactors	Staphylokinase -- Biovation, Prothera,
SKF 106160 -- GlaxoSmithKline	Thrombogenetics
S-nitroso-AR545C --	Streptococcal A vaccine -- M6
SNTP -- Active Biotech	Pharmaceuticals, North American Vaccine
somatomedin-1 -- GroPep, Mitsubishi-	Streptococcal B vaccine -- Microscience
Tokyo, NIH	Streptococcal B vaccine recombinant --
somatomedin-1 carrier protein -- Insmmed	Biochem Vaccines
somatostatin -- Ferring	Streptococcus pyogenes vaccine
Somatotropin/	STRL-33 -- NIH
Human Growth Hormone -- Bio-Tech.	Subalin -- SRC VB VECTOR
General, Eli Lilly	SUIS -- United Biomedical
somatropin -- Bio-Tech. General, Alkermes,	SUIS-LHRH -- United Biomedical
ProLease, Aventis Behring, Biovector,	SUN-E3001 -- Suntory
Cangene, Dong-A, Eli Lilly, Emisphere,	super high affinity monoclonal antibodies --
Enact, Genentech, Genzyme Transgenics,	YM BioSciences
Grandis/InfiMed, CSL, InfiMed, MacroMed,	Superoxide dismutase -- Chiron, Enzon,
Novartis, Novo Nordisk, Pharmacia	Ube Industries, Bio-Tech, Yeda
Serono, TranXenoGen	superoxide dismutase-2 -- OXIS
somatropin derivative -- Schering AG	suppressin -- UAB Research Foundation
somatropin, AIR -- Eli Lilly	SY-161-P5 -- ThromboGenics
Somatropin, inhaled -- Eli Lilly/Alkermes	SY-162 -- ThromboGenics
somatropin, Kabi -- Pharmacia	Systemic lupus erythematosus vaccine --
somatropin, Orasome -- Novo Nordisk	MedClone/VivoRx
Sonermin -- Dainippon Pharmaceutical	T cell receptor peptides -- Xoma
SP(V5.2)C -- Supertek	T cell receptor peptide vaccine
SPf66	T4N5 liposomes -- AGI Dermatics
sphingomyelinase -- Genzyme	TACI, soluble -- ZymoGenetics
SR 29001 -- Sanofi	targeted apoptosis -- Antisoma
SR 41476 -- Sanofi	tasonermin -- Boehringer Ingelheim
SR-29001 -- Sanofi	TASP
SS1(dsFV)-PE38 -- NeoPharm	TASP-V
β 2 microglobulin -- Avidex	Tat peptide analogues -- NIH
β 2-microglobulin fusion proteins -- NIH	TBP I -- Yeda
β -amyloid peptides -- CeNeS	TBP II
β -defensin -- Pharis	TBV25H -- NIH
Staphylococcus aureus infections --	Tc 99m ior cea1 -- Center of Molecular
Inhibitex/ZLB	Immunology
	Tc 99m P 748 -- Diatide

FIG. 28AA

58/497

Tc 99m votumumab -- Intracell	Tissue factor -- Genentech
Tc-99m rh-Annexin V -- Theseus Imaging	Tissue factor pathway inhibitor
teceleukin -- Biogen	TJN-135 -- Tsumura
tenecteplase -- Genentech	TM 27 -- Avant
Teriparatide -- Armour Pharmaceuticals, Asahi Kasei, Eli Lilly	TM 29 -- Avant
terlipressin -- Ferring	TMC-151 -- Tanabe Seiyaku
testisin -- AMRAD	TNF tumour necrosis factor -- Asahi Kasei
Tetra fibrin -- Roche	TNF Alpha -- Cytimmune
TFPI -- EntreMed	TNF antibody -- Johnson & Johnson
tgD-IL-2 -- Takeda	TNF binding protein -- Amgen
TGF-Alpha -- ZymoGenetics	TNF degradation product -- Oncotech
TGF- β -- Kolon	TNF receptor -- Immunex
TGF- β 2 -- Insmed	TNF receptor 1, soluble -- Amgen
TGF- β 3 -- OSI	TNF Tumour necrosis factor-alpha -- Asahi Kasei, Genentech, Mochida
Thalassaemia gene therapy -- Crucell	TNF-Alpha inhibitor -- Tripep
TheraCIM-h-R3 -- Center of Molecular Immunology, YM BioSciences	TNFR:Fc gene therapy -- Targeted Genetics
Theradigm-HBV -- Epimmune	TNF-SAM2
Theradigm-HPV -- Epimmune	Tolerimab -- Innogenetics
Theradigm-malaria -- Epimmune	Toxoplasma gondii vaccine -- GlaxoSmithKline
Theradigm-melanoma -- Epimmune	TP 9201 -- Telios
TheraFab -- Antisoma	TP10 -- Avant
ThGRF 1-29 -- Theratechnologies	TP20 -- Avant
ThGRF 1-44 -- Theratechnologies	tPA -- Centocor
Thrombin receptor activating peptide -- Abbott	trafermin -- Scios
thrombomodulin -- Iowa, Novocastra	TRAIL/Apo2L -- Immunex
Thrombopoietin -- Dragon Pharmaceuticals, Genentech	TRAIL-R1 MAb -- Cambridge Antibody Technologies
thrombopoietin, Pliva -- Recepton	transferrin-binding proteins -- CAMR
Thrombospondin 2 --	Transforming growth factor-beta-1 -- Genentech
thrombostatin -- Thromgen	transport protein -- Genesis
thymalfasin -- SciClone	Trastuzumab -- Genentech
thymocartin -- Gedeon Richter	TRH -- Ferring
thymosin Alpha1 -- NIH	Triabin -- Schering AG
thyroid stimulating hormone -- Genzyme	Triconal
tICAM-1 -- Bayer	Triflavin
Tick anticoagulant peptide -- Merck	troponin I -- Boston Life Sciences
TIF -- Xoma	TRP-2 ^A -- NIH
Tifacogin -- Chiron, NIS, Pharmacia	trypsin inhibitor -- Mochida



FIG. 28BB

59/497

TSP-1 gene therapy –	Vascular endothelial growth factors – R&D
TT-232	Systems
TTS-CD2 – Active Biotech	vascular targeting agents – Peregrine
Tuberculosis vaccine – Aventis Pasteur,	vasopermeation enhancement agents --
Genesis	Peregrine
Tumor Targeted Superantigens – Active	vasostatin – NIH
Biotech – Pharmacia	VCL – Bio-Tech. General
tumour vaccines – PhotoCure	VEGF – Genentech, Scios
tumour-activated prodrug antibody	VEGF inhibitor – Chugai
conjugates – Millennium/ImmunoGen	VEGF-2 – Human Genome Sciences
tumstatin – ILEX	VEGF-Trap – Regeneron
Tuvirumab – Novartis	viscumin, recombinant – Madaus
TV-4710 – Teva	Vitaxin
TWEAK receptor -- Immunex	Vitrax -- ISTA Pharmaceuticals
TXU-PAP	West Nile virus vaccine -- Bavarian Nordic
TY-10721 – TOA Eiyo	WP 652
Type I diabetes vaccine -- Research Corp	WT1 vaccine -- Corixa
Typhoid vaccine CVD 908	WX-293 – Willex BioTech.
U 143677 -- Pharmacia	WX-360 -- Willex BioTech.
U 81749 -- Pharmacia	WX-UK1 – Willex BioTech.
UA 1248 – Arizona	XMP-500 – XOMA
UGIF – Sheffield	XomaZyme-791 – XOMA
UIC 2	XTL 001 – XTL Biopharmaceuticals
UK 101	XTL 002 – XTL Biopharmaceuticals
UK-279276 – Corvas Intl.	yeast delivery system -- GlobelImmune
urodilatin – Pharis	Yersinia pestis vaccine
urofollitrophin – Serono	YIGSR-Stealth – Johnson & Johnson
Urokinase – Abbott	Yisum Project No. D-0460 -- Yisum
uteroferrin-- Pepgen	YM 207 – Yamanouchi
V 20 – GLYCODesign	YM 337 -- Protein Design Labs
V2 vasopressin receptor gene therapy	Yttrium-90 labelled biotin
vaccines -- Active Biotech	Yttrium-90-labeled anti-CEA MAb T84.66 –
Varicella zoster glycoprotein vaccine --	ZD 0490 – AstraZeneca
Research Corporation Technologies	ziconotide – Elan
Varicella zoster virus vaccine live – Cantab	ZK 157138 -- Berlex Laboratories
Pharmaceuticals	Zolimomab aritox
Vascular endothelial growth factor –	Zorcell – Immune Response
Genentech, University of California	ZRXL peptides -- Novartis

FIG. 28CC

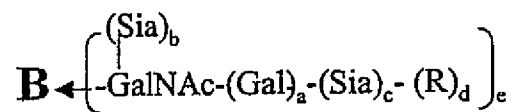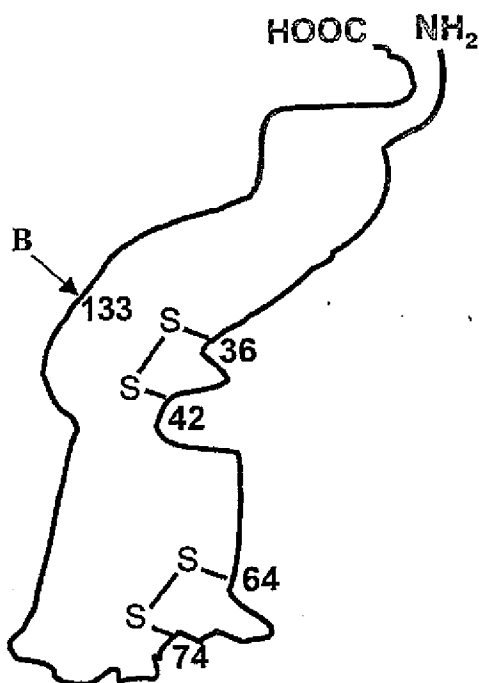60/497



a-c, e (independently selected) = 0 or 1;

d = 0;

R = modifying group, sialyl or
oligosialyl

FIG. 29A

61/497

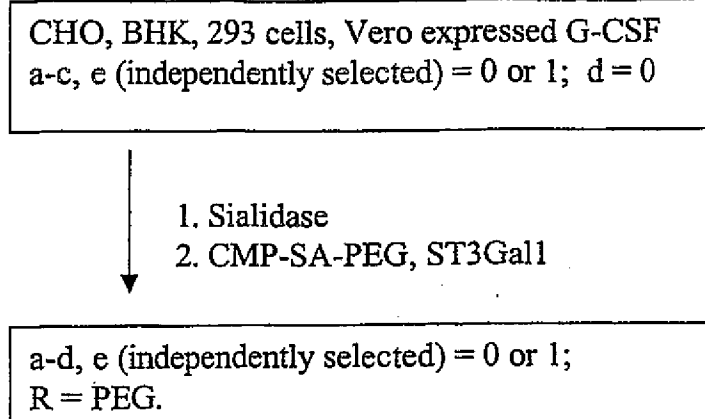


FIG. 29B

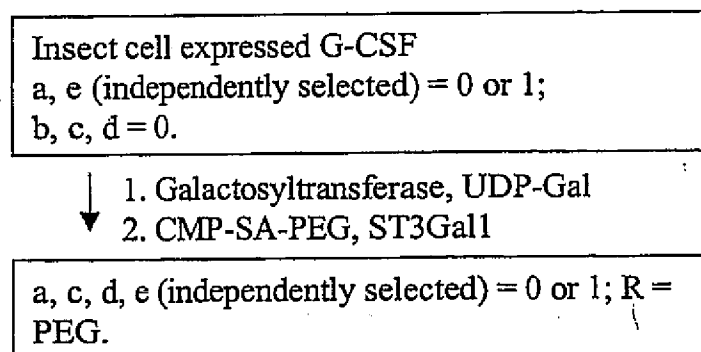


FIG. 29C

62/497

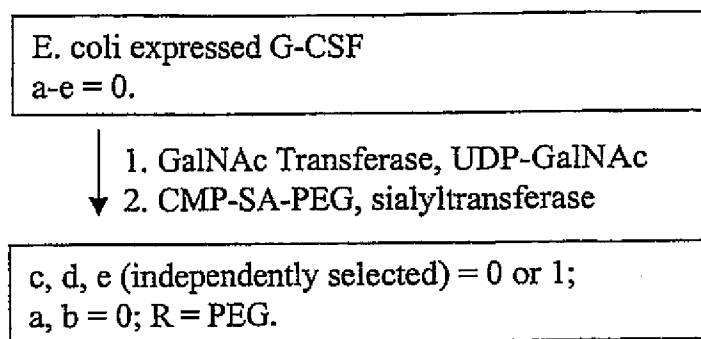


FIG. 29D

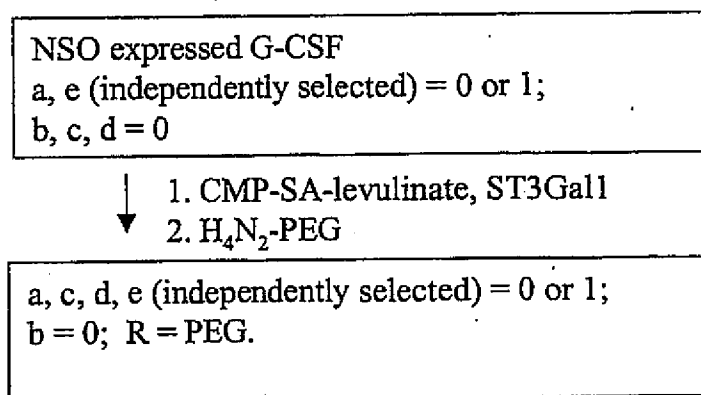


FIG. 29E

63/497

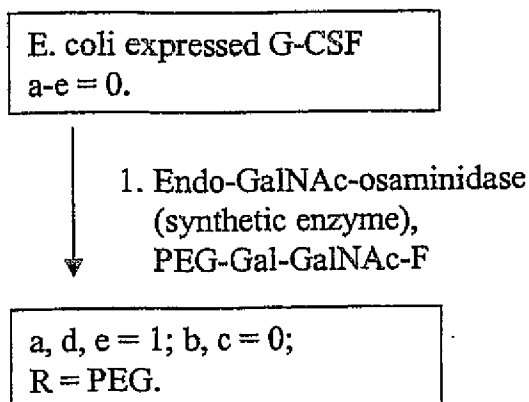


FIG. 29F


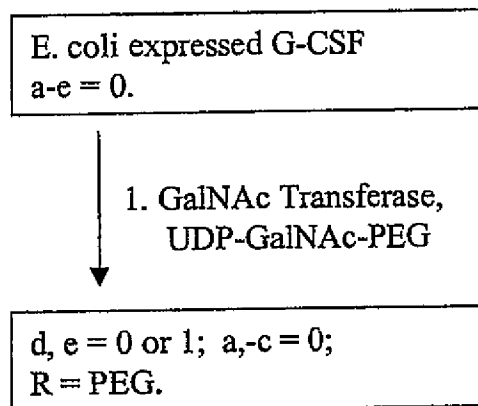


FIG. 29G

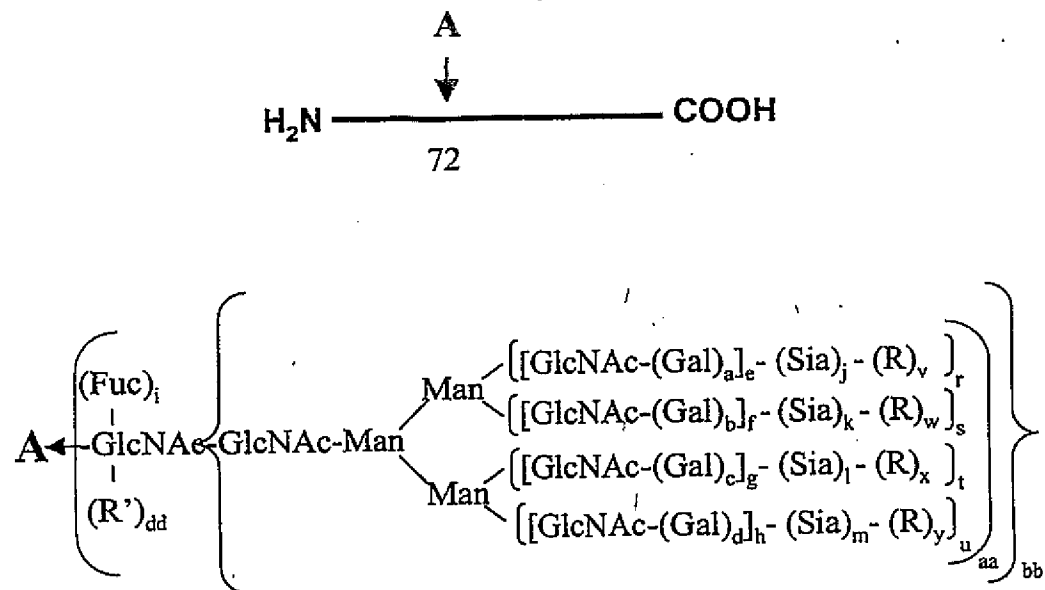64/497



a-d, i, n-u (independently selected) = 0 or 1.

aa, bb, cc, dd, ee (independently selected) = 0 or 1.

e-h (independently selected) = 0 to 6.

j-m (independently selected) = 0 to 20.

v-z = 0; R = modifying group, mannose, oligo-mannose.

R' = H, glycosyl residue, modifying group, glycoconjugate.

FIG. 30A

65/497

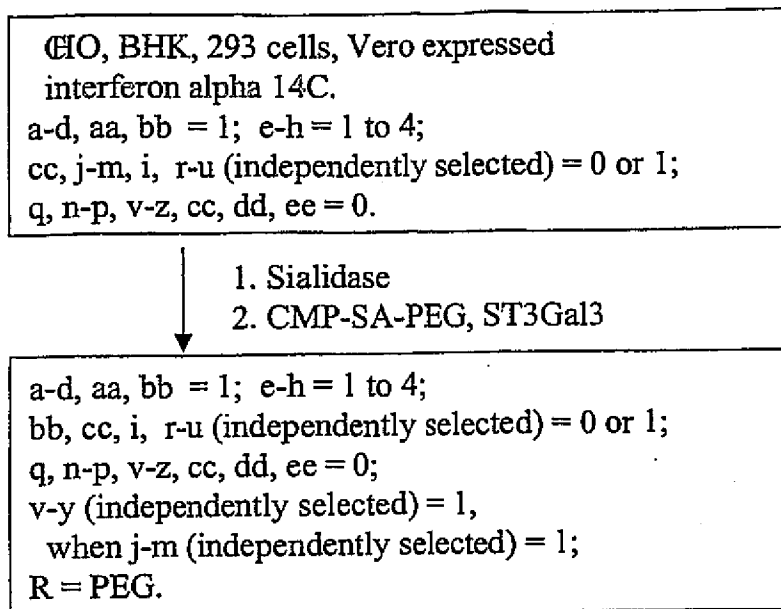


FIG. 30B

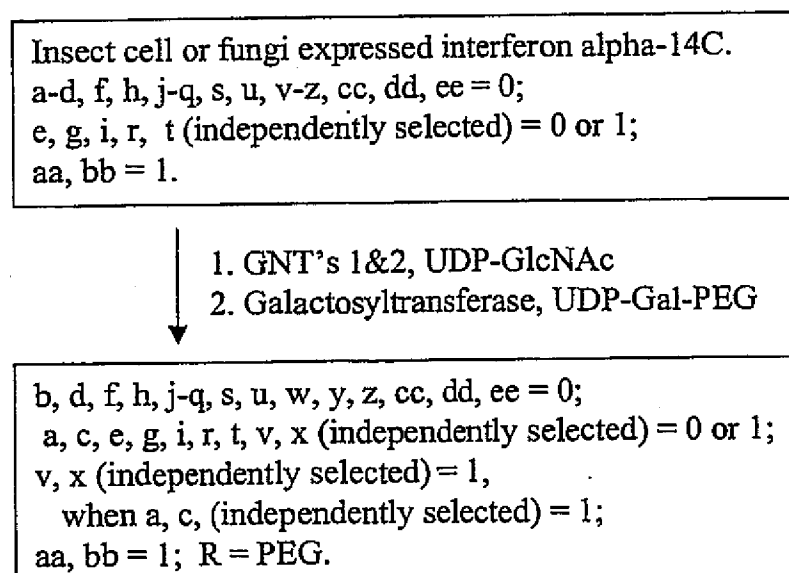


FIG. 30C

66/497


Yeast expressed interferon alpha-14C.
a-q, cc, dd, ee, v-z = 0;
r-y (independently selected) = 0 to 1;
aa, bb = 1;
R (branched or linear) = Man, oligomannose or
polysaccharide.

- ↓
1. Endo-H
2. Galactosyltransferase, UDP-Gal
3.. CMP-SA-PEG, ST3Gal3

a-z, bb = 0; aa = 1; R' = -Gal-Sia-PEG.

FIG. 30D

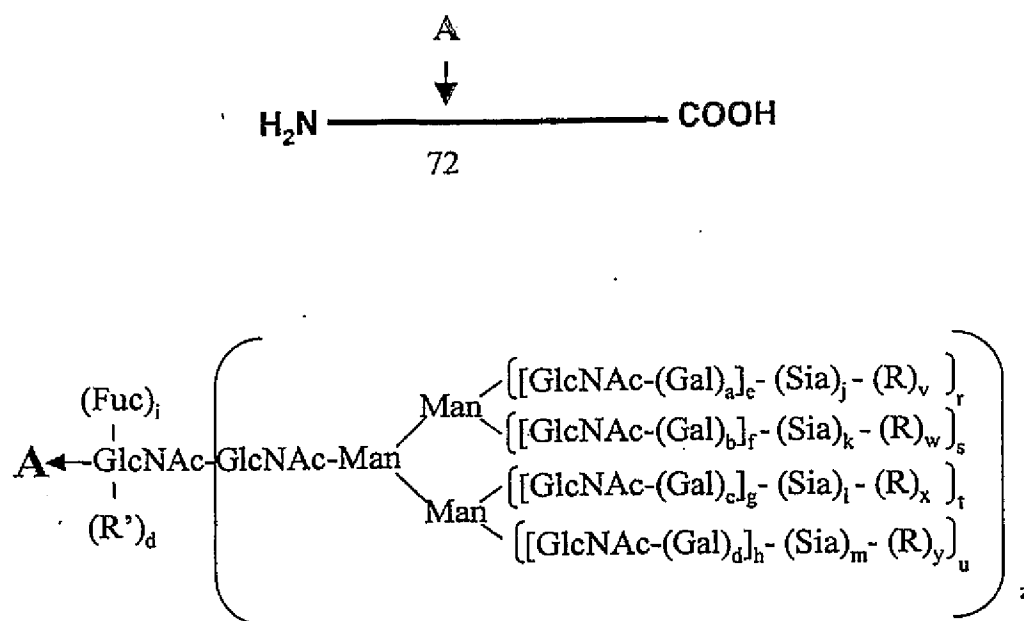67/497



a-d, i, r-u (independently selected) = 0 or 1.

e-h (independently selected) = 0 to 4.

j-m (independently selected) = 0 or 1.

n, v-y = 0; z = 0 or 1.

R = polymer; R' = sugar, glycoconjugate.

FIG. 30E

68/497

CHO, BHK, 293 cells, Vero expressed
interferon alpha-14C.

$h = 1$ to 3 ;

$a-g, j-m, i$ (independently selected) = 0 or 1 ;

$r-u$ (independently selected) = 0 or 1 ;

$n, v-y = 0$; $z = 1$.



1. CMP-SA-PEG, ST3Gal3

$h = 1$ to 3 ;

$a-g, i$ (independently selected) = 0 or 1 ;

$r-u$ (independently selected) = 0 or 1 ;

$j-m, v-y$ (independently selected) = 0 or 1 ;

$z = 1$; $n = 0$; $R = \text{PEG}$.

FIG. 30F

Insect cell or fungi expressed
interferon alpha-14C.

$a-d, f, h, j-n, s, u, v-y = 0$;

e, g, i, r, t (independently selected) = 0 or 1 ;

$z = 1$.



1. GNT's 1,2,4,5, UDP-GlcNAc

2. Galactosyltransferase, UDP-Gal

3. CMP-SA-PEG, ST3Gal3

$a-m, r-y$ (independently selected) = 0 or 1 ;

$z = 1$; $n = 0$; $R = \text{PEG}$.

FIG. 30G

69/497

Yeast expressed interferon alpha-14C.

a-n = 0; r-y (independently selected) = 0 to 1;
z = 1; R (branched or linear) = Man,
oligomannose.

1. mannosidases
2. GNT's 1,2,4,5, UDP-GlcNAc
3. Galactosyltransferase, UDP-Gal
4. CMP-SA-PEG, ST3Gal3

a-m, r-y (independently selected) = 0 or 1;
z = 1; n = 0; R = PEG.

FIG. 30H

NSO expressed interferon alpha 14C.

a-i, r-u (independently selected) = 0 or 1;
j-m, n, v-y = 0; z = 1.

1. CMP-SA-levulinate, ST3Gal3,
buffer, salt
2. H₄N₂-PEG

a-i, j-m, r-y (independently selected) = 0 or 1;
n = 0; z = 1; R = PEG.

FIG. 30I